

Phosphotyrosine HTRF® Assay

Catalog No. 17-10000

384-well

FOR RESEARCH USE ONLY Not for use in diagnostic procedures.

USA & Canada Phone: +1(800) 645-5476 In Europe, please contact Customer Service:

France: 0825.045.645; Spain: 901.516.645 Option 1 Germany: 01805.045.645

Italy: 848.845.645

United Kingdom: 0870.900.46.45

For other locations across the world please visit www.millipore.com/offices

Introduction

The phosphotyrosine HTRF[®] Assay (Cat. # 17-10000) is a generic homogenous method for measuring tyrosine kinase activities using the industry standard monoclonal phosphotyrosine antibody, 4G10[®], tyrosine peptide substrate and a universal detection system. An increase in the HTRF[®] signal is observed as a result of the kinase activity. Inhibitor compounds can be screened during the kinase reaction for an IC 50 value determination.


The assay has a wide dynamic range and high precision, making this assay a valuable new tool for the study of tyrosine phosphorylation.

Background

Cellular Signal Transduction through Tyrosine Phosphorylation Protein phosphorylation is often necessary for the initiation of biological processes such as cell growth, proliferation, and ubiquitin-mediated degradation. Tyrosine phosphorylation, in particular, often mediates receptor signaling events at the membrane by initiating protein-protein interactions in response to growth factors, hormones and cytokines. Research tools that detect tyrosine phosphorylation are central in understanding the significance of cellular signaling in many biological processes and disease origination. Millipore developed the Gold-Standard antibody (4G10) for Tyrosine Phosphorylation.

Tyrosine phosphorylation is one of the most common methods of research. There are many different assay formats available including fluorescence, ELISA, and IP. A number of homogeneous technologies (e.g. luminescence) are restricted by a reduction in signal intensity following assay measurement or prolonged incubations. Our patented HTRF® technology has many benefits including long signal stability. HTRF® fluorescence is not dampened by assay measurement, assay additives such as DMSO or extended incubations prior to reading. The advantages of the tremendous stability of HTRF® fluorescence include: time flexibility and assay security. In case of instrument failure, it is possible to read the microplates even after an extended period of time.

Assay Principle Flow Chart

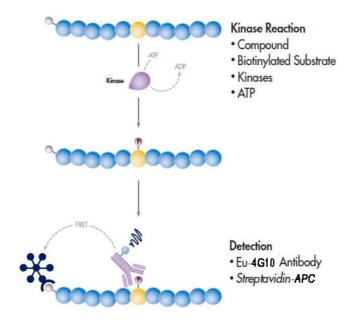

For Research Use Only; Not for use in diagnostic procedures

Kit Components

Store at 2-8° C

- 1. Black 384 well plate: (Part No. R5215) One 384-well black plate.
- 2. <u>Detection buffer</u>: (Part No. CS206478) Two bottles containing 5 mL each of detection buffer.
- 3. Plate Covers: Two plate covers.

Store at -20° C

- 4. <u>10X KinEASE buffer</u>: (Part No. CS206508). Two vials containing 1 mL each of kinase reaction buffer.
- 5. <u>Biotinylated TK1 Peptide</u>: (Part No. CS206504). Two vials containing 60 μ L each of Biotinylated TK1 Peptide Stock at 50 μ M.
- 6. Biotinylated TK2 Peptide: (Part No. CS206505). Two vials containing 60 μ L each of Biotinylated TK2 Peptide Stock at 50 μ M.
- 7. <u>Biotinylated TK3 Peptide</u>: (Part No. CS206506). Two vials containing 60 μ L each of Biotinylated TK3 Peptide Stock at 50 μ M.
- 8. <u>Biotinylated Phosphotyrosine Peptide Standard</u>: (Part No. 90192). One vial containing 10 μg lyophilized biotinylated phosphotyrosine peptide standard.
- 9. <u>Streptavidin-APC</u>: (Part No. CS206501). Two vials containing 360 μL each of streptavidin-APC stock at 1 μM.
- 10. <u>Eu-cryptate-conjugated 4G10[®] anti-phosphotyrosine</u>: (Part No. CS206507). Two vials containing 30 μL each Eu-cryptate conjugated 4G10 anti-phosphotyrosine antibody stock at 35 μg/mL.
- 11. ATP: (Part No. CS206509). Two vials containing 50 μL each of ATP stock at 10mM.
- 12. DTT: (Part No, 90499). Two vials containing 100 µL each of DTT stock at 1M.
- 13. MnCl₂: (Part No. CS206510). Two vials containing 100 μL each of MnCl₂ stock at 0.1M.
- 14. MgCl₂: (Part No. CS206511). Two vials containing 100 μL each of MgCl₂ stock at 1M.
- 15. <u>Sodium Orthovanadate:</u> (Part No. CS206465). Two vials containing 50 μL each of sodium orthovanadate stock at 50 mM.
- 16. <u>Staurosporine:</u> (Part No. CS201778). One vial containing 100 μL of staurosporine stock at 2.5 mM.

Materials Not Supplied

- 1. Multi-channel or repeating pipettes
- 2. Plate shaker (optional)
- 3. Pipettors and tips capable of accurately measuring 1-1000 μL

- 4. Graduated serological pipettes
- 5. 384-well microtiter plate reader with luminescence readout.
- 6. Graphing software for plotting data or graph paper for manual plotting of data
- 7. Microcentrifuge tubes for standard and sample dilutions
- 8. Mechanical vortex
- 9. 1 liter container
- 10. Distilled or deionized water
- 11. DMSO to dissolve inhibitor compounds.

Storage

Maintain the unopened kit at -20°C until expiration date. Once opened, please follow component storage instructions. Avoid repeated freeze/thaw cycles, aliquot if necessary.

Precautions

- The instructions provided have been designed to optimize the kit's performance. Deviation from the instructions may result in suboptimal performance of the kit and the failure to produce accurate data.
- Safety Warnings and Precautions: This kit is designed for research use only and not recommended for internal use in humans or animals. All chemicals should be considered potentially hazardous and principles of good laboratory practice should be followed.

Technical Notes

- All kit reagents should be brought to room temperature (20°C to 25°C) just prior to use.
- Do not use reagents beyond the expiration date of the kit.
- Do not mix or interchange reagent from various kit lots.
- Manual Plate Washing: Vigorous washing and complete removal of all liquid by aspiration at the end of each washing step is very important to obtain low background values. Gentle agitation during the wash steps or a 2-3 minutes soak may reduce background values.

Reagents Preparation

Prepare 1X working kinase buffer. Provided kinase buffer is 10X stock (Part No. CS206508).
 Prepare the required amount per assay setup just before use. Working solutions cannot be stored and must be used immediately, discard any unused.

For example, follow table 1 below to prepare 1 mL of the 1X working kinase buffer.

Table 1.

| Component | Part No. | Amount (µL) | Final concentration |
|------------------------------|----------|-------------|---------------------|
| 10X KinEASE Buffer | CS206508 | 100 | 1X |
| 1M MgCl ₂ | CS206511 | 5 | 5 mM |
| 0.1M MnCl ₂ | CS206510 | 10 | 1 mM |
| 1M DTT | 90499 | 1 | 1 mM |
| 50mM Sodium Orthovanadate | CS206465 | 2 | 0.1 mM |
| Distilled Water | N/A | 882 | N/A |

*Note: The appendix table A listed the recommended kinase reaction conditions. Please reference the recommended 1X reaction buffer column for the optimized buffer. MnCl₂ may not be needed for some kinases. Please adjust the preparation of 1X buffer accordingly.

- 2. Prepare working inhibitor compounds solution. Depends on the compound property, the end user will need to determine the suitable solvent to dissolve it to a concentrated stock. Typically, DMSO is a commonly used solvent. Then make further dilutions with 1X kinase buffer to give the final DMSO no more than 2% in the reaction to minimize any potential effects.
- 3. Prepare working biotinylated tyrosine kinase peptide solution as substrate. Total of 3 biotinylated tyrosine kinase (TK) peptides are included with the kit. All at 50 µM stock. Reference the appendix table attached in the end of this manual guide for a suitable TK substrate to the choice of kinase. Follow the optimization steps below to determine a proper final concentration for each kinase target. Then prepare a 5X concentrated working solution with 1X kinase buffer according to the determined proper final usage if inhibitor compounds screening will be applied to your target kinases.
- 4. Prepare working solution of biotinylated phosphotyrosine standard. Add 109.6 μ L Distilled Water to the component 90192 to dissolve completely. This is a 50 μ M solution.

Optimization of the Kinase Assay

Fix Biotin-substrate/Streptavidin-APC molar ratio to 8/1

Kinase titration

This allows the optimal kinase concentration for which the signal reaches around $80{\sim}90\%$ of the maximum to be determined. A fixed concentration of the biotin-TK substrate (1 μ M) and ATP (100 μ M) should be tested with a range of kinase

concentrations (0.5-20ng/reaction in a two fold serial dilutions). Setup the enzymatic reaction to run for 60 minutes at 30-37°C.

Kinetic study

Enzyme kinetic depends on the kinase and substrate concentrations. A time course is performed using the optimal kinase concentration determined at the kinase titration step, ATP (100 μ M) and biotin-TK substrate (1 μ M). The reaction is stopped at various end points by the addition of the detection reagents (1, 2, 5, 10, 15, 30, 60, 90 minutes). Determine the optimal reaction time when the signal reaches around 80-90% of the maximum. This optimal reaction time is used for the rest experiments.

Substrate titration

This allows the determination of substrate Km. Use the optimal kinase concentration determined from kinase titration step, fix concentration of ATP (100 μ M) and test a range of the biotin-TK substrate (0.03-1 μ M in a two fold serial dilutions). Adjust the streptavidin-APC concentration to have a fixed biotin/streptavidin molar ratio at 8/1. Reference to table 2 for making the adjustment. Let the kinase reaction run for 60 minutes at 30-37°C.

ATP titration

This allows the determination of ATP Km. Use the optimal kinase concentration, a fixed concentration of the biotin-TK substrate (1 μ M) to test a range of ATP concentrations (3.125-200 μ M in a two fold serial dilutions). The kinase reaction is stopped at 60 minutes at 30-37°C.

Inhibitor compounds IC 50 determination

The kinase activity is tested over a broad range of inhibitor concentrations to generate a dose response curve. The test is generally run using the previously determined optimal assay conditions. A proper ATP concentration (reference the Km value determined at the ATP titration step) should be considered for ATP-competitive compounds.

Suggested final concentration ranges of the assay components:

| Component | Concentration range | |
|---------------------|---------------------|--|
| Biotin-TK substrate | 0.03-1 μM | |
| Kinase | 0.5-20 ng/reaction | |
| ATP | 3.125-200 µM | |
| Streptavidin-APC | 0.04-125 nM | |
| Eu-4G10 | Ready to use | |

Assay Protocol


I. Setup a standard curve with provided biotinylated phosphotyrosine peptide.

Biotinylated Phosphotyrosine Peptide Standard serial dilutions should be setup along with various concentration of streptavidin-APC to keep Biotin/Streptavidin ratio fix at 8/1.

Use 3 columns to setup the standard in triplicates for each point. Cover the unused remaining wells with plate sealer.

Wells Layout:

STD (uM)

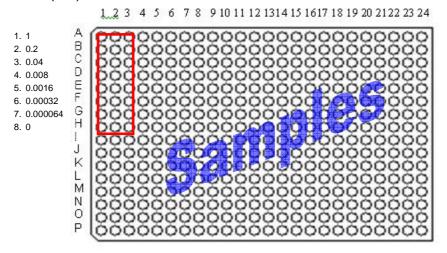

1. Label 7 test tubes #1-7 and "0 dose". Add 48 μ L of 1X kinase buffer to tubes #2-7 and "0 dose".

The biotinylated phosphotyrosine peptide standard solution prepared above is further used to make the #1 standard tube. It will be a 1 μ M final in the assay.

2.4 μ L of 50 μ M stock 15 μ L of streptavidin-APC (part No. CS206501) 42.6 μ L of Detection Buffer (part No. CS206478) Total 60 μ L.

- 2. Standards #2-7 are then prepared by performing a 5-fold serial dilution of the preceding standard. For example, to make Standard #2, remove 12 μ L of Standard #1 and add it to tube #2 and vortex, and so on. Do not add any Standard to the "0 dose" tube.
- 3. Table 2 below showed the final biotinylated-phosphotyrosine (as biotin-p-Tyr) and Streptavidin-APC (as Strep-APC) concentrations in the standard solutions.

Table 2.

| Biotin-p-Tyr | Standard (µM) | Streptavidin- | APC (nM) |
|--------------------|---------------|--------------------|---------------|
| Conc. In standard | Final assay | Conc. in standard | Final assay |
| solution (12.5 µL) | conc. (25 µL) | solution (12.5 µL) | conc. (25 μL) |
| 2 | 1 | 250 | 125 |
| 0.4 | 0.2 | 50 | 25 |
| 0.08 | 0.04 | 10 | 5 |
| 0.016 | 0.008 | 2 | 1 |
| 0.0032 | 0.0016 | 0.4 | 0.2 |
| 0.00064 | 0.00032 | 0.08 | 0.04 |
| 0.000128 | 0.000064 | 0.016 | 0.008 |
| 0 | 0 | 0 | 0 |

- 4. Dispense 12.5 μL/well of the standard solutions made above (reference table 3 below for the biotin-p-Tyr column.
- II. Setup kinase activity samples, compounds can be added to the kinase reactions in a dose-response curve to determine IC 50 value.
- 5. Setup kinase reactions by applying the optimal condition determined during the assay optimization. Reference table 3 for reaction setup examples.

<u>Note*.</u> It is recommended to have a buffer only and Eu-4G10 only control for the first time the assay is setup (as highlighted columns in table 3).

Table 3.

| Kinase Reaction | Negative | Positive | Inhibitor/ Compounds | Biotin-p-Tyr Standard | Eu-4G10 only | Buffer only |
|-----------------------------|----------|----------|-------------------------|--------------------------|-----------------|----------------|
| Inhibitor/Compounds | - | - | 2.5μL | - | - | - |
| Biotin TK Substrate (5µM) | 2.5μL | 2.5μL | 2.5μL | - | - | - |
| Kinase (4 ng/μL) | - | 2.5μL | 2.5µL | - | - | - |
| ATP (500μM) | 2.5μL | 2.5µL | 2.5μL | - | - | - |
| Biotin-p-Tyr Std as table 2 | - | - | - | 12.5μL | - | - |
| 1X Kinase Buffer | 7.5µL | 5μL | 2.5μL | - | 12.5μL | 12.5μL |

- 6. A total 12.5 µL reaction mix per well. Gently tap the plate to mix the solutions, cover with one plate sealer. Incubate for 60 minutes at 30-37°C on a plate shaker with gentle agitation.
- 7. Prepare the detection reagent:
 - 1). For biotin-p-Tyr and Eu-4G10 only wells, dilute Eu-4G10 conjugate (part No. CS206507) 1:100 with detection buffer (part No. CS206478), vortex briefly to mix well. Make sure to prepare some extra for pipetting.

- 2). For buffer only wells, directly add 12.5 μ L/well detection buffer (part No. CS206501).
- 3). For negative, positive and inhibitor/compounds wells, prepare the detection reagents as below, Make sure to prepare some extra for pipetting.

Eu-4G10 1:100

Streptavidin-APC 1:8 if 1 µM Biotin-TK substrate is applied in

the kinase reaction.

DM C Fill up the required volume.

- * Adjust the dilution factor of streptavidin-APC to have a fix biotin/streptavidin molar ratio at 8/1 if a different concentration of biotin-TK substrate is applied.
- 8. Pipet 12.5 µL detection reagent prepared in step 7 to each well accordingly. Gently tap to mix the solutions. Incubate for 60 minutes at room temperature on a plate shaker with gentle agitation. Protect from light.
- 9. Measure TR-FRET ratio on an appropriate reader according to the following parameters (these are guideline parameters based on Molecular Probes -Analyst, please also refer to parameters recommended in the instrument instruction manual):

Excitation 330 nm

Emission 660-50 nm and 620-35nm

HTRF Ratios are calculated as follows:

$$HTRF \ Ratio = \left(\frac{Emission \ at \ 660nm}{Emission \ at \ 620nm}\right) \times 10000$$

10. Data is interpreted by GraphPad PRISM nonlinear fit, Sigmoidal dose-response (variable slope).

Assay Results

Z' factor was calculated to be greater than 0.8 using 48 replicates for each of no biotin-p-Tyr detection mixture as background signal and biotin-p-Tyr detection mixture as positive signal.

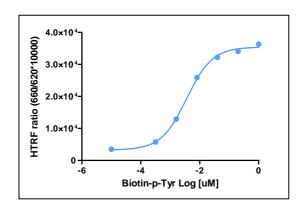


Figure 1. Dose response curve of Biotinylated-phosphotyrosine peptide.

Phosphotyrosine-containing peptide from $0\sim1~\mu\text{M}$ was incubated with Europium 4G10 and streptavidin APC. TR-FRET signal was determined, and a non linear curve fit with Sigmoidal Dose Response was applied.

NOTE: This data is presented for reference use only and should not be used to interpret actual assay results. A standard curve must be generated for each assay.

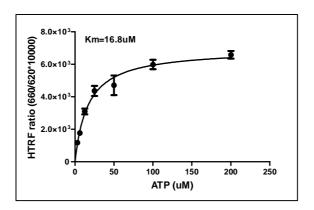


Figure 2. Determination of ATP dependence of active Src kinase by TR-FRET.10ng active Src kinase was applied with Biotin-TK3 substrate. The Km value may be used for

10ng active Src kinase was applied with Biotin-1K3 substrate. The Km value may be used to establishing ATP concentration when screening compounds.

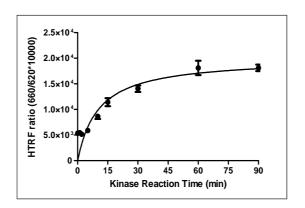


Figure 3. Kinetics of substrate phosphorylation by active FGFR1 kinase domain. The kinase reaction was complete by 60 minutes at 30°C.

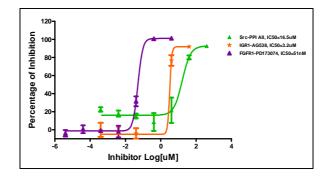
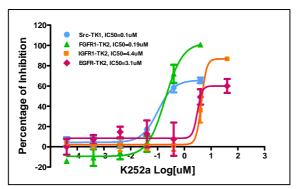


Figure 4. Determination of IC50 values for kinase-selective inhibitors.

The assay contained 10ng kinase, 100 μ M ATP, and 0.5 μ M biotinylated substrate, was incubated for 60 minutes at 30°C. Data obtained by TR-FRET was analyzed by non-linear curve fit with sigmoidal dose response with variable slope.



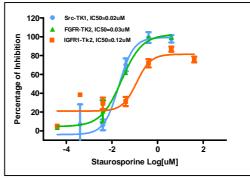


Figure 5. Determination of IC50 values of universal ATP-competitive inhibitors with a panel of kinases.

The assay contained 10ng kinase, 20 μ M ATP, and 1 μ M biotinylated substrate, and was incubated for 60 minutes at 30°C. Data obtained by TR-FRET was analyzed by non-linear curve fit with sigmoidal dose response with variable slope.

* K252a (alkaloid), a staurosporine analog, is a highly potent cell permeable inhibitor of CaM kinase and phosphorylase kinase. It acts by competition with the ATP binding site, non-selective protein kinase inhibitor.

References

- 1. Kanakura, Y., et al. (1991). J. Biol. Chem. 266:490.
- 2. Cohen, B., et al. (1990). Proc. Natl. Acad. Sci.USA. 87: 4458.
- 3. Druker, B. J., et al. (1989). New Eng. J. Med. 321: 1383.
- 4. O'Brian, C A and Ward, N E (1990). J Natl Cancer Inst 82: 1734-5.
- 5. Raychaudhuri, S.P., et al., *J. Invest. Dermatol.*, 122, 812-819 (2004).
- 6. 8. Winston, J.H., et al., J. Pain., 4, 329-337 (2003).

Appendix A. Recommended Kinase Conditions.

Table of recommended reaction conditions for Upstate kinases (h = human; m = mouse):

| Kinase | Millipore Part | Recommended Substrate | Recommended Kinase Reaction Condition |
|-----------------|-------------------|--------------------------|---------------------------------------|
| Abl (h) | 14-529 | Biotin-TK1 | 5mM MgCl ₂ |
| Abl (T315l) (h) | 14-522 | Biotin-TK1 | 5mM MgCl ₂ |
| Abl (m) | 14-459 | Biotin-TK1 | 5mM MgCl ₂ |
| Arg (h) | 14-521 | Biotin-TK1 | 5mM MgCl ₂ |

| Kinase | Millipore Part | Recommended Substrate | Recommended Kinase Reaction Condition |
|---------------|-------------------|--------------------------|---|
| Arg (m) | 14-460 | Biotin-TK1 | 5mM MgCl ₂ |
| Blk (m) | 14-316 | Biotin-TK1 | 5mM MgCl ₂ |
| Bmx (h) | 14-499 | Biotin-TK1 | 5mM MgCl ₂ |
| BTK (h) | 14-552 | Biotin-TK1 | 5mM MgCl ₂ |
| CSK (h) | 14-458 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Src (h) | 14-326 | Biotin-TK3 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| EGFR (h) | 14-531 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| EphA2 (h) | 14-560 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| EphB2 (h) | 14-553 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| EphB4 (h) | 14-554 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Fes/Fps (h) | 14-473 | Biotin-TK1 | 5mM MgCl ₂ |
| FGFR1 (h) | 14-582 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| FGFR4 (h) | 14-583 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| FGFR3 (h) | 14-464 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Flt1 (h) | 14-562 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Flt3 (h) | 14-500 | Biotin-TK1 | 5mM MgCl ₂ |
| Fms (h) | 14-551 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Fyn (h) | 14-441 | Biotin-TK3 | 5mM MgCl ₂ |
| IGF-1R (h) | 14-465 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Insulin R (h) | 14-466 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Lck (h) | 14-442 | Biotin-TK3 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Lyn (h) | 14-510 | Biotin-TK1 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Lyn (m) | 14-315 | Biotin-TK1 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Met (h) | 14-526 | Biotin-TK1 | 5mM MgCl ₂ |
| PDGFRα (h) | 14-467 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| PDGFRβ (h) | 14-463 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Ret (h) | 14-570 | Biotin-TK1 | 5mM MgCl ₂ |
| Ros (h) | 14-527 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Rse (h) | 14-535 | Biotin-TK3 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Syk (h) | 14-314 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |

| Kinase | Millipore Part | Recommended Substrate | Recommended Kinase Reaction Condition |
|------------|-------------------|--------------------------|---|
| Tie2 (h) | 14-540 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Trk A (h) | 14-571 | Biotin-TK1 | 5mM MgCl ₂ , 1mM MnCl ₂ |
| Trk B (h) | 14-507 | Biotin-TK2 | 5mM MgCl ₂ |
| Yes (h) | 14-478 | Biotin-TK2 | 5mM MgCl ₂ |
| ZAP-70 (h) | 14-404 | Biotin-TK2 | 5mM MgCl ₂ , 1mM MnCl ₂ |

Troubleshooting Guide

| Problem | Potential Cause | Experimental Suggestions |
|---|--|--|
| No signal or weak signal in all wells | Missing components or key steps | Check to make sure all components were added in the appropriate steps and amounts. |
| | Biotinylated-phosphotyrosine Peptide is no longer active or has reduced activity | Make sure all components are stored at the recommended temperature and minimize the freeze/thaw cycle as manual recommends. Make aliquots of components when first thawed if planning more than one assay. |
| | Plate reader or settings are not optimal | Verify the measurement, read time, and filter on the plate reader. |
| | Incorrect storage temperatures Incorrect assay temperature | Items are to be stored at the appropriate storage temperatures. Performance can be negatively affected if reagents are not stored and used in the appropriate time period. |
| No detectable signal in samples | Low phosphotyrosine level or missing key components. | Check the kinase reaction setup procedure, make sure ATP, biotin-TK substrate and MgCl ₂ , MnCl ₂ are added to the reaction buffer. |

Warranty

Millipore Corporation ("Millipore") warrants its products will meet their applicable published specifications when used in accordance with their applicable instructions for a period of one year from shipment of the products. MILLIPORE MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. The warranty provided herein and the data, specifications and descriptions of Millipore products appearing in Millipore's published catalogues and product literature may not be altered except by express written agreement signed by an officer of Millipore. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorized and if given, should not be relied upon.

In the event of a breach of the foregoing warranty, Millipore's sole obligation shall be to repair or replace, at its option, the applicable product or part thereof, provided the customer notifies Millipore promptly of any such breach. If after exercising reasonable efforts, Millipore is unable to repair or replace the product or part, then Millipore shall refund to the Company all monies paid for such applicable Product. MILLIPORE SHALL NOT BE LIABLE FOR CONSEQUENTIAL, INCIDENTAL, SPECIAL OR ANY OTHER DAMAGES RESULTING FROM ECONOMIC LOSS OR PROPERTY DAMAGE SUSTAINED BY ANY COMPANY CUSTOMER FROM THE USE OF ITS PRODUCTS.

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.

(c) 2010: Millipore Corporation. All rights reserved. No part of these works may be reproduced in any form without permission in writing