

sigma-aldrich.com

3050 Spruce Street, St. Louis, MO 63103 USA Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757 email: techservice@sial.com sigma-aldrich.com

Product Information

Membrane Scaffold Protein 1D1 BTN recombinant, expressed in *E. coli* buffered aqueous solution

Catalog Number **MSP13** Storage Temperature –20 °C

Synonyms: MSP1D1 BTN

Product Description

Nanodisc technology is a widely applicable approach to render membrane proteins soluble in aqueous solutions in a native-like bilayer environment, where the membrane proteins remain stable and active.

The Nanodisc concept is derived from high density lipoprotein (HDL) particles and their primary protein component, apolipoprotein. The Nanodisc is a non-covalent structure of phospholipid bilayer and membrane scaffold protein (MSP), a genetically engineered protein, which mimics the function of Apolipoprotein A-1 (ApoA-1).^{1,2} A soluble Nanodisc assembles as the phospholipid forms a bilayer, which is encircled by two amphipathic MSP molecules that cover the hydrophobic alkyl chains of the bilayer. The length of the MSP controls the size of the Nanodisc structure.³


Incorporation of a membrane protein target into a Nanodisc is accomplished by first solubilizing the target protein in detergent, and then mixing with cholatesolubilized phospholipid and MSP. Removal of the detergents results in the assembly of the Nanodisc structure around the target protein (see Figure 1).

The resulting structure renders the protein soluble in a model membrane system of defined phospholipid content and size with a native-like bilayer. This provides stability and functional requirements for the protein target. The size of the Nanodisc, determined by the length of the MSP, allows control of the oligomerization state of the membrane protein.

Some examples of membrane-bound proteins incorporated into Nanodiscs are GPCRs,⁴ cytochromes P450,⁵⁻⁷ bacteriorhodopsin,⁸ coagulation factors,⁹ cholera toxin,¹⁰ and TAR receptor.¹¹

Figure 1.

Example of a Nanodisc containing a 7-transmembrane protein

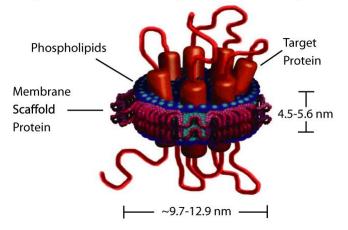

A critical component of Nanodiscs is the membrane scaffold protein (MSP), which forms the encircling amphipathic helical protein belt.^{12,13} The sequence similarities of MSPs to Apolipoprotein A-1 are presented in Figure 2.

Figure 2.

Protein Maps of Apolipoprotein A-1 and MSP1D1 Biotin-Labeled

		A	oolipo	proteir	ח A1				
Globular Re	gion H1	H2	Нз	H4 H	45 H	5 H7	H8	H9 H1	0
1	44	66	88 99	121	143	165	188 20	9 220	243
TEV pro	otease tion site	Si	gma M	SP13 =	= MSP1	D1 Biot	tin-Labe	eled	
Hist. tag	1	H2 H	3 H4	H5	H6	H7 I	H8 H9	H10 p	-Terminal Biotin
1 9	23 34	56	67	89 1	11 13	3 155	177 18	88 212	231

The first MSP, MSP1, was engineered with its sequence based on the sequence of A-1, but without the globular N-terminal domain of native A-1.¹⁴ The Membrane Scaffold Protein 1D1 (MSP1D1) variant of MSP1 deletes the first 11 amino acids in the Helix 1 portion (referred to as "H0.5" in Figure 2) of the original MSP1 sequence.³ The MSP1D1 protein is an N-terminal histidine-tagged protein with a TEV protease cleavage site between the histidine-tag and the protein sequence. The MSP1D1 BTN variant of MSP1D1 features an enzymatically biotinylated additional 20-amino acid sequence at the C-terminus of MSP1D1.

MSP1D1 BTN generates Nanodiscs that typically contain 120–160 lipid molecules and two MSPs per Nanodisc.³ The thickness of a Nanodisc is dependent on the type of phospholipid incorporated (typically 4.6–5.6 nm).³

This product is supplied as a buffered aqueous solution, at a protein concentration of ~5 mg/mL, and containing 20 mM Tris-HCl, pH 7.4, 100 mM NaCl, 0.5 mM EDTA, 5 mM sodium cholate, and 0.1% sodium azide.

MSP1D1 BTN sequence:¹⁵

GHHHHHHHDYDIPTTENLYFQGSTFSKLREQLGPVT QEFWDNLEKETEGLRQEMSKDLEEVKAKVQPYLDDF QKKWQEEMELYRQKVEPLRAELQEGARQKLHELQE KLSPLGEEMRDRARAHVDALRTHLAPYSDELRQRLA ARLEALKENGGARLAEYHAKATEHLSTLSEKAKPALE DLRQGLLPVLESFKVSFLSALEEYTKKLNTQSLLGGG LNDIFEAQKIEWHE + Biotin

Molecular mass: 27.123 Da (based on sequence, calculated with biotin)

Molecular mass: 26,900.1 Da (based on sequence,

calculated without biotin)

Molecular mass of histidine-tag cleaved variant: 24.282.4 (based on sequence, calculated without biotin)

The MSP1D1 BTN protein concentration can be determined spectrophotometrically using the following extinction coefficients:

MSP1D1, $\epsilon_{280} = 26,930 \text{ M}^{-1} \text{ cm}^{-1}$

Histidine-tag cleaved variant, $\varepsilon_{280} = 23,950 \text{ M}^{-1} \text{ cm}^{-1}$ (Solvent: 20 mM Tris, pH 7.4, with 0.1 M NaCl, 0.5 mM EDTA, and 0.01% NaN₃)

Reagents Required, but Not Provided

- Sodium Cholate (Catalog No. C6445)
- Amberlite[®] XAD[®]-2 (Catalog No. 20275)
- Phospholipid(s)
 - 2-Oleoyl-1-palmitoyl-*sn*-glycero-3-phosphocholine (POPC, Catalog No. P3017)
 - 1,2-Dipalmitoyl-*rac*-glycero-3-phosphocholine (DPPC, Catalog No. P5911)
 - 1,2-Dimyristoyl-*rac*-glycero-3-phosphocholine (DMPC, Catalog No. P7930)
- Triton[™] X-100 (optional, Catalog No. X100)
- Superdex[®] 200 (Catalog No. S6782) or

Superdex 200 Increase 10/300 GL column (Catalog No. GE28-9909-44)

Precautions and Disclaimer

This product is for R&D use only, not for drug, household, or other uses. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

Store this product at –20 °C. Working stock solutions of MSP1D1 can be stored at 2–8 °C for up to several days. 16

Procedure

The successful assembly of Nanodiscs is highly dependent on the properties of the target protein. The choice of the types and ratios of MSPs, phospholipids, and detergents might require trials and titrations to optimize disc assembly. A thorough review of the literature is essential for this aspect of Nanodisc construction.^{3,12,13,14}

Successful incorporation of the target membrane protein into Nanodiscs most likely will require optimization of the following important parameters:

- Choice and concentration of secondary detergent for preparation of target membrane protein
- Choice and preparation of phospholipid
- Temperature of reconstitution
- Lipid to MSP ratio
- MSP to target protein ratio

References

- Brouillette, C.G., et al., Biochim. Biophys. Acta, 1531(1-2), 4-46 (2001).
- 2. Jonas, A., *Methods Enzymol.*, **128**, 553-582 (1986).
- 3. Denisov, I.G., *et al.*, *J. Am. Chem. Soc.*, **126**(11), 3477–3487 (2004).
- Bayburt, T.H., et al., J. Biol. Chem., 282(20), 14875-14881 (2007).
- Denisov, I.G., et al., J. Biol. Chem., 282(10), 7066-7076 (2007).
- Grinkova, Y.V., et al., Biochem. Biophys. Res. Commun., 372(2), 379-382 (2008).
- 7. Denisov, I.G., *et al.*, *Biochim Biophys Acta*, **1814(1)**, 223-229 (2011).
- Bayburt, T.H., et al., Arch. Biochem. Biophys., 450(2), 215-222 (2006).
- Shaw, A.W., et al., J. Biol. Chem., 282(9), 6556-6563 (2007).
- 10. Borch, J., *et al.*, *Anal. Chem.*, **80(16)**, 6245-6252 (2008).
- Boldog, T., et al., Proc. Nat. Acad. Sci. USA, 103(31), 11509-11514 (2006).
- Nath, A., et al., Biochemistry, 46(8), 2059-2069 (2007).
- 13. Bayburt, T.H., and Sligar, S.G., *FEBS Lett.*, **584(9)**,1721-1727 (2010).
- 14. Bayburt, T.H., *et al., Nano Letters*, **2(8)**, 853-856 (2002).
- 15. Supplier information.

For a further and extensive list of citations and protocols, visit the website of Professor Stephen Sligar, University of Illinois:

http://sligarlab.life.uiuc.edu/nanodisc.html

Nanodisc technology, and many of its uses, are covered by the following patents held by the University of Illinois:

- 7,691,414: Membrane scaffold proteins
- 7,662,410: Membrane scaffold proteins and embedded membrane proteins
- 7,622,437: Tissue factor compositions and methods
- 7,592,008: Membrane scaffold proteins
- 7,575,763: Membrane scaffold proteins and tethered membrane proteins
- 7,083,958: Membrane scaffold proteins
- 7,048,949: Membrane scaffold proteins

Triton is a trademark of The Dow Chemical Company or an affiliated company of Dow.

Amberlite and XAD are registered trademarks of The Dow Chemical Company or an affiliated company of Dow.

Superdex is a registered trademark of GE Healthcare.

GCY,TD,RC,RBG,MAM 02/17-1