




# Material Matters

Chemistry Driving Performance

Vol. 1 No. 1

## Polymerization for Advanced Applications

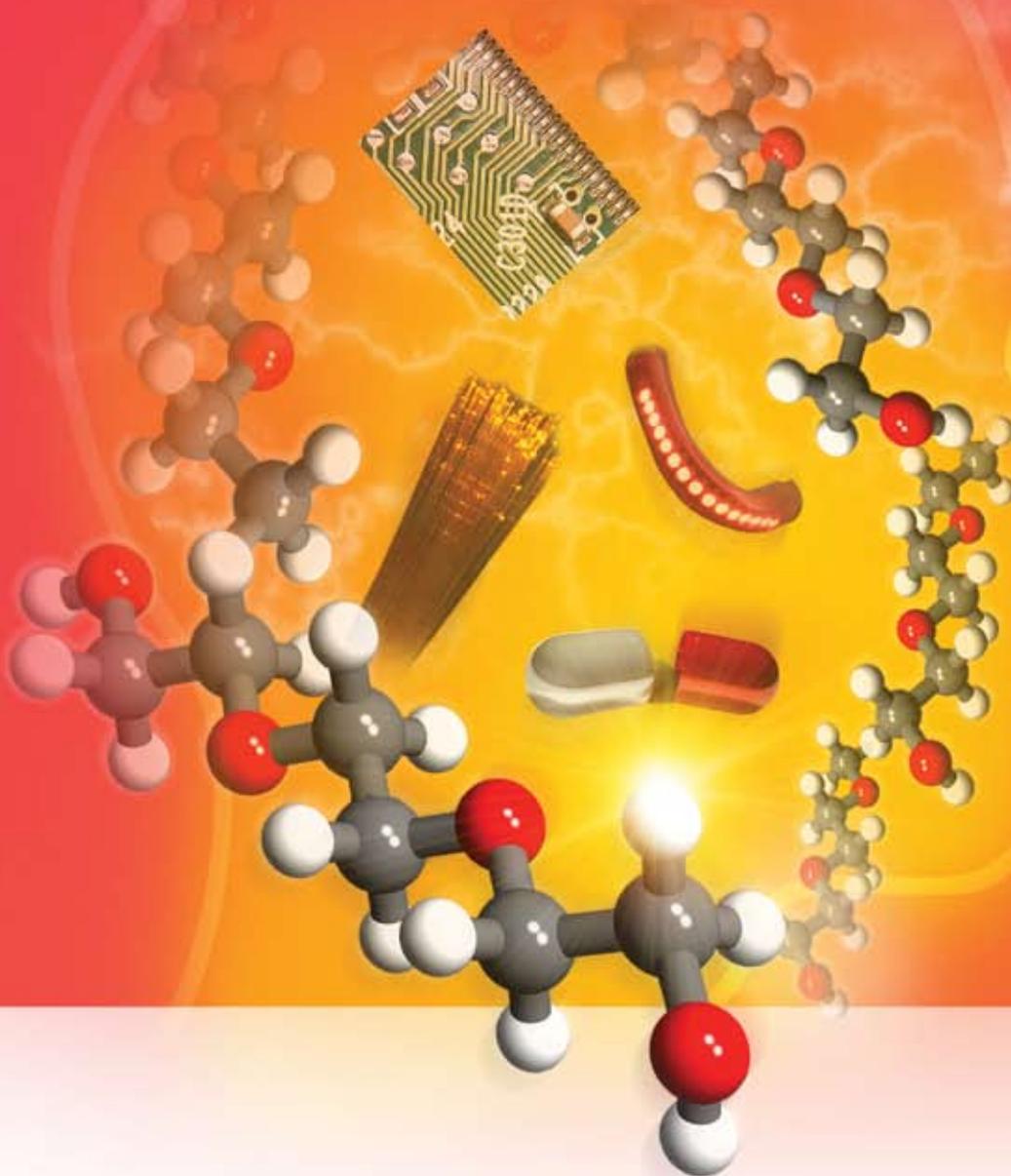

Initiator/Stabilizer FAQs

Polymer Analysis by NMR

Fluorinated Hyperbranched  
Polymers

Etch-Resistant Block  
Copolymers

Bioactive Hydrogels

Anionic Polymerization

Premiere Issue

[sigma-aldrich.com](http://sigma-aldrich.com)

 **ALDRICH**<sup>®</sup>  
Advancing Science

## Introduction

Welcome to the premiere issue of *Material Matters*, a technical guide from Sigma-Aldrich dedicated to addressing research needs in materials science and technology. The theme for this issue is polymerization. Innovative polymers are helping advance almost every field of materials science, from alternate energy to organic electronics. Included in this guide are reviews from researchers in relevant technical fields that discuss a subset of some of the tools available to scientists and engineers.

Our mission at Sigma-Aldrich is to inspire and advance your research. I hope that by highlighting the innovations and products featured in this technical guide, we will help generate the next ideas to significantly impact research in polymerization.

Sincerely,



Luke Grocholl, Ph.D.  
Materials Science Team  
Sigma-Aldrich Corporation



## The MatSci eMailbox

Contact us at [matsci@sial.com](mailto:matsci@sial.com) for:

- Comments on this technical guide
- Suggested topics for future technical guides
- New product suggestions
- Spectral and other data to facilitate your research
- Questions for inclusion in a FAQ
- Subscriptions to future issues of *Material Matters*

## The Polymerization Tools Web Site

Visit us at [sigma-aldrich.com/poly](http://sigma-aldrich.com/poly) for:

- Over 1000 monomers for applications from drug delivery to PLEDs
- Guide to thermal initiators with solvent-specific half-life temperature values
- Absorbance spectra of over 250 photoinitiators
- Functional polymers for the synthesis of advanced copolymers
- Comprehensive list of surfactants organized by HLB value
- Cross-linkers, chain transfer agents, plasticizers, and stabilizers for polymer modification

### About Our Cover

Central to advances in science and engineering are the dedicated researchers whose ideas drive tomorrow's technologies. This premiere issue of *Material Matters* features advanced and unique Sigma-Aldrich materials for polymerization such as monomers, cross-linkers, and functional polymers. The cover depicts how our products combine with your ideas to yield the next generation of nanolithography, fiber optics, flexible LEDs and bioactive polymers for drug delivery. *Material Matters—Chemistry Driving Performance.*

 **Material Matters**  
Chemistry Driving Performance

Vol. 1 No. 1

**Aldrich Chemical Co., Inc.**  
**Sigma-Aldrich Corporation**  
6000 N. Teutonia Ave.  
Milwaukee, WI 53209, USA

### To Place Orders

Telephone 800-325-3010 (USA)  
FAX 800-325-5052 (USA)

### Customer & Technical Services

Customer Inquiries	800-325-3010
Technical Service	800-231-8327
SAFC™	800-244-1173
Custom Synthesis	800-244-1173
Flavors & Fragrances	800-227-4563
International	414-438-3850
24-Hour Emergency	414-438-3850
Web Site	<a href="http://sigma-aldrich.com">sigma-aldrich.com</a>
Email	<a href="mailto:aldrich@sial.com">aldrich@sial.com</a>

### Subscriptions

To request your **FREE** subscription to *Material Matters*, please contact us by:

Phone: 800-325-3010 (USA)

Mail: **Attn: Marketing Communications**  
**Aldrich Chemical Co., Inc.**  
**Sigma-Aldrich Corporation**  
P.O. Box 355  
Milwaukee, WI 53201-9358

Email: [sams-usa@sial.com](mailto:sams-usa@sial.com)

International customers, please contact your local Sigma-Aldrich office. For worldwide contact information, please see back cover.

*Material Matters* is also available in PDF format on the Internet at [sigma-aldrich.com/matsci](http://sigma-aldrich.com/matsci).

Aldrich brand products are sold through Sigma-Aldrich, Inc. Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product for its particular use. See reverse side of invoice or packing slip for additional terms and conditions of sale.

*Material Matters* is a publication of Aldrich Chemical Co., Inc. Aldrich is a member of the Sigma-Aldrich Group. © 2006 Sigma-Aldrich Co.

## Initiator/Stabilizer FAQs



Dr. S. S. Newaz  
Polyorganix Inc.  
Houston, TX

### Q. How does one choose an appropriate initiator?

In a free-radical addition polymerization, the choice of polymerization initiator depends mainly on two factors: a) its solubility and b) its decomposition temperature. If the polymerization is performed in an organic solvent, then the initiator should be soluble in that solvent, and the decomposition temperature of the initiator must be at or below the boiling point of the solvent. Commonly, AIBN (2,2'-Azobis(2-methylpropanitrile)) (**441090**, dec. 102–104°C) and BPO (2-(4-Biphenyl)-5-phenyloxazole) (**216984**, mp 115–119°C) suit these requirements. If the desired polymerization occurs at or below 20 °C, then special, low-temperature free-radical initiators need to be used. Various azo-type initiators can be chosen to satisfy the decomposition temperature requirement.

For emulsion polymerization or polymerization in an aqueous system, a water-soluble initiator like  $K_2S_2O_8$  (**379824**) or an organic, water-soluble initiator (4,4'-Azobis(4-cyanovaleric acid)) (**118168**, dec. 118–125 °C) would be suitable.

### Q. How does one determine the reactivity of a monomer?

Determination of monomer reactivity is not always obvious or straightforward. Researchers rely on their experience and published data on individual monomers. In general, extent of conjugation in the molecular structure can be viewed as indicative of its tendency to form the initial free radical required for propagating a free-radical polymerization. Usually, a more conjugated system is more likely to undergo free-radical polymerization.

### Q. When is it necessary to remove a stabilizer prior to polymerization and how does one do so?

To inhibit polymerization during storage, many monomers are provided with a stabilizer as indicated by the label. Usually, it is not necessary to remove stabilizers. They are typically present in ppm level, and the use of a free radical initiator at the polymerization temperature will overwhelm the effect of the stabilizer. In worst-case scenarios, one may need to add extra amounts of initiator to sustain an acceptable polymerization rate. In most cases, once a polymerization initiates, the rate can be sustained without much difficulty. If, however, it is absolutely necessary to remove the stabilizer, column chromatography is the preferred method (for inhibitor removal columns, see products **306312**, **311332**, **306320**).

### Q. How does one remove residual initiator, stabilizer, and/or unreacted monomer after polymerization?

It is a common practice to dissolve the polymers in a solvent prior to end use, followed by precipitating the polymer using a cosolvent. Usually, the residual initiators and stabilizers will remain in solution and the polymers will separate out as a solid (powder, gum, or fibers). This process may be repeated until desirable polymer characteristics are obtained. This fractional precipitation is also effective in removing lower molecular weight polymers, resulting in narrower molecular weight distribution—of course accompanied with a loss of yield. Typical solvent/cosolvent pairs could be toluene/hexane, toluene/methanol, THF/water, etc., determined by the relative solubilities of the polymer versus the small-molecule component.

For initiator solubilities and decomposition temperatures, visit our Web site at [sigma-aldrich.com/poly](http://sigma-aldrich.com/poly).

## Polymer Analysis by NMR

### Sigma-Aldrich Quality Control Team

One of the challenges polymer scientists face is molecular weight (average chain length) determination of their materials. While membrane osmometry, gel permeation chromatography, viscosity analysis and mass spectrometry are typically used for molecular weight determination, the techniques can be time consuming, inaccurate for the molecular weight ranges involved, or require specialized instrumentation. End-group analysis by NMR offers an easy alternative method using an instrument commonly found in many analytical labs. In addition, NMR analysis can also be used to accurately determine monomer ratios for various copolymer.


Scientists at Sigma-Aldrich routinely determine number-average molecular weight ( $M_n$ ) by  $^1H$  NMR end-group analysis for polymers having  $M_n$  values under 3000. Sensitivity of the instrument to detect end-group protons will determine the upper limit that can be measured. In order to use this method, the following criteria must be met:

- Identifiable end-group protons distinguishable from repeating monomer-group protons by NMR
- Accurate integration of both end-group and monomer protons
- Knowledge of monomer formula weights

Once the ratio of protons on the end-groups to protons on the polymer chain is determined, using the NMR, simple math can be applied to generate the  $M_n$  value.

### This example illustrates this method:

#### 437441 Poly(ethylene glycol) diacrylate

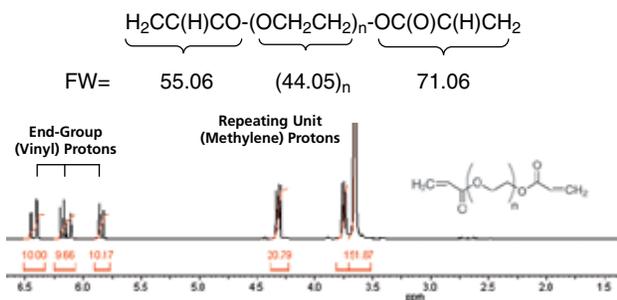

#### 1) Calculation, integral per proton:

$$\begin{aligned} &\text{Locate the end-group proton signals (ca. 5.8, 6.2 \& 6.4 ppm)} \\ &\text{integral per proton} = \frac{\text{sum of vinyl proton integrals}}{\text{\# of protons in the two vinyl end groups}} \\ &= \frac{10.00 + 9.66 + 10.17}{6} = \mathbf{4.97 \text{ per proton}} \end{aligned}$$

#### 2) Calculation, number of repeating monomer units, n:

$$\begin{aligned} &\text{Locate the OCH}_2\text{CH}_2 \text{ proton signals (ca. 3.6, 3.7 \& 4.3 ppm)} \\ &n = \frac{\text{(sum of CH}_2 \text{ proton integrals)/\# of CH}_2 \text{ protons}}{\{\text{integral per proton value}\}} \\ &= \frac{(20.79 + 151.87)/4}{4.97} = \mathbf{8.69 \text{ repeating units, n}} \end{aligned}$$

#### 3) Calculation, $M_n$ :

$$\begin{aligned} M_n &= (\text{FW end groups}) + (\text{FW repeating unit})(n) \\ &= (55.06 + 71.60) + (44.05)(8.69) = \mathbf{509} \end{aligned}$$

Therefore, the  $M_n$  of this polymer is approx. 509

For questions, product data, or new product suggestions,  
please contact the Materials Science team at [matsci@sial.com](mailto:matsci@sial.com).



Prof. Anja Mueller  
Department of Chemistry  
Central Michigan University,  
Mount Pleasant, MI

Fluorocarbon polymers, like small-molecule fluorocarbons, exhibit increased thermal stability, hydrophobicity, lipophobicity, improved chemical resistance, and decreased intermolecular attractive forces

in comparison to their hydrocarbon analogs.<sup>1</sup> These properties derive from the fundamental atomic properties of fluorine: high ionization potential, low polarizability, and high electronegativity. Due to the very high electronegativity, C–F bonds are always strongly polarized. The strength of the C–F bond is due to its highly ionic character, which accounts for the thermal stability of perfluorocarbons. The high ionization potential, combined with the low polarizability, leads to weak intermolecular interactions, which in turn leads to low surface energy and low refractive indices for perfluorocarbons. Therefore, perfluorocarbons have been used to create non-stick and non-wettable surfaces with low surface energies.

Linear fluorinated polymers, such as tetrafluoroethylene (Teflon®) exhibit high crystallinity, which increases the melting point even further. That often leads to inhibitive high processing temperatures. For applications such as mold releases or coatings, high crystallinity is often not needed or even unfavorable.


The superior chemical resistance, hydrophobicity, and low adhesive forces can be coupled with improved processibility (high solubility, low viscosity)<sup>2</sup> by making highly branched fluorocarbon polymers (**Figure 1**).<sup>3</sup> The glass transition temperature of these materials is up to 55 °C (depending on molecular weight), but they are thermally stable to 300 °C, which is sufficient for most applications. The contact angle with water for this hyperbranched fluoropolymer is just below 100° (tetrafluoroethylene: 105°), which can be increased to 120° by substituting one-third of the remaining *p*-fluorines of the structure with longer fluoroalkyl chains.

This material has improved lubricating properties and has been used as an imprinting mold release (**Figure 2**).<sup>4</sup> With a mold coated with the hyperbranched polymer, 250 nm circles and 50–60 nm lines can be imprinted without the pattern being destroyed by removing the imprinter.


For coatings applications, the hyperbranched fluorinated polymer has to be cross-linked to make it less brittle (fluorocarbon polymers have not only reduced adhesion, but also reduced cohesion). At the same time, the cross-linking molecule can be used to introduce other properties or additional functional groups.<sup>5</sup>

This family of materials thus combines the superior properties of fluorocarbon polymers with an easy synthesis and processibility, allowing for its use in a variety of applications.

**References:** (1) Chambers, R. D. *Fluorine in Organic Chemistry*; Olah, G.A., Ed.; Interscience Monograph on Organic Chemistry; John Wiley & Sons: New York, 1973. (2) Mourey, Y. H. et al. *Macromolecules* **1992**, *25*, 240. (3) Mueller, A. et al. *Macromolecules* **1998**, *31*, 776. (4) Mueller, A. *Hyperbranched Fluoropolymers: Synthesis, Characterization, Derivatization, and Applications*; Ph.D. Thesis 1998, Wooley, K.L., Adv.; Dept. of Chemistry, Washington University, St. Louis, MO. (5) Gan, D. et al. *J. Polym. Sci., Part A: Polym. Chem.* **2003**, *41*(22), 3531.

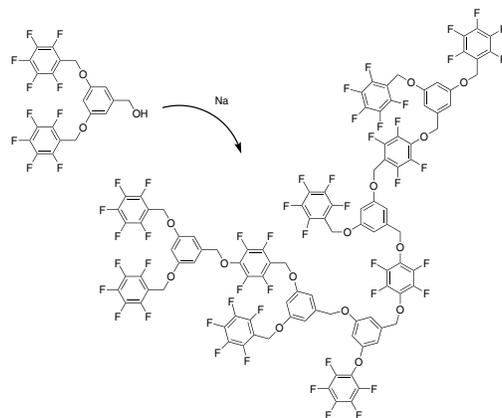

**Figure 1.** A highly cross-linked fluorinated polymer.

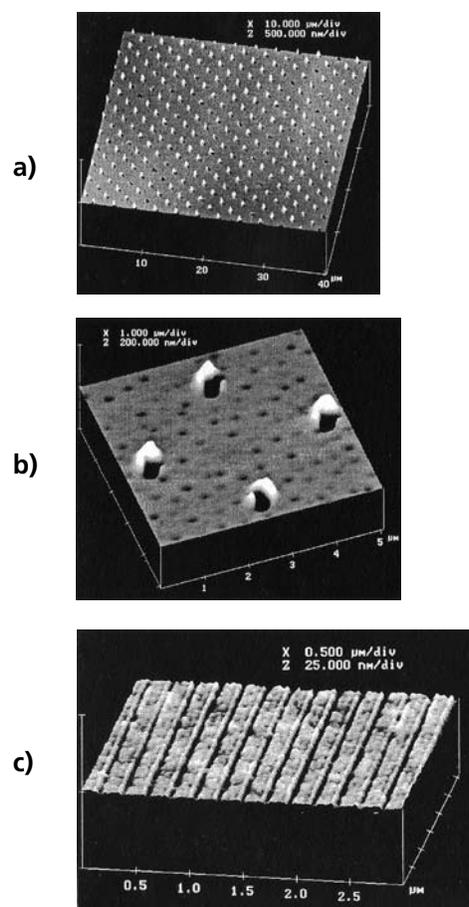

**Figure 2.** AFM images of a) 250nm punctures in ca. 100 nm thick film of fluoroalkyl-substituted HBFP; b) expelled material adjacent to punctures; c) imprint of 50–60 nm thick lines spaced 210 nm apart. Imprinting: Krchnavek, Dept. of Electrical Engineering, Washington University, St. Louis; AFM: Tomasz Kowalewski, Dept. of Chemistry, Carnegie Mellon University.

## Fluorinated Monomers

### 3,6-Difluorophthalic anhydride, 97%

$C_8H_2F_2O_3$   
MW: 184.1  
MP: 218–221 °C (lit.)



381128-100MG	100 mg
381128-500MG	500 mg


### 2-Fluorostyrene, 98%

$C_8H_7F$   
MW: 122.14  
BP: 29–30 °C (4 mm Hg) (lit.)

290505-1G	1 g
290505-5G	5 g

### 4-Fluorostyrene, 99%


$C_8H_7F$   
MW: 122.14  
BP: 67 °C (50 mm Hg) (lit.)

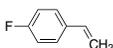

155799-1G	1 g
155799-10G	10 g

### Glycidyl 2,2,3,3-tetrafluoropropyl ether, 97%

$C_6H_8F_4O_2$   
MW: 188.12  
BP: 50 °C (4 mm Hg) (lit.)

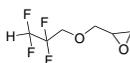

474150-5ML	5 mL
474150-25ML	25 mL

### Hexafluoroglutaric acid, 97%

$C_5H_2F_6O_4$   
MW: 240.06  
MP: 88–91 °C (lit.)  
BP: 134–138 °C (3 mm Hg) (lit.)

196908-5G	5 g
196908-25G	25 g


### 4,4'-(Hexafluoroisopropylidene)dianiline, 98%

$C_{15}H_{12}F_6N_2$   
MW: 334.26  
MP: 195–198 °C (lit.)

368148-1G	1 g
368148-5G	5 g

### 4,4'-(Hexafluoroisopropylidene)diphenol, 97%


$C_{15}H_{10}F_6O_2$   
336.23  
MP: 160–163 °C (lit.)

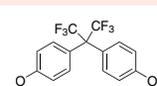

257591-5G	5 g
257591-25G	25 g
257591-100G	100 g

### 2,3,4,5,6-Pentafluorostyrene, 99%

$C_8H_3F_5$   
MW: 194.1  
BP: 62–63 °C (50 mm Hg) (lit.)

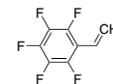

196916-5G	5 g
196916-25G	25 g

### Tetrafluorophthalic anhydride, 97%

$C_8F_4O_3$   
MW: 220.08  
MP: 94–96 °C (lit.)



339016-1G	1 g
339016-5G	5 g

### Tetrafluoroterephthalic acid, 97%

$C_8H_2F_4O_4$   
MW: 238.09  
MP: 275–277 °C (dec.) (lit.)

104418-1G	1 g
104418-5G	5 g


### 3-(Trifluoromethyl)styrene, 99%

$C_9H_7F_3$   
MW: 172.15  
BP: 64.5 °C (40 mm Hg) (lit.)

366692-1G	1 g
-----------	-----

For halogenated monomers that result in high and low refractive index polymers, see the *Advanced Polymers for Electronic/Optical Devices* technical guide. Request your free copy at [matsci@sial.com](mailto:matsci@sial.com); reference code GGE.

## Fluorinated Acrylates

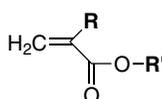

Monomer	R	R'	Quantity
2,2,3,3,3-Pentafluoropropyl acrylate, 98%	H	CF <sub>2</sub> CF <sub>2</sub> CF <sub>3</sub>	470961-5ML 470961-25ML
3,3,4,4,5,5,6,6,7,7,8,8,9,9,10,10-Heptafluorodecyl acrylate, 97%	H	CH <sub>2</sub> CH <sub>2</sub> (CF <sub>2</sub> ) <sub>7</sub> CF <sub>3</sub>	474487-5ML 474487-25ML
Heneicosfluorododecyl acrylate, 96%	H	CH <sub>2</sub> CH <sub>2</sub> (CF <sub>2</sub> ) <sub>9</sub> CF <sub>3</sub>	474355-5G
2,2,2-Trifluoroethyl methacrylate, 99%	CH <sub>3</sub>	CH <sub>2</sub> CF <sub>3</sub>	373761-5G 373761-25G
2-(Trifluoromethyl)acrylic acid, 98%	CF <sub>3</sub>	H	369144-1G 369144-5G
Zonyl® TM fluoromonomer	(CF <sub>2</sub> ) <sub>n</sub> CF <sub>3</sub> , n ~ 7–8	CH <sub>2</sub> CH <sub>3</sub>	421480-10ML 421480-50ML

For questions, product data, or new product suggestions, please contact the Materials Science team at [matsci@sial.com](mailto:matsci@sial.com).

## Etch-Resistant Block Copolymers



Prof. Padma Gopalan  
Department of Materials Science and Engineering,  
University of Wisconsin, Madison, WI

Prof. Shu Yang  
Department of Materials Science and Engineering  
University of Pennsylvania, Philadelphia, PA


### Introduction

Block copolymers offer a means of combining the desirable characteristics of different polymers in a new hybrid material. Polymers consisting of hydrophobic and hydrophilic blocks, for example, can be used to encapsulate organic molecules and deliver them into aqueous media. There has been tremendous interest in the self-assembly of block copolymers in nanoscale dimensions, especially in thin-film configuration.


Conventional lithography has its limitations when features of less than 30 nm are desired. Accessibility to a wide range of periodic structures with feature sizes less than 30 nm make block copolymers attractive as templates for nanopatterning.<sup>1-3</sup> Most of the literature approaches use selective ozonolysis or preferential staining of one block with heavy metals to increase etch selectivity between the blocks. Often, an intermediate silicon nitride (SiN) layer and selective etching of one block over another is required for successful pattern transfer. In general, the use of organic block copolymers is limited at high temperatures because of low thermal/mechanical stabilities. Thus, direct patterning of semiconductors that requires high growth temperature (>500 °C) using organic block copolymers as templates is nearly impossible.

It has been well established that incorporation of silicon (at least 10 wt %) in resist polymers provides improved oxygen-RIE (reactive ion etching) etch resistance. When exposed to oxygen plasma, the silicon-containing polymers are oxidized to silicon oxide that is stable in an O<sub>2</sub> environment. The high etch resistance to oxygen plasma compared to organic polymers makes silicon-containing polymers favorable as bilayer resists to pattern high-aspect ratio structures and to create nanoporous ceramic thin films in a variety of morphologies.<sup>4-8</sup> In addition, silicon oxide has high thermal and mechanical stability at a temperature greater than 500 °C, making it a long-time dielectric in microchip fabrication. Thus, the possibility of combining acid labile groups and silicon-containing groups in block copolymers offers a new route to directly pattern nanostructured semiconductors.

As the synthesis of silicon-containing block copolymers is quite challenging using traditional living anionic polymerization, post functionalization of polymers is often used to incorporate silicon. Recent advances in controlled living free-radical polymerization (LFRP),<sup>9-11</sup> including nitroxide-mediated radical polymerization (NMRP), atom transfer radical polymerization (ATRP), and reversible addition fragmentation chain transfer (RAFT), make it

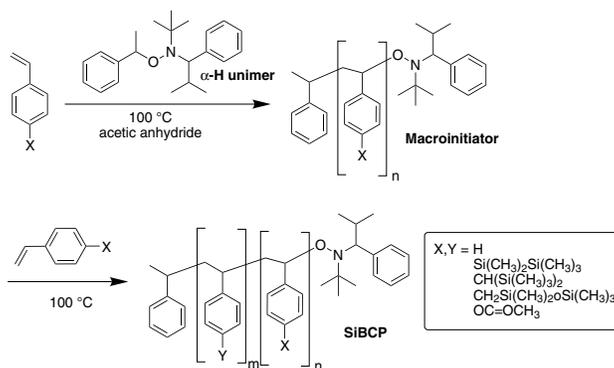

**Scheme 1.** Synthesis of macroinitiators and SiBCPs by LFRP at 100 °C.

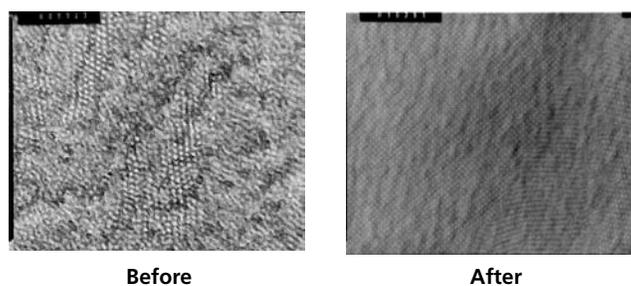

**Figure 1.** Transmission electron micrograph (TEM) images of PACOS-PSi<sub>2</sub>St (21/79 v/v) before and after O<sub>2</sub> plasma for 10 minutes, showing intact cylindrical morphology.

possible to design and synthesize a variety of block copolymers with novel functionalities. The LFRP procedures in general are easier to carry out as they are tolerant to a variety of functionalities and do not require stringent purification of the starting materials, unlike living anionic or cationic polymerization. We had recently applied NMRP towards, (i) the synthesis of narrow dispersed silicon-containing homopolymers from three kinds of silicon-containing styrenic monomers, including 4-(pentamethyldisilyl)-styrene (Si<sub>2</sub>St), 4-(bis(trimethylsilyl)-methyl)styrene (Si<sub>2</sub>-CSt), and 4-(pentamethyldisiloxy)methyl styrene (OSi<sub>2</sub>-St) (**Scheme 1**), each containing two silicon atoms to enhance the etch selectivity, and (ii) the synthesis of block copolymers from silicon-containing styrenic monomers with styrene and acid labile acetoxystyrene by sequential monomer addition using an nitroxide unimer initiator. By optimizing conditions such as solvent polarity, temperature of polymerization, and the monomer addition sequence, well-defined narrow dispersed silicon-containing block copolymers were synthesized from the above monomers. Both TEM (transmission electron microscopy) and SAXS (small angle X-ray scattering) data showed that these polymers formed cylindrical, lamellae, or disordered structures depending on the volume ratio between the blocks and their molecular weights. When the silicon-containing block was the major phase and silicon content was greater than 12 wt %, block copolymer morphology and its domain size were well maintained under exposure to oxygen plasma<sup>12</sup> (**Figure 1**).

Synthetic access to novel silicon-containing block copolymers via LFRP enables potential applications such as (1) growth of nanostructured semiconductor crystals at high temperatures, (2) formation of nanoporous ceramic films, or (3) creation of hierarchical hybrid nanostructures by combining photolithography and self-assembly of photosensitive silicon-containing block copolymers.


## Silicone-Containing Monomers

### Allyltriethoxysilane, 97%

$C_9H_{20}O_3Si$

MW: 204.34

BP: 78 °C (21 mm Hg) (lit.)

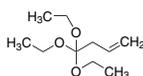

A36301-5G	5 g
A36301-25G	25 g

### Diphenylsilanediol, 95%

$C_{12}H_{12}O_2Si$

MW: 216.31

D213705-25G	25 g
D213705-100G	100 g

### Octadecyltrichlorosilane, 90+%

$C_{18}H_{37}Cl_3Si$

MW: 387.93

BP: 223 °C (10 mm Hg) (lit.)

104817-25G	25 g
104817-100G	100 g
104817-500G	500 g

### Poly(dimethylsiloxane), vinyl terminated, viscosity 1,000 centistokes

BP: >93 °C (lit.)

433012-100ML	100 mL
433012-500ML	500 mL

### Trichlorovinylsilane, 97%

$C_2H_3Cl_3Si$

MW: 161.49

MP: -95 °C (lit.)

BP: 90 °C (lit.)



104876-5G	5 g
104876-100G	100 g
104876-500G	500 g

## Styrene Monomers



Monomer	R	Quantity
Styrene, reagentplus, 99+%	H	240869-5ML 240869-100ML
4-Bromostyrene, 98%	Br	124141-1G 124141-10G 124141-25G


**References:** (1) Park, M. et al. *Science* **1997**, 276, 1401. (2) Black, C. T. et al. *Appl. Phys. Lett.* **2001**, 79, 409. (3) Kim, H. C. et al. *Adv. Mater.* **2001**, 13, 795. (4) Gabor, A. H. et al. *Chem. Mater.* **1994**, 6, 927. (5) Zharov, I. et al. *Chem. Mater.* **2002**, 14, 656. (6) Bowden, M. et al. *J. Photopolym. Sci. Technol.* **2003**, 16, 629. (7) Avgeropoulos, A. et al. *Chem. Mater.* **1998**, 10, 2109. (8) Chan, V. Z.-H. et al. *L. Science* **1999**, 286, 1716. (9) Matyjaszewski, K. *Advances in Controlled/Living Radical Polymerization*; American Chemical Society: Washington, DC, 2003; Vol. 854, p 2. (10) Hawker, C. J. *J. Am. Chem. Soc.* **1994**, 116, 1185. (11) Benoit, D. et al. *J. Am. Chem. Soc.* **1999**, 121, 3904. (12) Fukukawa, K. et al. *Macromolecules* **2005**, 38, 263.

### Trimethylsilyl methacrylate, 98%

$C_7H_{14}O_2Si$

MW: 158.27

BP: 51–51.5 °C (20 mm Hg) (lit.)

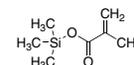

347493-25G	25 g
347493-100G	100 g

### Vinyltrimethoxysilane, 98%

$C_5H_{12}O_3Si$

MW: 148.23

BP: 123 °C (lit.)




235768-5ML	5 mL
235768-100ML	100 mL
235768-500ML	500 mL

### 1,4-bis(dimethyl(2-(5-norbornen-2-yl)ethyl)silyl)benzene, mixture of endo and exo

$C_{28}H_{42}Si_2$

MW: 434.8

BP: 314 °C (lit.)

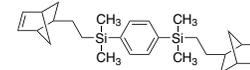


523607-5ML	5 mL
------------	------

### Hexamethyl-1,5-bis(2-(5-norbornen-2-yl)ethyl)trisiloxane, mixture of endo and exo

$C_{24}H_{44}O_2Si_3$

MW: 448.86

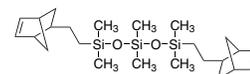

523593-1ML	1 mL
523593-5ML	5 mL

For an up-to-date list of silsesquioxane monomers, visit our nanomaterials Web site at [sigma-aldrich.com/nano](http://sigma-aldrich.com/nano).

For styrenic and other vinyl monomers, use our advanced sub-structure search: [sigma-aldrich.com](http://sigma-aldrich.com).

For questions, product data, or new product suggestions, please contact the Materials Science team at [matsci@sial.com](mailto:matsci@sial.com).

## Bioactive Hydrogels



**Anthony Guiseppi-Elie**  
ABTECH Scientific, Inc.  
Richmond, VA


The first report of the synthesis of composite materials comprising conducting polymers and hydrogels was in 1994 by Wallace et al.<sup>1</sup> Their objective was to enhance the porosity and ion-transport properties of hydrogels for controlled drug

delivery through electrochemically stimulated release of analytes. Since then, the electrochemical and oxidative polymerization of pyrrole, aniline, and thiophene and their derivatives within hydrogel hosts, such as polyacrylamide, poly(acrylic acid), chitosan, and poly(HEMA) have been routinely accomplished for biosensor applications and to achieve voltage-stimulated or controlled release. There have also been several studies conducted on the fabrication of polymer blends of hydrogels, such as poly(methyl methacrylate) (PMMA), poly(vinyl methyl ether) (PVME), poly(4-vinylpyridine),<sup>2</sup> and poly(2-hydroxyethyl methacrylate) (p(HEMA)), primarily for the construction of artificial muscles.<sup>3,4</sup> With hydrogels, high degrees of hydration (ca. 90 %) could be reversibly achieved along with biocompatibility, good refractive index matching with water, and relative ease of molecular engineering. In general, the conducting polymer component of these composites retains their electroactive properties.


Brahim et al.<sup>5</sup> have fabricated bioactive polypyrrole-p(HEMA) composites to function as sensing membranes for clinically important amperometric biosensors. A monomer cocktail containing, among other components, the relevant methacrylate monomers, pyrrole or aniline, and photoinitiator was spin-cast onto microfabricated electrodes and first irradiated by UV to effect polymerization of the hydrogel components. This was immediately followed by potentiostatic electropolymerization of the pyrrole/aniline monomer in a phosphate-buffered potassium chloride solution saturated with further monomer. Amperometric enzyme biosensors for the detection of glucose, cholesterol, and galactose were demonstrated, each possessing extensive linear dynamic response ranges, high sensitivities, and prolonged storage stabilities.<sup>6</sup>

Of particular recent interest is the development of **bioactive** (containing biologically active moieties such as bioactive peptides, growth factors, enzymes and the like) and **biosmart** (responsive to biologically derived external stimuli) electroconductive hydrogels for implant biocompatibility. A novel polymer composite material consisting of a water-dispersed complex of polypyrrole doped with polystyrenesulfonate and embedded in polyacrylamide hydrogel was prepared and evaluated as a matrix for enzyme immobilization.<sup>7</sup> The enzyme glucose oxidase was physically entrapped in the polymer by inclusion in the aqueous phase during emulsion polymerization. The resulting bioactive microparticles (3.5–7.0 μm diameter) were cast onto platinum electrodes and the polymer-modified electrodes used as amperometric glucose biosensors. This configuration displayed rapid response times and efficient screening of interferents.

For implantable biosensing applications, the synthesis of hydrogel composite polymers consisting of cross-linked p(HEMA) with incorporated polypyrrole and/or polyaniline chains are rendered “bioactive” by the covalent immobilization of oxidoreductase enzymes. Enzymes were first “monomerized” by hetero-bifunctional coupling of the amines of the lysine residues of

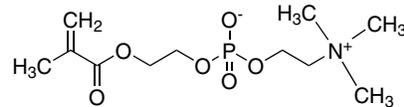

**Figure 1.** Structure of 2-methacryloyloxyethyl phosphorylcholine

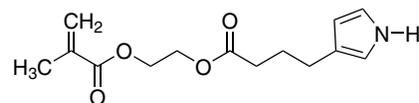

**Figure 2.** Structure of 2-methacryloyloxyethyl pyrrolylbutyrate

the enzyme (typically 1:2) with acryloyl (polyethylene glycol)<sub>110</sub> *N*-hydroxy succinamide ester (Acryl-PEG-NHS). This allowed the covalent immobilization of the tethered oxidoreductase enzyme within the hydrogel milieu. To provide for stabilization of the immobilized enzymes, poly(ethylene glycol)<sub>200</sub> monomethacrylate (PEGMA) was also included in the monomer cocktail at 0.5 mol %. Together these components allowed photolithographically defined, spin-cast membranes formed on microlithographically defined electrodes to recognize and amperometrically respond to the enzyme's substrate and achieve approximately one year of retained enzyme activity (ca. 80%). For implant biocompatibility, the synthesis of the monomer 2-methacryloyloxyethyl phosphorylcholine (MPC) (**Figure 1**) was accomplished by the coupling of HEMA with 2-chloro-2-oxo-1,3,2-dioxaphospholane (COP).<sup>8</sup> When incorporated into the hydrogels at the level of ca. 5–10 mol %, this monomer conferred nonthrombogenicity, reduced protein adsorption, and supported cell viability. Mimicking the zwitterionic head group of the outer leaflet of cell membranes, phosphatidyl choline, the phosphorylcholine moiety confers the molecular equivalent of “stelt” to the polymer when the biosensor is implanted. When cultured with muscle fibroblasts and endothelial cells, these highly porous hydrogel composites support the migration and mobility of cells within its 3-D network<sup>9</sup>; a property that will render these materials appropriate for the fabrication of nerve electrodes and cellular interfaces as the polymer mimics the biological structures that enable cells to grow.


One challenge faced by the use of electroactive polymers as components of hydrogels for mammalian implantation is the potential toxicity. Early work has established polypyrrole as bio-benign. However, to address this issue, bi-functional monomers of pyrrole such as 2-methacryloyloxyethyl pyrrolylbutyrate (MPB) (**Figure 2**) and aniline that may be UV-polymerized, and hence, covalently coupled into the hydrogel network and also oxidatively polymerized with imbibed free pyrrole or aniline monomer to form the electroactive polymer component were developed. In this way, an interpenetrating network of the electroactive polymer is formed within the preformed hydrogels network that serves as the reactor. Studies are ongoing to evaluate the potential cytotoxicity and biocompatibility of these polymers.

**References:** (1) Small, C. J. et al. *Polymer Gels and Networks* **1997**, 5, 251. (2) Asberg, P.; Ingana, O. *Biosens. Bioelectron.* **2003**, 19, 199. (3) Pich, A. et al. *Polymer* **2002**, 43, 5723. (4) Douglass, P. M, et al. *Soc. Automotive Eng.* **2000**, 1. (5) Brahim, S. et al. *Biosens. Bioelectron.* **2002**, 17, 53. (6) Brahim, S. et al. *Electroanalysis*. **2002**, 14(9), 627. (7) Rubio Retama, J. et al. *Biosens. Bioelectron.* **2004**, 20, 1111. (8) Brahim, S. et al. *Microchimica Acta* **2003**, 143, 123. (9) Abraham, S et al. *Biomaterials* **2005**, 26(23), 4767. (10) Abraham, S.; Guiseppi-Elie, A. *Biomaterials* **2006**, submitted for publication.

## Monomers for Bioactive Polymers

### Benzyl 2-ethyl acrylate, 99% NEW


C<sub>12</sub>H<sub>14</sub>O<sub>2</sub>  
MW: 190.24

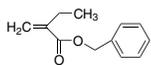

589136-250MG 250 mg

### Benzyl 2-propylacrylate, 99% NEW

C<sub>13</sub>H<sub>16</sub>O<sub>2</sub>  
MW: 204.26

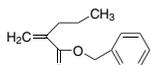

590126-250MG 250 mg


### 1,6-Bis(p-acetoxycarbonylphenoxy) hexane, 97% NEW

C<sub>24</sub>H<sub>26</sub>O<sub>8</sub>  
MW: 442.46

657174-1G 1 g

### 1,3-Bis(4-carboxyphenoxy)propane, 97% NEW

C<sub>17</sub>H<sub>16</sub>O<sub>6</sub>  
MW: 316.31  
MP: 310 °C(lit.)

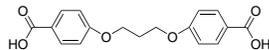

655538-5G 5 g

### 1,6-Bis(p-carboxyphenoxy)hexane, 90% NEW

C<sub>20</sub>H<sub>22</sub>O<sub>6</sub>  
MW: 358.39

655546-5G 5 g

### 2-Chloro-1,3,2-dioxaphospholane-2-oxide NEW

C<sub>2</sub>H<sub>4</sub>ClO<sub>3</sub>P  
MW: 142.48  
MP: 12–14 °C (lit.)  
BP: 89–91 at 0.8 mm Hg (lit.)



377953-1G 1 g

377953-5G 5 g

### 2-Ethylacrylic acid, 98% NEW

C<sub>5</sub>H<sub>8</sub>O<sub>2</sub>  
MW: 100.12  
BP: 176 °C (lit.)



589128-250MG 250 mg


### Ethyl 2-propyl acrylate, 99%

C<sub>8</sub>H<sub>14</sub>O<sub>2</sub>  
MW: 142.2  
BP: 141 °C (lit.)

590118-250MG 250 mg

### Glycosyloxyethyl methacrylate, 5% (w/v) solution in ethanol NEW


C<sub>12</sub>H<sub>20</sub>O<sub>8</sub>  
MW: 292.28

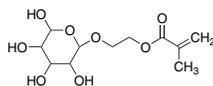

659576-25ML 25 mL

### 2-Propylacrylic acid, 99% NEW

C<sub>6</sub>H<sub>10</sub>O<sub>2</sub>  
114.14  
165–188 °C (lit.)

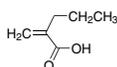

591009-100MG 100 mg

591009-1G 1 g

### Pyrrole, 98%

C<sub>4</sub>H<sub>5</sub>N  
MW: 67.09  
MP: –23 °C (lit.)  
BP: 131 °C (lit.)



131709-25ML 25 mL

131709-100ML 100 mL

131709-500ML 500 mL

### 2-Vinylpyridine, 97%

C<sub>7</sub>H<sub>7</sub>N  
MW: 105.14  
BP: 79–82 °C (29 mm Hg) (lit.)



132292-5ML 5 mL

132292-100ML 100 mL

132292-500ML 500 mL

### 4-Vinylpyridine, 95%

C<sub>7</sub>H<sub>7</sub>N  
MW: 105.14  
BP: 62–65 °C (15 mm Hg) (lit.)



V3204-5ML 5 mL

V3204-100ML 100 mL

V3204-500ML 500 mL

### 1-Vinyl-2-pyrrolidinone, 99+% NEW

C<sub>6</sub>H<sub>9</sub>NO  
MW: 111.14  
BP: 92–95 °C (11 mm Hg) (lit.)



V3409-5G 5 g

V3409-250G 250 g

V3409-1KG 1 kg

V3409-18KG 18 kg

Three grades of HEMA (2-hydroxyethyl methacrylate) are available: 97% (128635), 98% (525464), 99+% (477028).

For a comprehensive list of hydrogel hosts such as polyacrylamide, poly(acrylic acid), chitosan and poly(HEMA), as well as functionalized PEGs (linear, 4-arm and 6-arm), visit [sigma-aldrich.com/biocomp](http://sigma-aldrich.com/biocomp).

## Cross-linkers



DVB, 85% (535583)

Cross-linking is the formation of chemical links between molecular chains to form a three-dimensional network of connected molecules. The strategy of covalent cross-linking is key to the formation of hydrogels. It is also used in several other technologies of commercial and scientific interest to control and enhance the properties of the resulting polymer system or interface, such as thermosets and coatings.

For a complete list of cross-linkers, visit us at [sigma-aldrich.com/biocomp](http://sigma-aldrich.com/biocomp) and scroll down to cross-linkers.

For questions, product data, or new product suggestions, please contact the Materials Science team at [matsci@sial.com](mailto:matsci@sial.com).

## Anionic Polymerization




Prof. Roderic P. Quirk  
and  
Ms. Manuela Ocampo  
Maurice Morton  
Institute of  
Polymer Science  
The University of  
Akron, Akron, OH

Living anionic polymerization, especially using alkyl lithium initiators, has been demonstrated to be a convenient and useful method to make well-defined polymers with low degrees of compositional heterogeneity and with control of the major structural variables that affect polymer properties.<sup>1,2</sup> Living polymerizations are chain-reaction polymerizations that proceed in the absence of the kinetic steps of chain termination and chain transfer. For a living polymerization, one initiator molecule generates one polymer molecule; thus, it is possible to calculate and control the number average molecular weight ( $M_n$ ) of the final polymer via the stoichiometry of the reaction using the following relationship.

$$M_n = \text{g of monomer consumed / moles of initiator}$$

Given a comparable or faster rate of initiation relative to propagation, it is possible to obtain narrow molecular weight distribution polymers, i.e.,  $M_w/M_n \leq 1.1$ .<sup>3</sup> Due to the absence of termination and transfer steps, the product after complete monomer consumption is a reactive, polymeric organolithium compound. The living nature of alkyl lithium-initiated anionic polymerizations using suitable monomers provides versatile methods for the preparation of well-defined block copolymers by sequential addition of monomers,<sup>4</sup> chain-end functionalized polymers by reaction of the living chain ends with appropriate monomers and/or electrophilic terminating agents<sup>5,6</sup> and branched polymers by linking reactions with multi-functional linking agents.<sup>7</sup>

The monomers that can be polymerized anionically are classified into two categories: (a) unsaturated monomers with one or more double bonds, such as vinyl (e.g., styrenes, vinylpyridines, alkyl methacrylates), dienes (e.g., isoprene, 1,3-butadiene) and carbonyl-type monomers (e.g., formaldehyde); and (b) heterocyclic monomers (e.g., epoxides, thiranes, lactones, lactams, and siloxanes). In the case of vinyl monomers, the presence of electron-withdrawing substituents (e.g., X, Y) in the double bond is generally required to stabilize the negative charge that develops in the transition state as shown below.

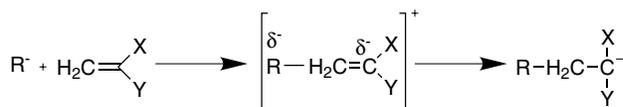

### Organolithium Initiators

Of all alkali metals, lithium is unique in that it exhibits the highest electronegativity, the smallest covalent and ionic bond radii, along with low-lying, unoccupied p-orbitals available for bonding.<sup>8,9</sup> Organolithium compounds are unique among organoalkali compounds in exhibiting properties characteristic of both covalent and ionic compounds. Thus, they are aggregated in solution, in the solid state and in the gas phase, and they are generally soluble in hydrocarbon solution. In general, the initiation of anionic polymerization of styrene and diene monomers is effected with alkyl lithium compounds such as *sec*-butyllithium (**195596**) and

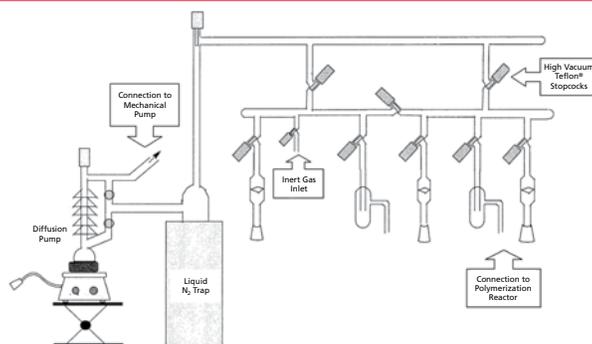


Figure 1. The typical construction of a glass high vacuum line for anionic polymerization.<sup>15-17</sup>

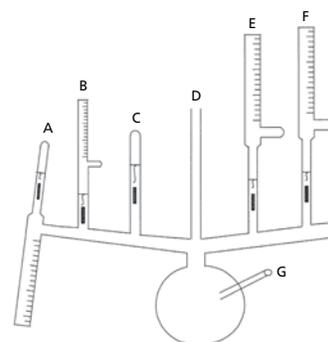


Figure 2. General set-up for a glass polymerization reactor.

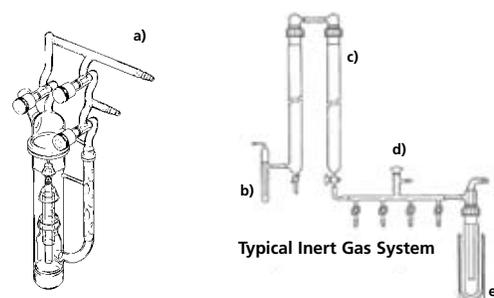


Figure 3. Aldrich products a) Z544787 b) Z173053 c) Z174254 d) Z174432 e) Z220418

*n*-butyllithium (**186171**, **230707**, **230715**, **302104**, and **302120**) in hydrocarbon solution. Under these conditions, the unique ability of organolithium compounds to effect 1,4-enchainment of 1,3-dienes is achieved.<sup>1,10</sup> The concentrations of solutions of active organolithium compounds can be determined using the Gilman double titration method.<sup>11</sup>

### Experimental Methods

Due to the reactivity of organolithium compounds and other carbanionic species toward impurities such as oxygen, moisture or carbon dioxide,<sup>12</sup> it is necessary to exclude these contaminants from the reaction environment by the use of an inert gas atmosphere<sup>13,14</sup> or high vacuum techniques.<sup>15-17</sup>

**High Vacuum Techniques** The use of high vacuum techniques provides the most effective experimental method to exclude impurities from the reaction system.<sup>15-17</sup> In order to attain high vacuum, the combination of a mechanical pump and an oil diffusion pump (**Z220418**) is used in conjunction with a two-stage glass manifold as shown in **Figure 1**.

In order to achieve the desired levels of purity for controlled anionic polymerization, all monomers, reactants, and solvents should be purified, dried, and degassed, preferably on the vacuum line. Solvents are distilled directly into the requisite glass reactors (**Figure 2**) via **D** followed by flame sealing from the vacuum line. Ampules **B**, **E**, and **F** contain monomers or functionalizing agents. Ampule **C** contains a terminating agent such as degassed methanol. Ampule **A** is equipped with a degassed methanol tube, and it is used to remove a base sample of the living polymer.

**Schlenk Line and Glove Box Techniques** are often suitable for carrying out many living anionic polymerization procedures. Alkylolithium-initiated polymerizations are somewhat forgiving in the sense that one can add a calculated excess of initiator to clean the reactor/solvent/monomer system of reactive impurities. See **Figure 3** for representative Schlenk Line Glassware. For representative and the Sigma-Aldrich glass center, visit [sigma-aldrich.com/glass](http://sigma-aldrich.com/glass).

## Safety Considerations

Vacuum traps should be vented while warming because of the possibility of trapped, liquefied gases. Hydrocarbon solutions of alkylolithium compounds are air- and moisture-sensitive; they should be either handled under an inert atmosphere or using syringes and recommended procedures for handling air-sensitive compounds.<sup>13</sup> Carbon dioxide extinguishers should not be used because Li

compounds and many other organometallic compounds react with carbon dioxide exothermically. An all-purpose fire extinguisher, or one designed specifically for combustible metals, should be available when working with these organometallic compounds and alkali metals.<sup>18</sup>

**References:** (1) Hsieh, H. L.; Quirk, R. P. *Anionic Polymerization: Principles and Practical Applications*; Dekker: New York, 1996. (2) Quirk, R.P. Anionic Polymerization. In *Encyclopedia of Polymer Science and Technology* Kroschwitz, J. I., Ed.; 3rd ed.; Wiley-Interscience: New York, 2003; Vol. 5, p 111. (3) Fetters, L. J. Monodisperse Polymers. In *Encyclopedia of Polymer Science and Engineering* Kroschwitz, J. I., Ed.; 2nd ed.; Wiley-Interscience: New York, 1985; Vol. 2, p 478. (4) Hadjichristidis, N.; Pispas, S.; Floudas, G. A. *Block Copolymers: Synthetic Strategies, Physical Properties, and Applications*; Wiley-Interscience: New York, 2003. (5) Quirk, R. P. Anionic Synthesis of Polymers with Functional Groups. In *Comprehensive Polymer Science, First Supplement*; Aggarwal, S. L., Russo, S., Eds.; Pergamon Press: Oxford, 1992; p 83. (6) Hirao, A.; Hayashi, M. *Acta Polym.* **1999**, *50*, 219. (7) Hadjichristidis, N. et al. *Chem. Rev.* **2001**, *101*, 3747. (8) Wardell, J. L. Alkali Metals. In *Comprehensive Organometallic Chemistry: The Synthesis, Reactions and Structures of Organometallic Compounds*; Wilkinson, G., Gordon, F., Stone, A. Abel, E. W., Eds.; Pergamon Press: Oxford; 1982; Vol. 1, p 43. (9) Sanderson, R. T. *Chemical Periodicity*; Reinhold: New York; 1960. (10) Bywater, S. In *Comprehensive Polymer Science*; Eastmond, G. C., Ledwith, A., Russo, S., Sigwalt, P., Eds.; Chain Polymerization I; Pergamon Press: Elmsford, New York, 1989, Vol. 3, p. 433. (11) Gilman, H.; Cartledge, F. K. *J. Organomet. Chem.* **1964**, *2*, 447. (12) Wakefield, B. J. *The Chemistry of Organolithium Compounds*; Pergamon Press: New York; 1974. (13) Shriver, D. F.; Drezdson, M. A. *The manipulation of Air-Sensitive Compounds*; Wiley: New York, 1986. (Cat. No. Z558486) (14) Ndoni, S.; et al. *Rev. Sci. Instrum.* **1995**, *66*, 1090. (15) Morton, M.; Fetters, L. J. *Rubber Chem. Technol.* **1975**, *48*, 359. (16) Hadjichristidis, N. et al. *J. Polym. Sci., Part A: Polym. Chem.* **2000**, *38*, 3211. (17) Uhrig, D; Mays, J. W. *J. Polym. Sci., Part A: Polym. Chem.* **2005**, *43*, 6179. (18) Wietelmann, U.; Bauer, R. J. In *Ullmann's Encyclopedia. Industrial Inorganic Chemicals and Products*; Wiley-VCH Verlag: Weinheim, Germany; 1998; Vol. 4, p. 2899.

For the manipulation of air-sensitive compounds, see also Aldrich Technical Bulletins **AL-134**, **AL-136**, **AL-164**, and **AL-166**.

## Monomers for Anionic Polymerization



	X	Y
<b>Vinylidene chloride, 99%</b> 163023-100G 163023-500G	H	Cl
<b>Acrylic acid, 99%</b> 147230-5G 147230-100G 147230-500G	H	
<b>Methyl acrylate, 99%</b> M27301-5ML M27301-250ML	H	
<b>Ethyl acrylate, 99%</b> E9706-5ML E9706-100ML	H	
<b>Butyl acrylate, 99+%</b> 234923-5ML 234923-100ML	H	
<b>2-Hydroxyethyl acrylate, 96%</b> 292818-5ML 292818-250ML	H	
<b>2-Hydroxy-3-phenoxypropyl acrylate</b> 407364-100ML 407364-500ML	H	
<b>N,N-Dimethylacrylamide, 99%</b> 274135-5ML 274135-100ML 274135-500ML	H	
<b>N-Isopropylacrylamide, 97%</b> 415324-10G 415324-50G	H	

	X	Y
<b>Methacrylic acid, 99%</b> 155721-5G 155721-100G 155721-500G	CH <sub>3</sub>	
<b>Methyl methacrylate, 99%</b> M55909-25ML M55909-500ML	CH <sub>3</sub>	
<b>Ethyl methacrylate, 99%</b> 234893-5ML 234893-100ML 234893-500ML	CH <sub>3</sub>	
<b>2-Hydroxyethyl methacrylate, 99+%</b> 128635-5G	CH <sub>3</sub>	
<b>Glycidyl methacrylate, 97%</b> 151238-5G 151238-100G 151238-500G	CH <sub>3</sub>	
<b>2-Isocyanatoethyl methacrylate, 98%</b> 477060-5ML 477060-50ML	CH <sub>3</sub>	
<b>Methyl 2-(bromomethyl)acrylate, 97%</b> 302546-250MG 302546-1G 302546-5G	CH <sub>2</sub> Br	

For a comprehensive list of acrylate-methacrylate monomers visit [sigma-aldrich.com/polymer](http://sigma-aldrich.com/polymer).

**Argentina**

SIGMA-ALDRICH DE ARGENTINA, S.A.  
Tel: 54 11 4556 1472  
Fax: 54 11 4552 1698

**Australia**

SIGMA-ALDRICH PTY., LIMITED  
Free Tel: 1800 800 097  
Free Fax: 1800 800 096  
Tel: 612 9841 0555  
Fax: 612 9841 0500

**Austria**

SIGMA-ALDRICH HANDELS GmbH  
Tel: 43 1 605 81 10  
Fax: 43 1 605 81 20

**Belgium**

SIGMA-ALDRICH NV/SA.  
Free Tel: 0800-14747  
Free Fax: 0800-14745  
Tel: 03 899 13 01  
Fax: 03 899 13 11

**Brazil**

SIGMA-ALDRICH BRASIL LTDA.  
Tel: 55 11 3732-3100  
Fax: 55 11 3733-5151

**Canada**

SIGMA-ALDRICH CANADA LTD.  
Free Tel: 800-565-1400  
Free Fax: 800-265-3858  
Tel: 905-829-9500  
Fax: 905-829-9292

**China**

SIGMA-ALDRICH CHINA INC.  
Tel: 86-21-6386 2766  
Fax: 86-21-6386 3966

**Czech Republic**

SIGMA-ALDRICH S.R.O.  
Tel: +420 246 003 200  
Fax: +420 246 003 291

**Denmark**

SIGMA-ALDRICH DENMARK A/S  
Tel: 43 56 59 10  
Fax: 43 56 59 05

**Finland**

SIGMA-ALDRICH FINLAND  
Tel: (09) 350 9250  
Fax: (09) 350 92555

**France**

SIGMA-ALDRICH CHIMIE S.à.r.l.  
Tel appel gratuit: 0800 211 408  
Fax appel gratuit: 0800 031 052

**Germany**

SIGMA-ALDRICH CHEMIE GmbH  
Free Tel: 0800-51 55 000  
Free Fax: 0800-649 00 00

**Greece**

SIGMA-ALDRICH (O.M.) LTD  
Tel: 30 210 9948010  
Fax: 30 210 9943831

**Hungary**

SIGMA-ALDRICH Kft  
Tel: 06-1-235-9054  
Fax: 06-1-269-6470  
Ingyenes zöld telefon: 06-80-355-355  
Ingyenes zöld fax: 06-80-344-344

**India**

SIGMA-ALDRICH CHEMICALS  
PRIVATE LIMITED  
Telephone  
Bangalore: 91-80-4112-7272  
New Delhi: 91-11-4165-4255  
Mumbai: 91-22-2570-2364  
Hyderabad: 91-40-5584-5488  
Fax  
Bangalore: 91-80-4112-7473  
New Delhi: 91-11-4165-4266  
Mumbai: 91-22-2579-7589  
Hyderabad: 91-40-5584-5466

**Ireland**

SIGMA-ALDRICH IRELAND LTD.  
Free Tel: 1800 200 888  
Free Fax: 1800 600 222  
Tel: 353 1 4041900  
Fax: 353 1 4041910

**Israel**

SIGMA-ALDRICH ISRAEL LTD.  
Free Tel: 1-800-70-2222  
Tel: 08-948-4100  
Fax: 08-948-4200

**Italy**

SIGMA-ALDRICH S.r.l.  
Telefono: 02 33417310  
Fax: 02 38010737  
Numero Verde: 800-827018

**Japan**

SIGMA-ALDRICH JAPAN K.K.  
Tokyo Tel: 03 5796 7300  
Tokyo Fax: 03 5796 7315

**Korea**

SIGMA-ALDRICH KOREA  
Tel: 031-329-9000  
Fax: 031-329-9090

**Malaysia**

SIGMA-ALDRICH (M) SDN. BHD  
Tel: 603-56353321  
Fax: 603-56354116

**Mexico**

SIGMA-ALDRICH QUÍMICA, S.A. de C.V.  
Free Tel: 01-800-007-5300  
Free Fax: 01-800-712-9920

**The Netherlands**

SIGMA-ALDRICH CHEMIE BV  
Tel Gratis: 0800-0229088  
Fax Gratis: 0800-0229089  
Tel: 078-6205411  
Fax: 078-6205421

**New Zealand**

SIGMA-ALDRICH PTY., LIMITED  
Free Tel: 0800 936 666  
Free Fax: 0800 937 777  
Tel: 61 2 9841 0500  
Fax: 61 2 9841 0500

**Norway**

SIGMA-ALDRICH NORWAY AS  
Tel: 23 17 60 60  
Fax: 23 17 60 50

**Poland**

SIGMA-ALDRICH Sp. z o.o.  
Tel: 061 829 01 00  
Fax: 061 829 01 20

**Portugal**

SIGMA-ALDRICH QUÍMICA, S.A.  
Free Tel: 800 202180  
Free Fax: 800 202178  
Tel: 21 9242555  
Fax: 21 9242610

**Russia**

SIGMA-ALDRICH RUS, LLC  
Tel: +7 (095) 621-5828/6037  
Fax: +7 (095) 975-4792

**Singapore**

SIGMA-ALDRICH PTE. LTD.  
Tel: 65-67791200  
Fax: 65-67791822

**South Africa**

SIGMA-ALDRICH  
SOUTH AFRICA (PTY) LTD.  
Free Tel: 0800 1100 75  
Free Fax: 0800 1100 79  
Tel: 27 11 979 1188  
Fax: 27 11 979 1119

**Spain**

SIGMA-ALDRICH QUÍMICA S.A.  
Free Tel: 900 101376  
Free Fax: 900 102028  
Tel: 91 661 99 77  
Fax: 91 661 96 42

**Sweden**

SIGMA-ALDRICH SWEDEN AB  
Tel: 020-350510  
Fax: 020-352522  
Outside Sweden Tel: +46 8 7424200  
Outside Sweden Fax: +46 8 7424243

**Switzerland**

SIGMA-ALDRICH CHEMIE GmbH  
Swiss Free Call: 0800 80 00 80  
Tel: +41 81 755 2828  
Fax: +41 81 755 2815

**United Kingdom**

SIGMA-ALDRICH COMPANY LTD.  
Free Tel: 0800 717181  
Free Fax: 0800 378785  
Tel: 01747 833000  
Fax: 01747 833313  
SAFC (UK): 01202 712305

**United States**

SIGMA-ALDRICH  
P.O. Box 14508  
St. Louis, Missouri 63178  
Toll-free: 800-325-3010  
Call Collect: 314-771-5750  
Toll-Free Fax: 800-325-5052  
Tel: 314-771-5765  
Fax: 314-771-5757

**Internet**

sigma-aldrich.com



**ALDRICH®**

P.O. Box 14508  
St. Louis, MO 63178  
USA

**Order/Customer Service** 1-800-325-3010 • Fax 1-800-325-5052  
**Technical Service** 1-800-325-5832 • sigma-aldrich.com/techservice

**Development/Bulk Manufacturing Inquiries** SAFC™ 1-800-244-1173

**World Headquarters** • 3050 Spruce St., St. Louis, MO 63103 • (314) 771-5765

*Accelerating Customers' Success through Leadership in Life Science, High Technology and Service*

*The SIGMA-ALDRICH Group* SIGMA ALDRICH Fluka Riedel-deHaën SAFC SUPELCO ISOTEC

©2006 Sigma-Aldrich Co. All rights reserved.  
SIGMA, SAFC, SAFC, SIGMA-ALDRICH, ISOTEC, ALDRICH, FLUKA, and SUPELCO are trademarks belonging to Sigma-Aldrich Co. and its affiliate Sigma-Aldrich Biotechnology LP. Riedel-de Haën\*: trademark under license from Riedel-de Haën GmbH. Sigma products are sold through Sigma-Aldrich, Inc. Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip. Teflon is a registered trademark of DuPont. Prices subject to change.



IRM  
00589-503403  
0036



**SIGMA-ALDRICH**