


YOUR **DIGITAL** QC COMPANION FOR STERILITY TESTING & OTHER REGULATED PROCESSES

M-Trace™ Software for Trusted and Efficient Data Recording

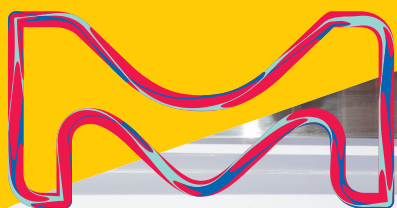

The life science business
of Merck operates as
MilliporeSigma in the
U.S. and Canada.



Millipore®

Preparation, Separation,
Filtration & Monitoring Products



STERILISATION

ECLAIRAGE

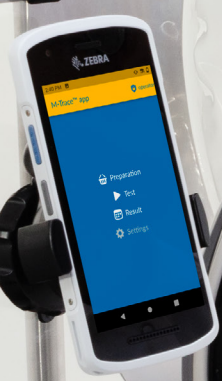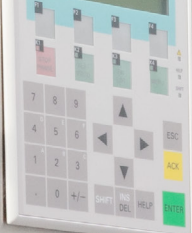



HOW DOES IT WORK?

This is how a fully traceable, reliable and efficient sterility test performed with the assistance of the M-Trace™ software could look like at your site.



Get started

Transfer and digitalize your SOPs in the M-Trace™ software, optionally connect the system with your QC-IT environment such as LIMS and finally connect the M-Trace™ software with your Steritest® Symbio pump to be all set.

Workflow preparation

Scan LIMS label of sample with the M-Trace™ mobile device to register sample information in the M-Trace™ software.

Collect and register consumables or equipment needed for testing based on item list provided by the M-Trace™ mobile device to make sure no item is forgotten.

**DESIGNED
TO ENSURE
COMPLIANCE**

**DIGITALIZED
TO BOOST
PRODUCTIVITY**

**GUIDANCE
TO AVOID
ERRORS**

**A PRE-
VALIDATED
SETUP**



Report

A comprehensive digital test report is automatically generated for convenient and reliable test record review.



Examination

Your personalized SOPs allow you to scan the test sample anywhere at your facility, so you can continuously document the results with the M-Trace™ mobile app. This app also lets you organize your sample management and reference a sample ID to a specific incubator.



Testing

Scan consumables for your work steps and be alerted if an expired or invalid consumable was selected. Get work instructions either by voice control or on the Steritest® Symbio pump and confirm proper execution real-time. In all this, the pump is controlled semi-automatically.

BOOST PRODUCTIVITY ENSURE COMPLIANCE

with the M-Trace™ solution



FULL TRACEABILITY

for your sterility testing data with the M-Trace™ solution

Ordering information



Product		Ordering No.
M-Trace™ All-in-One Computer	All-in-One Computer for running the M-Trace™ Electronic Test Record Software	MTRACEPC1
M-Trace™ Touch Mobile Computer	Touch Mobile Computer for running the M-Trace™ Mobile App	MTRACEMC1
M-Trace™ Charging Cradle	Charging Cradle for M-Trace™ Touch Mobile Computer	MTRACECR1
M-Trace™ Headset	Headset for voice control of the M-Trace™ Electronic Test Record Software	MTRACEHS1
M-Trace™ Mobile Computer Mount	Mount for the M-Trace™ Mobile Computer	MTRACEMM1
M-Trace™ Battery	Replacement battery of the M-Trace™ Touch Mobile Computer	MTRACEBT1
M-Trace™ Backup Drive	Drive to create backups out of the M-Trace™ Electronic Test Record Software	MTRACEBD1
M-Trace™ Software - Annual License Fee	Annual license to run the M-Trace™ Electronic Test Record Software	MTRASOFTW
M-Trace™ On-Site Installation Service	Service for on-site installation of the M-Trace™ hardware and software	MTRAINSTAL
M-Trace™ On-Site IQ OQ Service	Service for on-site IQ/OQ of the M-Trace™ hardware and software	MTRAVAL11
M-Trace™ IQ OQ e-Protocol	Electronic protocol for IQ/OQ of the M-Trace™ hardware and software	MTRAIQOQEP
M-Trace™ On-Site Software User Training	On-site end-user training of the M-Trace™ Electronic Test Record Software	MTRATRAIN
M-Trace™ Integration Service	Customized integration service of the M-Trace™ hardware and software into end-users IT environment	MTRAINTEG
M-Trace™ Master Test Record Creation	Service to create customized Master Test Records in the M-Trace™ Electronic Test Record Software	MTRAMTRC

SHAPE YOUR FUTURE

We help you bring robotics and automation together

Our mission is to transform microbiology quality control processes to accelerate the release of safe pharmaceutical products. We collaborate with the global scientific community to develop integrated automation and robotic solutions to improve QC testing reliability and data integrity.

Learn more:
SigmaAldrich.com/shapeyourfuture

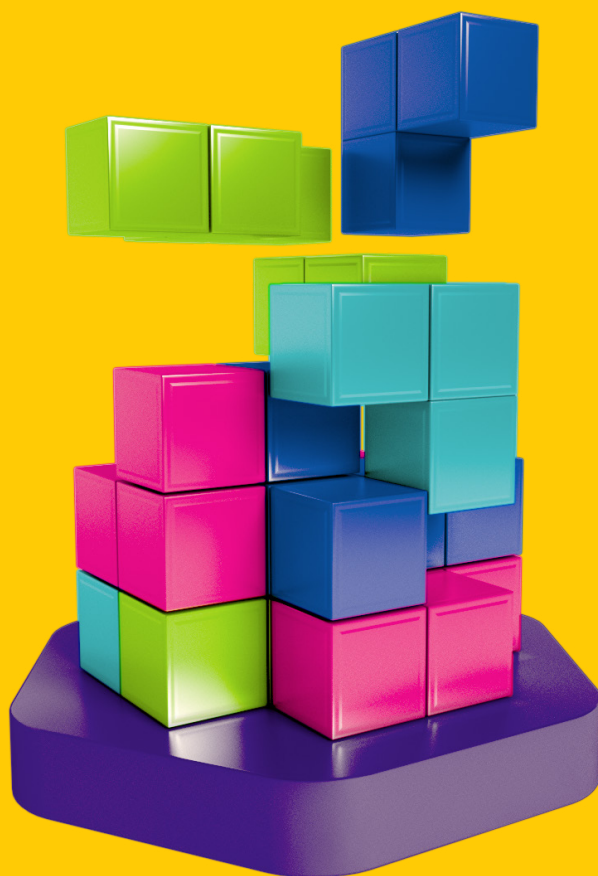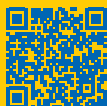

Learn more at **SigmaAldrich.com/M-Trace**

To Place an Order or Receive Technical Assistance

Order/Customer Service: **SigmaAldrich.com/order**

Technical Service: **SigmaAldrich.com/techservice**

Merck KGaA
Frankfurter Strasse 250
64293 Darmstadt, Germany

