

BioTracker™ 640 Red C2M(FM5-95) Synaptic Dye

Cat. # SCT129

pack size: 5 mg

FOR RESEARCH USE ONLY.
NOT FOR USE IN DIAGNOSTIC PROCEDURES.
NOT FOR HUMAN OR ANIMAL CONSUMPTION.

Store at 2-8°C



Data Sheet

page 1 of 2

Background

BioTracker™ synaptic dyes are a series of fluorescent cationic styryl dyes developed to follow synaptic activities at neuromuscular junctions or synapses. These dyes typically have a lipophilic tail at one end and a highly hydrophilic at the other end. In the presence of cells or tissue preparations, the dyes partition between the aqueous phase, where the dyes are virtually non-fluorescent, and the outer leaflet of the cell surface membranes, where the dyes insert the lipophilic end into the membranes and become intensely fluorescent. During endocytosis following nerve stimulation, the dyes become trapped inside the vesicles. Thus, after washing off the dyes on the cell surface, the fluorescent signal is proportional to the number of newly formed vesicles. On the other hand, during exocytosis, the dyes are released from the vesicles along with neurotransmitters, causing a decrease in fluorescent signal. As a result, the change in fluorescent intensity reflects the amount of endocytosis/exocytosis or synaptic activity.

BioTracker™ 640 Red C2M(FM5-95) Synaptic Dye is a red fluorescent probe used to follow synaptic vesicles. The dye is slightly more water soluble than BioTracker™ 640 Red C2(FM4-64) Synaptic Dye.

Storage

Store BioTracker™ 640 Red C2M(FM5-95) Synaptic Dye at 2-8°C, protected from light. Stock solutions can be prepared at 10 mM and stored at 4°C or -20°C for six months or longer.

Spectral Properties

Absorbance (Membranes): 510 nm

Emission (Membranes): 750 nm

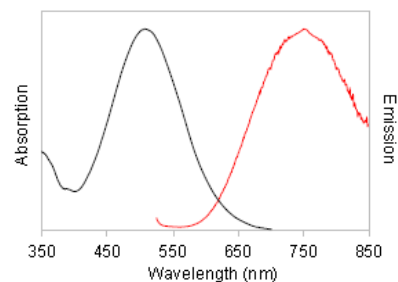


Figure 1. Absorption and emission spectra of BioTracker™ 640 Red C2M(FM5-95) Synaptic Dye in liposomes.

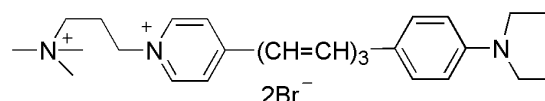


Figure 2. Chemical structure of the BioTracker™ 640 Red C2M(FM5-95) Synaptic Dye.

Please visit www.millipore.com for additional product information and references.

Submit your published journal article, and earn credit toward future purchases. Visit www.millipore.com/publicationrewards to learn more!

Assay Protocol

The following is an example of a protocol for nerve terminal staining of cultured neurons on coverslips. Nerve terminal dyes also can be used to label endocytic vesicles in non-neuronal cell types. Staining can be performed at 4°C for selective labeling of the plasma membrane; at room temperature or 37° C, endocytosis of the dye generally occurs within 10 minutes. Buffers other than Tyrode solution may be used. The addition of the sodium channel blocker tetrodotoxin (TTX) is optional, its purpose is to block action potentials and prevent synaptic vesicle release after staining. Optimal protocols for specific applications may need to be determined by the user.

1. Dilute BioTracker™ synaptic dye to a final concentration of 4 µM in 50 mM Tyrode solution. Place the coverslip with your cells in this solution for 1 minute at room temperature. Use enough solution to completely submerge the cells.
2. Transfer the coverslip to Tyrode + 0.5 µM tetrodotoxin (TTX) solution for 1 minute at room temperature.
3. Wash the coverslip several times in Tyrode + TTX at room temperature.

Note: to reduce background, 1 mM ADVASEP-7 can be added to the wash solution. Alternatively, SCAS can be used to quench background without repeated washes. Incubate the coverslip for 4 minutes at room temperature in Tyrode + TTX + 0.5 mM SCAS.

4. Mount the coverslip in Tyrode + TTX and image.

Note: 50 µM sulforhodamine 101 can be included during mounting to quench extracellular fluorescence.

Note: BioTracker™ synaptic dyes are not fixable.

BioTracker™ is a registered trademark of Merck KGaA

■ antibodies ■ Multiplex products ■ biotools ■ cell culture ■ enzymes ■ kits ■ proteins/peptides ■ siRNA/cDNA products

Please visit www.millipore.com for additional product information, test data and references

EMD Millipore Corporation, 28820 Single Oak Drive, Temecula, CA 92590, USA 1-800-437-7500

Technical Support: T: 1-800-MILLIPORE (1-800-645-5476) • F: 1-800-437-7502

FOR RESEARCH USE ONLY. Not for use in diagnostic procedures. Not for human or animal consumption. Purchase of this Product does not include any right to resell or transfer, either as a stand-alone product or as a component of another product. Any use of this Product for purposes other than research is strictly prohibited.

EMD Millipore®, the M mark, Upstate®, Chemicon®, Linco® and all other registered trademarks, unless specifically identified above in the text as belonging to a third party, are owned by Merck KGaA, Darmstadt, Germany. Copyright ©2008-2019 Merck KGaA, Darmstadt, Germany. All rights reserved.



We Buy 100% Certified
Renewable Energy