

PhotoCol®-IRG, Methacrylated Collagen Hydrogel Kit

3D CC Hydrogel

Cat. # CC320

pack size: 1 Kit

FOR RESEARCH USE ONLY.
NOT FOR USE IN DIAGNOSTIC PROCEDURES.
NOT FOR HUMAN OR ANIMAL CONSUMPTION.



Data Sheet

page 1 of 3

Background

3D cell culture, including bioprinting, allows for the creation of more physiological cell models by allowing cells to simultaneously interact with integrins on all cell surfaces, resulting in the activation of specific signaling pathways not activated in traditional 2D cell culture methods. Hydrogels are water swollen polymers that allow for the culture of cells in 3-dimensions and can have profound effects on cellular development, differentiation, migration, and function. New areas of tissue engineering such as 3D bioprinting, have utilized UV photocrosslinked methacrylated hydrogel biomaterials (PEGMA, GelMA, HAMA and ColMA etc.) to encapsulate cells to make printable bioinks.

The PhotoCol®-IRG, Methacrylated Collagen Hydrogel Kit is based upon purified type I bovine collagen methacrylate (ColMA), which when photocrosslinked provides a native-like 3D environment for cells. In addition to type I bovine collagen methacrylate, the kit includes the photoinitiator Irgacure 2959 and other reagents for users to easily fine tune their photocrosslinking experiments (i.e. altering hydrogel stiffness or gelling speeds). The methacrylated Type I collagen is produced from telo-peptide intact bovine collagen that has been modified by reacting with the protein's free amines, primarily the ϵ -amines groups of the lysine residues as well as the α -amines groups on the N-termini. > 20% of the total lysine residues of the collagen molecule have been methacrylated. The collagen is extracted from bovine hide and contains a high monomer content that was isolated from a closed herd and purified using controlled manufacturing processes.

Kit Components

The PhotoCol®-IRG, Methacrylated Collagen Hydrogel Kit (CC320) contains:

- 1) CC320-1 (Store at Room Temp): Neutralization Solution, 1 X 10 mL (CS226415)
- 2) CC320-2 (Store at 2-8°C): Methacrylated Collagen, 1 X 100 mg (CS226418), 20 mM Acetic Acid, 1 X 50 mL (CS226417), Irgacure Photoinitiator, 1 X 100 mg (CS226427).

Note: Do not freeze reagents.

Quality Control

Sterility (USP modified): No Growth
Endotoxin LAL: ≤ 10.0 EU/ml
Gel Formation: ≤ 40 min
Kinetic Gel Test: ≤ 40 min
Purity (SDS Page Silver Stain): $\geq 99\%$
Degree of Methacrylation: $\geq 20\%$
pH: 3.0-4.0
Osmolality: ≤ 40 mOsm/Kg H₂O
Young's Modulus: Characteristic
Electrophoretic Pattern: Characteristic

References

- 1) Isaacson, A, Swioklo, S. & Connon, C. J. 3D bioprinting of a corneal stroma equivalent. *Experimental Eye Research* 173, 188–193 (2018).
- 2) Maloney, E, et al. Immersion Bioprinting of Tumor Organoids in Multi-Well Plates for Increasing Chemotherapy Screening Throughput. *Micromachines* 11, 208 (2020).
- 3) Lelièvre, S, et al. Cell Culture and Coculture for Oncological Research in Appropriate Microenvironments. *Current Protocols Chem Biol.* 2019 Jun;11(2):e65.

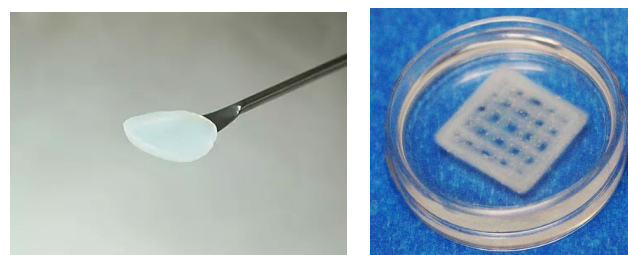


Figure 1. 3D printing of PhotoCol® Methacrylated Collagen Hydrogels can be used as native bioinks for tissue engineering bioprinting applications.

Instructions for Use

Note: Employ aseptic practices to maintain the sterility of the product throughout the preparation and handling of the collagen and other solutions. It is recommended that the collagen and other working solutions be chilled and kept on ice during the preparation of the collagen. Vortexing is not recommended at any step.

1. Add volume of 20 mM acetic acid to the lyophilized methacrylated collagen to achieve desired concentration. Recommend concentration(s) range from 3 to 8 mg/mL.
2. Mix on a shaker table or rotator plate at 2-8°C until fully solubilized or overnight. Avoid formation of air bubbles as possible. *Note: The higher concentrations of collagen will take longer to solubilize.*
3. Determine the volume of the neutralization solution (NS) to mix with the collagen. To achieve a final pH of 7.0 to 7.4, follow the guidelines below.

Collagen Concentration	Collagen Weight	Neutralization Solution Volume
3 mg/mL	1.0 g	100 µL
4 mg/mL	1.0 g	114 µL
6 mg/mL	1.0 g	120 µL
8 mg/mL	1.0 g	128 µL

4. Transfer the required volume of the neutralization solution (NS) into a sterile vessel or tube and chill.
5. If crosslinking is desired, first add 1 mL of neat methanol to the amber vial containing 100 mg of Irgacure, and vortex. (Irgacure only has a 2-week shelf life upon solubilizing. If you need the Irgacure to last longer, remove required amount and solubilize (10% solution). Calculate the volume of the photoinitiator required by multiplying the total volume required (collagen and neutralization solution) by 0.01.
6. Add the calculated volume of chilled photoinitiator to the volume of chilled neutralization solution (NS) and mix thoroughly.
7. Transfer the total volume of the chilled collagen into the chilled neutralization solution (NS)/photoinitiator. Mix quickly and thoroughly by pipetting or rotating a vessel or tube. *Note: Keep the collagen mixture chilled throughout this process. Check to ensure the pH is neutral. The high viscosity of this material can make it harder to mix.*
8. If desired, add dispersed chilled cells to the collagen mixture. Mix quickly and thoroughly by pipetting or rotating a vessel or tube. *Note: If air bubbles are a concern, allow to sit on ice until the bubbles come to the surface.*
9. Dispense the collagen mixture in the desired sterile plates or culture vessels.
10. Incubate at 37°C for > 30 minutes for gel formation.
11. If crosslinking is desired, place directly under a 365 nm UV light crosslinking source. Tunability of PhotoCol® product exposed to UV for 45, 90 and 600 seconds results respectively in a 22, 53, and 75% increase in gel stiffness while maintaining a stable gel stiffness over time. Longer exposure allows more crosslinking, though each cell type withstands different degrees of UV light and the free radicals (generated by the photoinitiator) that mediate crosslinking. *Note: The consistency and fidelity of crosslinking is improved by plating gels on glass-bottom substrates with good optical properties that produce minimal light scattering.*

■ antibodies ■ Multiplex products ■ biotools ■ cell culture ■ enzymes ■ kits ■ proteins/peptides ■ siRNA/cDNA products

Please visit www.milliporesigma.com for additional product information, test data and references

MilliporeSigma, 28820 Single Oak Drive, Temecula, CA 92590, USA 1-800-437-7500

Technical Support: T: 1-800-MILLIPORE (1-800-645-5476) • F: 1-800-437-7502

FOR RESEARCH USE ONLY. Not for use in diagnostic procedures. Not for human or animal consumption. Purchase of this Product does not include any right to resell or transfer, either as a stand-alone product or as a component of another product. Any use of this Product for purposes other than research is strictly prohibited.

MilliporeSigma®, the M mark, Upstate®, Chemicon®, Linco® and all other registered trademarks, unless specifically identified above in the text as belonging to a third party, are owned by Merck KGaA, Darmstadt, Germany. Copyright ©2008-2020 Merck KGaA, Darmstadt, Germany. All rights reserved.



We Buy 100% Certified
Renewable Energy

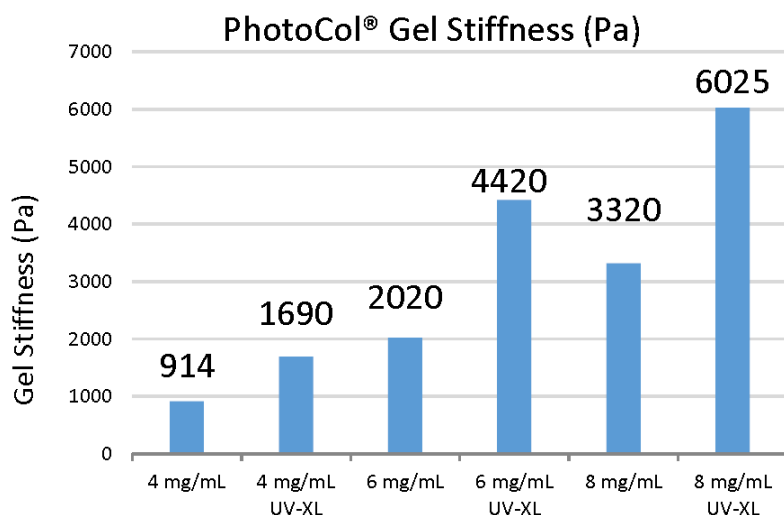


Figure 2. Gel stiffness of PhotoCol® at varying concentrations with and without UV crosslinking. Rheology testing was done on a Bohlin CVO-100 rheometer. Crosslinked collagen was exposed to 365 nm UV light for 5 minutes.

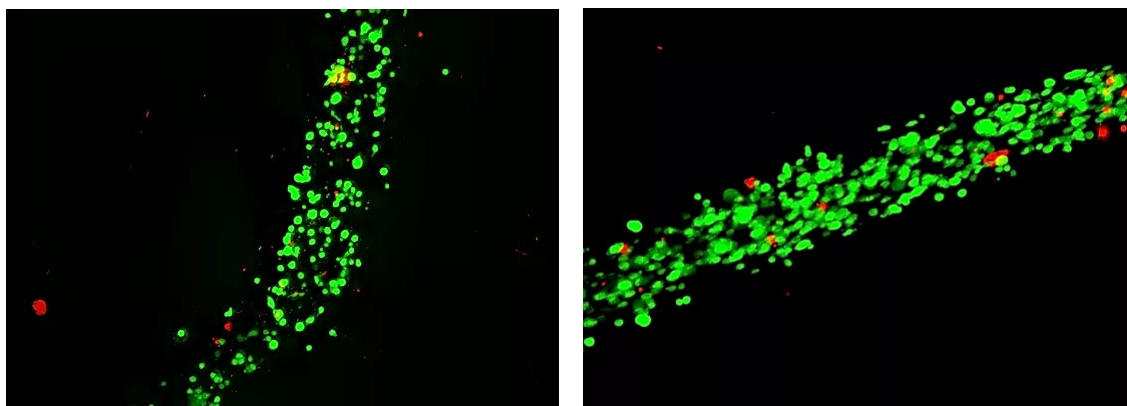


Figure 3. Bioprinting of human mesenchymal stem cells using PhotoCol®. PhotoCol® was used as a rapidly self-assembling type I collagen to form cross-linked hydrogels for tissue engineering. Human adipose derived mesenchymal stem cells (MSCs) (10×10^6 cells/ml) 1-day post bioprinting using PhotoCol® as the bioink material. Cell health was assayed using a live/dead assay (green= viable, red=dead). The majority (>90%) of the stem cells are viable when bioprinted using PhotoCol®.

PhotoCol® is a registered trademark of Advanced BioMatrix, Inc.

antibodies Multiplex products biotools cell culture enzymes kits proteins/peptides siRNA/cDNA products

Please visit www.milliporesigma.com for additional product information, test data and references

MilliporeSigma, 28820 Single Oak Drive, Temecula, CA 92590, USA 1-800-437-7500

Technical Support: T: 1-800-MILLIPORE (1-800-645-5476) • F: 1-800-437-7502

FOR RESEARCH USE ONLY. Not for use in diagnostic procedures. Not for human or animal consumption. Purchase of this Product does not include any right to resell or transfer, either as a stand-alone product or as a component of another product. Any use of this Product for purposes other than research is strictly prohibited.

MilliporeSigma®, the M mark, Upstate®, Chemicon®, Linco® and all other registered trademarks, unless specifically identified above in the text as belonging to a third party, are owned by Merck KGaA, Darmstadt, Germany. Copyright ©2008-2020 Merck KGaA, Darmstadt, Germany. All rights reserved.



We Buy 100% Certified Renewable Energy