

User Guide

MILLIPLEX[®] Feline Cytokine/Chemokine Magnetic Bead Panel

96-Well Plate Assay

FCYTMAG–20K-PMX FCYTMAG20KPX19BK

Introduction	3	Preparation of Serum Matrix.....	14
Principle.....	5	Preparation of Feline Cytokine/Chemokine Standard	15
Storage Conditions Upon Receipt	6	Immunoassay Procedure	18
Reagents Supplied.....	7	Plate Washing.....	20
Materials Required (not included)	9	Equipment Settings	21
Reagents.....	9	Quality Controls	22
Instrumentation/Materials.....	9	Assay Characteristics	23
Safety Precautions.....	10	Cross-Reactivity	23
Symbol Definitions.....	10	Assay Sensitivities	23
Technical Guidelines	11	Precision	24
Sample Collection and Storage	13	Accuracy	25
Preparation of Serum Samples	13	Troubleshooting	26
Preparation of Plasma Samples.....	13	Product Ordering.....	30
Preparation of Tissue Culture Supernatant	13	Well Map.....	31
Preparation of Reagents for Immunoassay	14	Notice	32
Preparation of Antibody-Immobilized Beads	14	Contact Information.....	32
Preparation of Quality Controls	14	Technical Assistance	32
Preparation of Wash Buffer	14	Standard Warranty	32
		Safety Data Sheets (SDS).....	32

For research use only. Not for use in diagnostic procedures.

00004429MAN Rev 03/25

2

Introduction

Cytokines are immunomodulatory polypeptides that play key roles in both adaptive and innate immune responses. A generic term, “cytokines” includes myokines, lymphokines, adipokines, etc. and have been divided into major classes: interleukins (IL), interferons (IFN), chemokines, colony stimulating factors (CSF) and tumor necrosis factors (TNF). Acting as regulators under both normal and pathological conditions, cytokines modulate the functional activities of individual cells and tissues. Consequently, research involving cytokines plays a significant role in achieving a deeper understanding of the immune system and its multi-faceted response to most antigens, especially those responses that make up the inflammatory process.

The MILLIPLEX® portfolio offers the broadest selection of analytes across a wide range of disease states and species. Once the analytes of interest have been identified, you can rely on the quality that we build into each kit to produce results you can trust. In addition to the assay characteristics listed in the protocol, other performance criteria evaluated during the verification process include: cross-reactivity, dilution linearity, kit stability, and sample behavior (for example, detectability and stability).

Each MILLIPLEX® panel and kit includes:

- Quality controls (QCs) provided to qualify assay performance
- Comparison of standard (calibrator) and QC lots to a reference lot to ensure lot-to-lot consistency
- Optimized serum matrix to mimic native analyte environment
- Detection antibody cocktails designed to yield consistent analyte profiles within panel

In addition, each panel and kit meets stringent manufacturing criteria to ensure batch-to-batch reproducibility. The MILLIPLEX® Feline Cytokine/Chemokine Magnetic Bead Panel thus enables you to focus on the therapeutic potential of the study of inflammatory and immune response in feline models. Coupled with the Luminex xMAP® platform in a magnetic bead format, you receive the advantage of ideal speed and sensitivity, allowing quantitative multiplex detection of dozens of analytes simultaneously, which can dramatically improve productivity.

The MILLIPLEX® Feline Cytokine/Chemokine Magnetic Bead Panel is part of the most versatile system available for inflammatory and immune response research. From our single to multiplex biomarker solutions, we partner with you to design, develop, analytically verify and build the most comprehensive library available for protein detection and quantitation.

MILLIPLEX® products offer you:

- A 19-plex premixed option
- A convenient “all-in-one” box format that gives you the assurance that you will have all the necessary reagents you need to run your assay.

The MILLIPLEX® Cytokine/Chemokine Magnetic Bead Panel is a 19-plex premixed kit to be used for the simultaneous quantification of all the following analytes in tissue/cell lysate and culture supernatant samples and serum or plasma samples: Fas, Flt-3L, GM-CSF, IFN γ , IL-1 β , IL-2, IL-4, IL-6, IL-8, IL-12 (p40), IL-13, IL-18, KC, MCP-1, PDGF-BB, RANTES, SCF, SDF-1, and TNF α .

For research use only. Not for use in diagnostic procedures.

Please read entire protocol before use.

It is important to use same assay incubation conditions throughout your study.

Principle

MILLIPLEX® products are based on the Luminex xMAP® technology — one of the fastest growing and most respected multiplex technologies offering applications throughout the life-sciences and capable of performing a variety of bioassays including immunoassays on the surface of fluorescent-coded magnetic beads known as MagPlex®-C microspheres.

- Luminex products use proprietary techniques to internally color-code microspheres with two fluorescent dyes. Through precise concentrations of these dyes, distinctly colored bead sets of 500-5.6 µm polystyrene microspheres or 80-6.45 µm magnetic microspheres can be created, each of which is coated with a specific capture antibody.
- After an analyte from a test sample is captured by the bead, a biotinylated detection antibody is introduced.
- The reaction mixture is then incubated with Streptavidin-PE conjugate, the reporter molecule, to complete the reaction on the surface of each microsphere.
- The following Luminex® instruments can be used to acquire and analyze data using two detection methods:
 - The Luminex® analyzers, Luminex® 200™, FLEXMAP 3D®, and xMAP® INTELLIFLEX, are flow cytometry-based instruments that integrate key xMAP® detection components, such as lasers, optics, advanced fluidics and high-speed digital signal processors.
 - The Luminex analyzer (MAGPIX®), a CCD-based instrument that integrates key xMAP® capture and detection components with the speed and efficiency of magnetic beads.
- Each individual microsphere is identified, and the result of its bioassay is quantified based on fluorescent reporter signals. We combine the streamlined data acquisition power of Luminex xPONENT® acquisition software with sophisticated analysis capabilities of the new MILLIPLEX® Analyst 5.1, integrating data acquisition and analysis seamlessly with all Luminex instruments.
- xMAP® INTELLIFLEX runs on INTELLIFLEX software for instrument control, run setup and generating high quality data with flexible output options. Data can be exported in xPONENT® style CSV files for compatibility with many existing analytical applications, or in the new, customizable INTELLIFLEX file format. The INTELLIFLEX file format is intended for flexibility and simplicity, allowing the user to freely select which data points to include and to reduce the time to analysis.

The capability of adding multiple conjugated beads to each sample results in the ability to obtain multiple results from each sample. Open-architecture xMAP® technology enables multiplexing of many types of bioassays reducing time, labor and costs over traditional methods.

Storage Conditions Upon Receipt

- Recommended storage for kit components is 2–8 °C.
- For long-term storage, freeze reconstituted standards and controls at ≤ -20 °C. Avoid multiple (> 2) freeze/thaw cycles.
- **DO NOT FREEZE** Antibody-Immobilized Beads, Detection Antibody, and Streptavidin-Phycoerythrin.

Reagents Supplied

Store all reagents at 2–8 °C.

Reagents	Volume	Quantity	Cat. No.
Feline Cytokine/Chemokine Standard	Lyophilized	1 vial	FCY-8020
Feline Cytokine/Chemokine Quality Controls 1 and 2	Lyophilized	2 vials	FCY-6020
Serum Matrix*	Lyophilized	2 vials	LMC-SD
Set of one 96-Well Plate with 2 sealers	-	1 set	-
Assay Buffer	30 mL	1 bottle	L-AB
10X Wash Buffer**	60 mL	1 bottle	L-WB
Feline Cytokine/Chemokine Detection Antibodies	3.2 mL	1 bottle	FCY-1020
Streptavidin-Phycoerythrin	3.2 mL	1 bottle	L-SAPE6
Mixing Bottle (not provided with premixed panel)	-	1 bottle	-

* Contains 0.08% Sodium azide

** Contains 0.05% ProClin™

Feline Cytokine/Chemokine Antibody-Immobilized Premixed Magnetic Beads

	Volume	Quantity	Cat. No.
Premixed 19-plex Beads	3.5 mL	1 bottle	FCYPMX19-MAG

The following Feline Cytokine/Chemokine Antibody-Immobilized Magnetic Beads are included within the Premixed 19-plex Beads FCYPMX19-MAG.

Feline Cytokine/Chemokine Antibody-Immobilized Magnetic Beads

Bead/Analyte Name	Luminex Magnetic Bead Region	19-Plex
		Magnetic Premixed Beads Available
Fas	13	✓
FIt-3L	15	✓
GM-CSF	18	✓
IFN γ	21	✓
IL-1 β	30	✓
IL-2	33	✓
PDGF-BB	36	✓
IL-12 (p40)	42	✓
IL-13	44	✓
IL-4	47	✓
IL-6	51	✓
IL-8	55	✓
KC	61	✓
SDF-1	64	✓
RANTES	67	✓
SCF	72	✓
MCP-1	74	✓
TNF α	76	✓
IL-18	78	✓

Materials Required (not included)

Reagents

MAGPIX® Drive Fluid PLUS (Cat. No. 40-50030), xMAP® Sheath Fluid PLUS (Cat. No. 40-50021), or xMAP® Sheath Concentrate PLUS (Cat. No. 40-50023)

Instrumentation/Materials

- Adjustable pipettes with tips capable of delivering 25 µL to 1000 µL
- Multichannel pipettes capable of delivering 5 µL to 50 µL, or 25 µL to 200 µL
- Reagent reservoirs
- Polypropylene microfuge tubes
- Rubber bands
- Aluminum foil
- Absorbent pads
- Laboratory vortex mixer
- Sonicator (Branson Ultrasonic Cleaner Model No. B200 or equivalent)
- Titer plate shaker (VWR® Microplate Shaker Cat. No. 12620-926 or equivalent)
- Luminex® 200™, HTS, FLEXMAP 3D®, MAGPIX® instrument with xPONENT® software, or xMAP® INTELLIFLEX instrument with INTELLIFLEX software by Luminex® Corporation
- Automatic Plate Washer for magnetic beads (BioTek® 405 LS and 405 TS, Cat. No. 40-094, 40-095, 40-096, 40-097 or equivalent) or Handheld Magnetic Separation Block (Cat. No. 40-285 or equivalent).

Note: If a plate washer or handheld magnetic separation block for magnetic beads is not available, one can use a microtiter filter plate (Cat. No. MX-PLATE) to run the assay using a vacuum filtration unit (Vacuum Manifold, Cat. No. MSMHTS00 or equivalent with Vacuum Pump, Cat. No. WP6111560 or equivalent).

Safety Precautions

- All blood components and biological materials should be handled as potentially hazardous. Follow universal precautions as established by the Centers for Disease Control and Prevention and by the Occupational Safety and Health Administration when handling and disposing of infectious agents.
- Sodium azide or ProClin™ has been added to some reagents as a preservative. Although the concentrations are low, Sodium azide and ProClin™ may react with lead and copper plumbing to form highly explosive metal azides. Dispose of unused contents and waste in accordance with international, federal, state, and local regulations.

Symbol Definitions

Ingredient	Cat. No.	Label	
Feline Cytokine Quality Control 1 & 2	FCY-6020	 	Warning. Harmful if swallowed. Toxic to aquatic life with long lasting effects. Avoid release to the environment.
Feline Cytokine Standard	FCY-8020	 	Warning. Harmful if swallowed. Toxic to aquatic life with long lasting effects. Avoid release to the environment.
Serum Matrix	LMC-SD	No Symbol Required	Harmful to aquatic life with long lasting effects. Avoid release to the environment.
Streptavidin-Phycoerythrin	L-SAPE6		Warning. Causes serious eye irritation. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.
10X Wash Buffer	L-WB		Warning. May cause an allergic skin reaction. Wear protective gloves. IF ON SKIN: Wash with plenty of soap and water.

Technical Guidelines

To obtain reliable and reproducible results, the operator should carefully read this entire manual and fully understand all aspects of each assay step before running the assay. The following notes should be reviewed and understood before the assay is set up.

- FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- Do not use beyond the expiration date on the label.
- Do not mix or substitute reagents with those from other lots or sources.
- The Antibody-Immobilized Beads are light sensitive and must be protected from light at all times. Cover the assay plate containing beads with opaque plate lid or aluminum foil during all incubation steps.
- It is important to allow all reagents to warm to room temperature (20-25 °C) before use in the assay.
- Incomplete washing can adversely affect the assay outcome. All washing must be performed with the Wash Buffer provided.
- The standards prepared by serial dilution must be used within 1 hour of preparation. Discard any unused standards except the standard stock which may be stored at ≤ -20 °C for 1 month and at ≤ -80 °C for greater than one month.
- If samples fall outside the dynamic range of the assay, further dilute the samples with the appropriate diluent and repeat the assay.
- Any unused mixed Antibody-Immobilized Beads may be stored in the Mixing Bottle at 2-8 °C for up to one month.
- During the preparation of the standard curve, make certain to mix the higher concentration well before making the next dilution. Use a new tip with each dilution.
- The plate should be read immediately after the assay is finished. If, however, the plate cannot be read immediately, seal the plate, cover with aluminum foil or an opaque lid, and store the plate at 2-8 °C for up to 24 hours. Prior to reading, agitate the plate on the plate shaker at room temperature for 10 minutes. Delay in reading a plate may result in decreased sensitivity for some analytes.
- The titer plate shaker should be set at a speed to provide maximum orbital mixing without splashing of liquid outside the wells. For the recommended plate shaker, this would be a setting of 5-7 which is approximately 500-800 rpm.
- Ensure that the needle probe is clean. This may be achieved by sonication and/or alcohol flushes.

-
- When reading the assay on the Luminex 200™ instrument, adjust probe height according to the protocols recommended by Luminex to the kit solid plate or to the recommended filter plates using 3 alignment discs. When reading the assay on the MAGPIX® instrument, adjust probe height according to the protocols recommended by Luminex to the kit solid plate or to the recommended filter plates using 2 alignment discs. When reading the assay on the FLEXMAP 3D® instrument, adjust probe height according to the protocols recommended by Luminex to the kit solid plate using 1 alignment disc.
 - For the FLEXMAP 3D® instrument, when using the solid plate in the kit, the final resuspension should be with 150 µL Sheath Fluid PLUS in each well and 75 µL should be aspirated.
 - For the xMAP® INTELLIFLEX instrument, adjust probe height based on the type of plate you are using, place an alignment disk or an alignment sphere in the well according to the protocol recommended by Luminex®.
 - For cell culture supernatants or tissue extraction, use the culture or extraction medium as the matrix solution in background, standard curve and control wells. If samples are diluted in assay buffer, use the assay buffer as matrix.
 - For serum/plasma samples that require dilution, use the serum matrix provided in the kit.
 - For cell/tissue homogenate, the final cell or tissue homogenate should be prepared in a buffer that has a neutral pH, contains minimal detergents or strong denaturing detergents, and has an ionic strength close to physiological concentration. Avoid debris, lipids, and cell/tissue aggregates. Centrifuge samples before use.
 - Vortex all reagents well before adding to plate.

Sample Collection and Storage

Preparation of Serum Samples

- Allow the blood to clot for at least 30 minutes before centrifugation for 10 minutes at 1000 x g. Remove serum and assay immediately or aliquot and store samples at ≤ -20 °C.
- Avoid multiple (> 2) freeze/thaw cycles.
- When using frozen samples, it is recommended to thaw the samples completely, mix well by vortexing and centrifuge prior to use in the assay to remove particulates.
- Neat serum samples are used. If dilution is required, use Serum Matrix provided as diluent.

Preparation of Plasma Samples

- Plasma collection using EDTA as an anti-coagulant is recommended. Centrifuge for 10 minutes at 1000 x g within 30 minutes of blood collection. Remove plasma and assay immediately or aliquot and store samples at ≤ -20 °C.
- Avoid multiple (> 2) freeze/thaw cycles.
- When using frozen samples, it is recommended to thaw the samples completely, mix well by vortexing and centrifuge prior to use in the assay to remove particulates.
- Neat plasma samples are used. If dilution is required, use Serum Matrix provided as diluent.

Preparation of Tissue Culture Supernatant

- Centrifuge the sample to remove debris and assay immediately or aliquot and store samples at ≤ -20 °C.
- Avoid multiple (> 2) freeze/thaw cycles.
- Tissue culture supernatant may require a dilution with an appropriate control medium prior to assay. Tissue/cell extracts should be done in neutral buffers containing reagents and conditions that do not interfere with assay performance. Excess concentrations of detergent, salt, denaturants, high or low pH, etc. will negatively affect the assay. Organic solvents should be avoided. The tissue/cell extract samples should be free of particles such as cells or tissue debris.

NOTE:

- A maximum of 10 μ L per well of neat serum or plasma can be used. A maximum of 25 μ L of tissue culture or other media may also be used and dilution must be determined by the user.
- All samples must be stored in polypropylene tubes. **DO NOT STORE SAMPLES IN GLASS.**
- Avoid debris, lipids and cells when using samples with gross hemolysis or lipemia.
- Care must be taken when using heparin as an anti-coagulant since an excess of heparin will provide falsely high values. Use no more than 10 IU heparin per mL of blood collected.

Preparation of Reagents for Immunoassay

Preparation of Antibody-Immobilized Beads

For the premixed beads, sonicate the premixed bead bottle 30 seconds and then vortex for 1 minute before use.

Preparation of Quality Controls

Before use, reconstitute Quality Control 1 and Quality Control 2 with 250 μ L deionized water. Invert the vial several times to mix and vortex. Allow the vial to sit for 5-10 minutes. Transfer the reconstituted Quality Control 1 and Quality Control 2 into two polypropylene microfuge tubes and set in an ice bath. These should be added to the plate within 1 hour of reconstitution. Unused portion may be stored at ≤ -20 °C for up to one month.

Preparation of Wash Buffer

Bring the 10X Wash Buffer to room temperature and mix to bring all salts into solution. Dilute 60 mL of 10X Wash Buffer with 540 mL deionized water. Store the unused portion at 2-8 °C for up to one month.

Preparation of Serum Matrix

This step is required for serum or plasma samples only.

Add 1 mL deionized water to each of the bottles containing lyophilized Serum Matrix. Mix well. Allow at least 5 minutes for complete reconstitution. Leftover reconstituted Serum Matrix should be stored at ≤ -20 °C for up to one month.


Preparation of Feline Cytokine/Chemokine Standard

1. Prior to use, reconstitute the Feline Cytokine/Chemokine Standard with 250 μL deionized water. Allow the vial to sit for 5-10 minutes. Transfer the reconstituted standard to a polypropylene microfuge tube. Invert the vial several times to mix. Vortex the vial for 10 seconds. This will be used as "Standard 7"; This reconstituted standard and the serially diluted standards in the following steps should be set in an ice bath, during the assay procedure. These need to be added to the plate within 1 hour of preparation. The unused portion of the "Standard 7" may be stored at $\leq -20\text{ }^{\circ}\text{C}$ for up to one month.
2. Preparation of Working Standards
Label 6 polypropylene microfuge tubes Standard 6, Standard 5, Standard 4, Standard 3, Standard 2, and Standard 1. Add 200 μL of Assay Buffer to each of the 6 tubes. Prepare serial dilutions by adding 100 μL of the reconstituted standard 7 to the Standard 6 tube, mix well and transfer 100 μL of the standard 6 to the Standard 5 tube, mix well and transfer 100 μL of the standard 5 to the Standard 4 tube, mix well and transfer 100 μL of the standard 4 to Standard 3 tube, mix well and transfer 100 μL of the standard 3 to the Standard 2 tube, mix well and transfer 100 μL of the standard 2 to the Standard 1 tube and mix well. The 0 ng/mL standard (Background) will be Assay Buffer.

Standard No.	Add Deionized Water (μL)	Add Standard (volume)
Original (Std 7)	250	0

Standard No.	Add Assay Buffer (μL)	Add Standard (volume)
Standard 6	200	100 μL of Std 7
Standard 5	200	100 μL of Std 6
Standard 4	200	100 μL of Std 5
Standard 3	200	100 μL of Std 4
Standard 2	200	100 μL of Std 3
Standard 1	200	100 μL of Std 2

Preparation of Standards

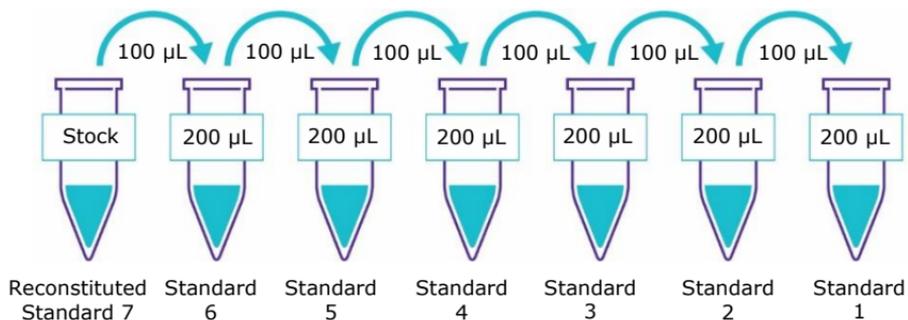

Standard	KC, RANTES (pg/mL)	Fit-3L (pg/mL)	Fas, GM-CSF, IL-2, IL-8, IL- 12p40, IL-13 (pg/mL)	IL-1β (pg/mL)
Standard 1	2.7	6.9	13.7	34.3
Standard 2	8.2	20.6	41.2	102.9
Standard 3	24.7	61.7	123.5	308.6
Standard 4	74.1	185.2	370.4	925.9
Standard 5	222.2	555.6	1,111	2,778
Standard 6	666.7	1,667	3,333	8,333
Standard 7	2,000	5,000	10,000	25,000

Standard	IFNγ, IL-4, IL-6, SCF, SDF-1 (pg/mL)	IL-18, MCP-1 (pg/mL)	PDGF-BB (pg/mL)	TNFα (pg/mL)
Standard 1	68.6	137.2	274.3	Refer to QC analysis sheet for exact concentration
Standard 2	205.8	411.5	823	
Standard 3	617.3	1,235	2,469	
Standard 4	1,852	3,704	7,407	
Standard 5	5,556	11,111	22,222	
Standard 6	16,667	33,333	66,667	
Standard 7	50,000	100,000	200,000	

Immunoassay Procedure

- Prior to beginning this assay, it is imperative to read this protocol completely and to thoroughly understand the Technical Guidelines.
- Allow all reagents to warm to room temperature (20-25 °C) before use in the assay.
- Diagram the placement of Standards [0 (Background), Standard 1, 2, 3, 4, 5, 6, and 7], Controls 1 and 2, and Samples on Well Map Worksheet in a vertical configuration.
(Note: Most instruments will only read the 96-well plate vertically by default.) It is recommended to run the assay in duplicate.
- If using a filter plate, set the filter plate on a plate holder at all times during reagent dispensing and incubation steps so that the bottom of the plate does not touch any surface.

1. Add 200 μL of Assay Buffer into each well of the plate. Seal and mix on a plate shaker for 10 minutes at room temperature (20-25 °C).
2. Decant Assay Buffer and remove the residual amount from all wells by inverting the plate and tapping it smartly onto absorbent towels several times.
3. Add 25 μL of each Standard or Control into the appropriate wells. Assay Buffer should be used for 0 ng/mL standard (Background).
4. Add 25 μL of Assay Buffer to the sample wells.
5. Add 25 μL of appropriate matrix solution to the background, standards, and control wells. When assaying serum or plasma, use the Serum Matrix. When assaying tissue culture or other supernatant, use proper control culture medium as the matrix solution.
6. If assaying serum/plasma samples, add 15 μL of Serum Matrix to the serum/plasma sample wells.
7. If assaying serum/plasma, add 10 μL of Sample (Neat) into the appropriate wells or 25 μL of tissue culture/supernatant to sample wells.
8. Vortex Mixing Bottle and add 25 μL of the Mixed or Premixed Beads to each well. (Note: During addition of Beads, shake bead bottle intermittently to avoid settling.)

Add 200 μL Assay Buffer per well



Shake 10 min, RT

Decant

- Add 25 μL Standard or Control to appropriate wells
- Add 25 μL Assay Buffer to background and sample wells
- Add 25 μL appropriate matrix solution to background, standards, and control wells
- Add 15 μL of Serum Matrix to Serum/Plasma sample wells
- Add 10 μL Neat Samples to sample wells or 25 μL tissue culture/supernatant to sample wells
- Add 25 μL Beads to each well

9. Seal the plate with a plate sealer. Wrap the plate with foil and incubate with agitation on a plate shaker for 16-18 hours at 4 °C.
10. Gently remove well contents and wash plate 2 times following instructions listed in the Plate Washing section.
11. Add 25 µL of Detection Antibodies into each well.
(Note: Allow the Detection Antibodies to warm to room temperature prior to addition.)
12. Seal, cover with foil and incubate with agitation on a plate shaker for 1 hour at room temperature (20-25 °C). DO NOT ASPIRATE AFTER INCUBATION.
13. Add 25 µL Streptavidin-Phycoerythrin to each well containing the 25 µL of Detection Antibodies.
14. Seal, cover with foil and incubate with agitation on a plate shaker for 30 minutes at room temperature (20-25 °C).
15. Gently remove well contents and wash plate 2 times following instructions listed in the Plate Washing section.
16. Add 100 µL of Sheath Fluid PLUS (or Drive Fluid PLUS if using MAGPIX®) to all wells. Resuspend the beads on a plate shaker for 5 minutes.
17. Run plate on Luminex® 200™, HTS, FLEXMAP 3D®, MAGPIX® instrument with xPONENT® software or xMAP® INTELLIFLEX instrument with INTELLIFLEX software.
18. Save and analyze the Median Fluorescent Intensity (MFI) data using a 5-parameter logistic or spline curve-fitting method for calculating analyte concentrations in samples.
(Note: For serum/plasma samples, multiply by a 2.5 dilution factor. For diluted samples, multiply the calculated concentration by the dilution factor.)



Incubate 16-18 hours at 4 °C



Remove well contents and wash 2X with 200 µL Wash Buffer

Add 25 µL Detection Antibodies per well



Incubate 1 hour at RT

Do Not Aspirate

Add 25 µL Streptavidin-Phycoerythrin per well



Incubate for 30 minutes at RT

Remove well contents and wash 2X with 200 µL Wash Buffer

Add 100 µL Sheath Fluid PLUS or Drive Fluid PLUS per well

Read on Luminex (50 µL, 50 beads per bead set)

Plate Washing

Solid Plate

If using a solid plate, use either a handheld magnet or magnetic plate washer.

- Handheld magnet (Cat. No. 40-285)
Rest plate on magnet for 60 seconds to allow complete settling of magnetic beads. Remove well contents by gently decanting the plate in an appropriate waste receptacle and gently tapping on absorbent pads to remove residual liquid. Wash plate with 200 μ L of Wash Buffer by removing plate from magnet, adding Wash Buffer, shaking for 30 seconds, reattaching to magnet, letting beads settle for 60 seconds and removing well contents as previously described after each wash. Repeat wash steps as recommended in Assay Procedure.
- Magnetic plate washer (Cat. No. 40-094, 40-095, 40-096 and 40-097)
Please refer to specific automatic plate washer manual for appropriate equipment settings. Please note that after the final aspiration, there will be approximately 25 μ L of residual wash buffer in each well. This is expected when using the BioTek[®]
• plate washer and this volume does not need to be aspirated from the plate.

If using an automatic plate washer other than BioTek[®] 405 LS or 405 TS, please refer to the manufacturer's recommendations for programming instructions.

Filter Plate (Cat. No. MX-PLATE)

If using a filter plate, use a vacuum filtration manifold to remove well contents. Wash plate with 200 μ L/well of Wash Buffer, removing Wash Buffer by vacuum filtration after each wash. Repeat wash steps as recommended in the Assay Procedure.

Equipment Settings

Luminex® 200™, HTS, FLEXMAP 3D®, MAGPIX® instruments with xPONENT® software and xMAP® INTELLIFLEX instrument with INTELLIFLEX software:

These specifications are for the above listed instruments and software. Luminex® instruments with other software (for example, MasterPlex®, StarStation, LiquiChip, Bio-Plex® Manager™, LABScan™100) would need to follow instrument instructions for gate settings and additional specifications from the vendors for reading Luminex® magnetic beads.

For magnetic bead assays, each instrument must be calibrated and performance verified with the indicated calibration and verification kits.

Instrument	Calibration Kit	Verification Kit
Luminex® 200™ and HTS	xPONENT® 3.1 compatible Calibration Kit (Cat. No. LX2R-CAL-K25)	Performance Verification Kit (Cat. No. LX2R-PVER-K25)
FLEXMAP 3D®	FLEXMAP 3D® Calibrator Kit (Cat. No. F3D-CAL-K25)	FLEXMAP 3D® Performance Verification Kit (Cat. No. F3D-PVER-K25)
xMAP® INTELLIFLEX	xMAP® INTELLIFLEX Calibration Kit (Cat. No. IFX-CAL-K20)	xMAP® INTELLIFLEX Performance Verification Kit (Cat. No. IFX-PVER-K20)
MAGPIX®	MAGPIX® Calibration Kit (Cat. No. MPX-CAL-K25)	MAGPIX® Performance Verification Kit (Cat. No. MPX-PVER-K25)

NOTE: When setting up a Protocol using the xPONENT® software, you must select MagPlex® as the Bead Type in the Acquisition settings.

NOTE: These assays cannot be run on any instruments using Luminex IS 2.3 or Luminex 1.7 software.

The Luminex probe height must be adjusted to the plate provided in the kit. Please use Cat. No. MAG-PLATE, if additional plates are required for this purpose.

Events	50, per bead
Sample Size	50 μ L
Gate Settings	8,000 to 15,000
Reporter Gain	Default (low PMT)
Time Out	60 seconds
Bead Set	19-Plex Beads
	<hr/>
Fas	13
Flt-3L	15
GM-CSF	18
IFN γ	21
IL-1 β	30
IL-2	33
PDGF-BB	36
IL-12p40	42
IL-13	44
IL-4	47
IL-6	51
IL-8	55
KC	61
SDF-1	64
RANTES	67
SCF	72
MCP-1	74
TNF α	76
IL-18	78

Quality Controls

The ranges for each analyte in Quality Control 1 and 2 are provided on the card insert or can be located at our website SigmaAldrich.com using the catalogue number as the keyword.

Assay Characteristics

Cross-Reactivity

There was no or negligible cross-reactivity between the antibodies for an analyte and any of the other analytes in this panel.

Assay Sensitivities (minimum detectable concentrations, pg/mL)

Minimum Detectable Concentration (MinDC) is calculated using MILLIPLEX® Analyst 5.1. It measures the true limits of detection for an assay by mathematically determining what the empirical MinDC would be if an infinite number of standard concentrations were run for the assay under the same conditions.

Overnight Protocol (n = 8 Assays)

Analyte	MinDC (pg/mL)	MinDC+2SD (pg/mL)
Fas	4	5
Flt-3L	4	6
GM-CSF	6	11
IFN γ	43	69
IL-1 β	14	23
IL-2	4	9
PDGF-BB	198	296
IL-12p40	9	15
IL-13	7	12
IL-4	30	66
IL-6	25	45
IL-8	7	13
KC	1	3
SDF-1	97	209
RANTES	1	3
SCF	48	76
MCP-1	164	494
TNF α	7	12
IL-18	30	46

Precision

Intra-assay precision is generated from the mean of the %CVs from 8 reportable results across two different concentrations of analytes in a single assay. Inter-assay precision is generated from the mean of the %CVs across two different concentrations of analytes across 5 different assays.

Analyte	Intra-assay %CV	Inter-assay %CV
Fas	< 10	< 15
Flt-3L	< 10	< 15
GM-CSF	< 10	< 15
IFN γ	< 10	< 15
IL-1 β	< 10	< 15
IL-2	< 10	< 15
PDGF-BB	< 10	< 15
IL-12p40	< 10	< 15
IL-13	< 10	< 15
IL-4	< 10	< 15
IL-6	< 10	< 15
IL-8	< 10	< 15
KC	< 10	< 15
SDF-1	< 15	< 15
RANTES	< 10	< 15
SCF	< 10	< 15
MCP-1	< 10	< 15
TNF α	< 10	< 17
IL-18	< 10	< 15

Accuracy

Spike Recovery: The data represent mean percent recovery of spiked standards ranging from low, medium, and high concentration in serum matrices (n=4).

Analyte	% Recovery in Serum Matrix
Fas	104
Flt-3L	104
GM-CSF	105
IFN γ	113
IL-1 β	105
IL-2	98
PDGF-BB	114
IL-12p40	105
IL-13	98
IL-4	108
IL-6	106
IL-8	109
KC	99
SDF-1	105
RANTES	104
SCF	104
MCP-1	96
TNF α	100
IL-18	96

Troubleshooting

Problem	Probable Cause	Solution
	Plate washer aspirate height set too low	Adjust aspiration height according to manufacturers' instructions.
	Bead mix prepared inappropriately	Sonicate bead vials and vortex just prior to adding to bead mix bottle according to protocol. Agitate bead mix intermittently in reservoir while pipetting this into the plate.
	Samples cause interference due to particulate matter or viscosity	See above. Also sample probe may need to be cleaned with alcohol flushes, back flushes and washes; or if needed, probe should be removed and sonicated.
Insufficient bead count		When reading the assay on the Luminex 200™ instrument, adjust probe height to the kit solid plate or to the recommended filter plates using 3 alignment discs. When reading the assay on the MAGPIX® instrument, adjust probe height to the kit solid plate or to the recommended filter plates using 2 alignment discs. When reading the assay on the FLEXMAP 3D® instrument, adjust probe height to the kit solid plate using 1 alignment disc. For the FLEXMAP 3D® instrument, when using the solid plate in the kit, the final resuspension should be with 150 µL Sheath Fluid in each well and 75 µL should be aspirated.
	Probe height not adjusted correctly	When reading the assay on the xMAP® INTELLIFLEX instrument, adjust probe height based on the type of plate you are using, place an alignment disk or an alignment sphere in the well according to the protocol recommended by Luminex®.
Background is too high	Background wells were contaminated	Avoid cross-well contamination by using sealer appropriately and pipetting with multichannel pipettes without touching reagent in plate.
	Matrix used has endogenous analyte or interference	Check matrix ingredients for cross-reacting components (for example, interleukin modified tissue culture medium).
	Insufficient washes	Increase number of washes.

Problem	Probable Cause	Solution
Beads not in region or gate	Luminex instrument not calibrated correctly or recently	Calibrate Luminex instrument based on manufacturer's instructions, at least once a week or if temperature has changed by >3 °C.
	Gate settings not adjusted correctly	Some Luminex instruments (for example, Bio-Plex®) require different gate settings than those described in the kit protocol. Use instrument default settings.
	Wrong bead regions in protocol template	Check kit protocol for correct bead regions or analyze selection.
	Incorrect sample type used	Samples containing organic solvents or if highly viscous should be diluted or dialyzed as required.
	Instrument not washed or primed	Prime the Luminex instrument 4 times to rid it of air bubbles, wash 4 times with sheath fluid or water if there is any remnant alcohol or sanitizing liquid.
Signal for whole plate is same as background	Incorrect or no Detection Antibody was added	Add appropriate Detection Antibody and continue.
	Streptavidin-Phycoerythrin was not added	Add Streptavidin-Phycoerythrin according to protocol. If Detection Antibody has already been removed, sensitivity may be low.
Low signal for standard curve	Detection Antibody may have been removed prior to adding Streptavidin-Phycoerythrin	May need to repeat assay if desired sensitivity not achieved.
	Incubations done at inappropriate temperatures, timings or agitation	Assay conditions need to be checked.
Signals too high, standard curves are saturated	Calibration target value set too high	With some Luminex instruments (for example, Bio-Plex®) default target setting for RP1 calibrator is set at high PMT. Use low target value for calibration and reanalyze plate.
	Plate incubation was too long with standard curve and samples	Use shorter incubation time.

Problem	Probable Cause	Solution
Sample readings are out of range	Samples contain no or below detectable levels of analyte	If below detectable levels, it may be possible to use higher sample volume. Check with technical support for appropriate protocol modifications.
	Samples contain analyte concentrations higher than highest standard point	Samples may require dilution and reanalysis for just that particular analyte.
	Standard curve was saturated at higher end of curve	See above.
High variation in samples and/or standards	Multichannel pipette may not be calibrated	Calibrate pipettes.
	Plate washing was not uniform	Confirm all reagents are removed completely in all wash steps.
	Samples may have high particulate matter or other interfering substances	See above.
	Plate agitation was insufficient	Plate should be agitated during all incubation steps using an orbital plate shaker at a speed where beads are in constant motion without causing splashing.
	Cross-well contamination	Check when reusing plate sealer that no reagent has touched sealer. Care should be taken when using same pipette tips that are used for reagent additions and that pipette tip does not touch reagent in plate.

FOR FILTER PLATES ONLY

Problem	Probable Cause	Solution
Filter plate will not vacuum	Vacuum pressure is insufficient	Increase vacuum pressure such that 0.2 mL buffer can be suctioned in 3-5 seconds.
	Samples have insoluble particles	Centrifuge samples just prior to assay set-up and use supernatant.
	High lipid concentration	After centrifugation, remove lipid layer and use supernatant.
Plate leaked	Vacuum pressure too high	Adjust vacuum pressure such that 0.2 mL buffer can be suctioned in 3-5 seconds. May need to transfer contents to a new (blocked) plate and continue.
	Plate set directly on table or absorbent towels during incubations or reagent additions	Set plate on plate holder or raised edge so bottom of filter is not touching any surface.
	Insufficient blotting of filter plate bottom causing wicking	Blot the bottom of the filter plate well with absorbent towels after each wash step.
	Pipette touching plate filter during additions	Pipette to the side of plate.
	Probe height not adjusted correctly	Adjust probe to 3 alignment discs in well H6.
	Sample too viscous	May need to dilute sample.

Product Ordering

Replacement Reagents	Cat. No.
Feline Cytokine/Chemokine Premix 19 plex Beads	FCYPMX19-MAG
Feline Cytokine/Chemokine Panel Standard	FCY-8020
Feline Cytokine/Chemokine Panel Quality Controls	FCY-6020
Serum Matrix	LMC-SD
Feline Cytokine/Chemokine Panel Detection Antibodies	FCY-1020
Assay Buffer	L-AB
Set of two 96-Well plates with sealers	MAG-PLATE
10X Wash Buffer	L-WB
Streptavidin-Phycoerythrin	L-SAPE6
Feline Cytokine/Chemokine 19 Plex Premixed Magnetic Bead Panel – BULK PACKAGING	FCYTMAG20KPX19BK

Well Map

	1	2	3	4	5	6	7	8	9	10	11	12
A	0 ng/mL Standard (Background)	Standard 4	QC-1 Control	Etc.								
B	0 ng/mL Standard (Background)	Standard 4	QC-1 Control									
C	Standard 1	Standard 5	QC-2 Control									
D	Standard 1	Standard 5	QC-2 Control									
E	Standard 2	Standard 6	Sample 1									
F	Standard 2	Standard 6	Sample 1									
G	Standard 3	Standard 7	Sample 2									
H	Standard 3	Standard 7	Sample 2									

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

By purchasing this product, which contains fluorescently labeled microsphere beads authorized by Luminex® Corporation ("Luminex®"), you, the customer, acquire the right under the Luminex® patent rights, if any, to use this product or any portion of this product, including without limitation the microsphere beads contained herein, only with the Luminex® laser based fluorescent analytical test instrumentation marketed under the name of Luminex® 200™, HTS, FLEXMAP 3D®, xMAP® INTELLIFLEX, MAGPIX®.

Contact Information

For the location of the office nearest you, go to [SigmaAldrich.com/offices](https://www.SigmaAldrich.com/offices).

Technical Assistance

Visit the tech service page at [SigmaAldrich.com/techservice](https://www.SigmaAldrich.com/techservice).

Standard Warranty

The applicable warranty for the products listed in this publication may be found at [SigmaAldrich.com/terms](https://www.SigmaAldrich.com/terms).

Safety Data Sheets (SDS)

Safety Data Sheets are available on the product page at [SigmaAldrich.com](https://www.SigmaAldrich.com).

The life science business of Merck operates
as MilliporeSigma in the U.S. and Canada.

Merck, Millipore, Milliplex and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2022-2025 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

For research use only. Not for use in diagnostic procedures.

The Merck logo is displayed in a bold, blue, sans-serif font. The letters are closely spaced and have a slightly irregular, hand-drawn appearance.