

Product Information

α -Amanitin

From *Amanita phalloides***A2263**Synonym: α -Amatoxin

CAS Number: 23109-05-9

Physical Description

Appearance

White to light yellow powder

Molecular formula

 $C_{39}H_{54}N_{10}O_{14}S$

Molecular weight

919.0

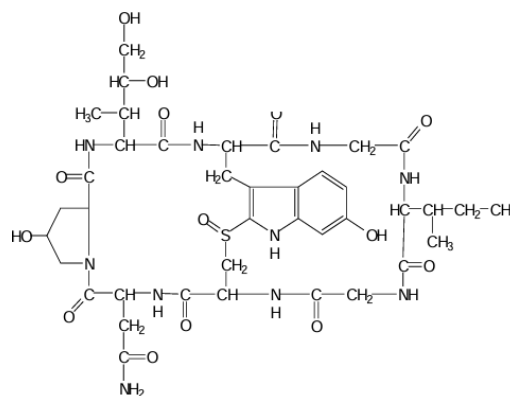
Melting point

254-255 °C

 ϵ^M (310 nm) = 13,500 (water)¹

Structure

α -amanitin differs structurally from β -amanitin in that it has an amide rather than a carboxyl group.

**Figure 1:** α -Amanitin.

Storage/Stability as Supplied

When stored frozen in the dark, retained samples were found to be greater than 99% by HPLC after two years.

Solubility/Stability of Solutions

We do not routinely test α -amanitin for solubility due to the hazardous nature of the material. It does dissolve at 1 mg/mL in water when tested by HPLC. It is reportedly soluble in water, methanol and ethanol.²

It is suggested to prepare a stock solution at 1 mg of α -amanitin per mL of water; this stock solution should be kept frozen in aliquots (-20 °C or below). Dilution with the appropriate buffer should be done immediately before use. The solution can be checked for intactness of α -amanitin by UV spectrometry, where it should show the typical absorption spectrum³ (ϵ^M = 13,500 at 310 nm in water).¹

α -Amanitin is destroyed by concentrated acids (sulfuric acid or hydrochloric acid) and by concentrated sodium hydroxide. Glassware may be quickly decontaminated by washing in sulfuric acid/chromic oxide cleaning solution.³

This product is extremely toxic. Be aware of the risks and be familiar with safety procedures before you use this product.

General Remarks

Among the various families of mushroom, the "true poisonous mushrooms" are of the genus *Amanita*. Not all the members of this genus are highly poisonous, but *Amanita phalloides*, the "green death cap" mushroom, contains a variety of compounds that are toxic for mammalian organisms. In this group, amatoxins and phallotoxins are two families of cyclic peptides that act on cytoplasmic or nuclear cell components. The amatoxins, one of which is α -amanitin, inhibit transcription of eukaryotic cells when present in nanomolar concentrations.⁴

In molecular biology applications, α -amanitin has become a tool in probing RNA-synthesis, since it has been shown to be a specific inhibitor of RNA Polymerase II.⁵⁻⁹

References

1. Wieland, T. and Faulstich, H., Critical Reviews in Biochemistry, 5, 185-260 (1978).
2. Data for Biochemical Research, 3rd Ed., Dawson, R.M.C. et al., eds. (Oxford Press, 1986), 257258.
3. Supplier information.
4. Wieland, T., Science, 159, 946-952 (1968). "Poisonous Principles of Mushrooms of the Genus *Amanita*."
5. Keding, C. et al., Biochem. Biophys. Res. Commun., 38, 165-171 (1970). " α -Amanitin: A Specific Inhibitor..."
6. Seifart, K.H. et al., Arch. Biochem. Biophys., 151, 519-532 (1972). "Multiple RNA Polymerase Species from Rat Liver Tissue..."
7. Huet, J. et al., Proc. Natl. Acad. Sci., 72, 3034-3038 (1975). "Dissociation of two polypeptide chains from yeast RNA polymerase A."
8. Schultz, L.D. and Hall, B.D., Proc. Natl. Acad. Sci., 73, 1029-1033 (1976). " α -Amanitin sensitivity...distinguish between RNA polymerases I and III."
9. Austoker, J.L., Cell, 3, 227-234 (1974). "DNA-Dependent RNA Polymerase...sensitive to... α Amanitin."

Review Articles

1. Faulstich, H., Progress in Molecular & Subcellular Biology, 7, 88-134 (1980). "The Amatoxins."
2. Zahler, A.M. and Prescott, D.M., Nucleic Acids Research, 17, 6299-6317 (1989). "DNA Primase and the replication of the telomeres in *Oxytricha nova*."
3. Wieland, T., Int. J. Peptide Protein Res., 22, 257-276 (1983). "The toxic peptides from *Amanita* mushrooms."
4. Culliton, B.J., Science, 185, 600-601 (1974). "The Destroying Angel: A Story of a Search for an Antidote."

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Technical Assistance

Visit the tech service page at SigmaAldrich.com/techservice.

Terms and Conditions of Sale

Warranty, use restrictions, and other conditions of sale may be found at SigmaAldrich.com/terms.

Contact Information

For the location of the office nearest you, go to SigmaAldrich.com/offices.

The life science business of Merck operates
as MilliporeSigma in the U.S. and Canada.

Merck and Sigma-Aldrich are trademarks of Merck or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 1996-2025 Merck and/or its affiliates. All Rights Reserved.
23263892 Rev 09/25

