

**Rat Adipocyte**

**96 Well Plate Assay**

**Cat. # RADPCYT-82K**

# MILLIPLEX<sup>®</sup> MAP

## RAT ADIPOCYTE Panel 96 Well Plate Assay

### #RADPCYT-82K

<b><u>TABLE OF CONTENTS</u></b>	<b><u>PAGE</u></b>
Introduction	2
Principle	3
Storage Conditions Upon Receipt	3
Reagents Supplied	4
Materials Required But Not Provided	5
Safety Precautions	5
Technical Guidelines	5
Sample Collection And Storage	7
Preparation of Reagents for Immunoassay	8
Immunoassay Procedure	10
Equipment Settings	12
Quality Controls	12
Assay Characteristics	13
Troubleshooting Guide	16
Replacement Reagents	18
Ordering Information	19
Well Map	20

### **For Research Use Only. Not for Use in Diagnostic Procedures.**

By purchasing this product, which contains fluorescently labeled microsphere beads authorized by Luminex Corporation ("Luminex"), you, the customer, acquire the right under Luminex's patent rights, if any, to use this product or any portion of this product, including without limitation the microsphere beads contained herein, only with Luminex's laser based fluorescent analytical test instrumentation marketed under the name of Luminex 100<sup>™</sup> IS, 200<sup>™</sup>, HTS.

## INTRODUCTION

In the past the role of white adipose tissue was thought to be limited to energy storage and internal organ protection. The discovery of the secretion of leptin from adipocytes in 1994 led to the recognition that white adipose tissue is involved in a variety of metabolic and physiological processes. Adipocytes secrete a number of hormones called adipokines with functions that include appetite and energy balance, insulin sensitivity, and lipid metabolism.

EMD Millipore recognizes the need to better understand the role that adipocytes and adipokines play in the disorders of metabolic syndrome and is therefore proud to announce that the former LINCOp<sup>lex</sup> Rat Adipocyte Panel now has the MILLIPLEX<sup>®</sup> MAP optimized format. While you will immediately recognize the quality and reproducibility that you have always trusted, you will also enjoy the enhancements that we have built into MILLIPLEX<sup>®</sup> MAP.

EMD Millipore's MILLIPLEX<sup>®</sup> MAP Rat Adipocyte Panel is to be used for the simultaneous quantification of adiponectin, IL-6, leptin, MCP-1, PAI-1 (active or total) and TNF $\alpha$  in rat tissue extract or cell/tissue culture supernatant samples. The panel provides biomedical researchers using rat models with quality tools for the study of obesity related disorders of metabolic disease.

***For Research Use Only. Not for Use in Diagnostic Procedures.***

***Please read entire protocol before use.***

***It is important to use same assay incubation conditions throughout your study.***

## PRINCIPLE

MILLIPLEX<sup>®</sup> MAP is based on the Luminex<sup>®</sup> xMAP<sup>®</sup> technology — one of the fastest growing and most respected multiplex technologies offering applications throughout the life sciences and capable of performing a variety of bioassays including immunoassays on the surface of fluorescent-coded beads known as microspheres.

- Luminex<sup>®</sup> uses proprietary techniques to internally color-code microspheres with two fluorescent dyes. Through precise concentrations of these dyes, 100 distinctly colored bead sets can be created, each of which is coated with a specific capture antibody.
- After an analyte from a test sample is captured by the bead, a biotinylated detection antibody is introduced.
- The reaction mixture is then incubated with Streptavidin-PE conjugate, the reporter molecule, to complete the reaction on the surface of each microsphere.
- The microspheres are allowed to pass rapidly through a laser which excites the internal dyes marking the microsphere set. A second laser excites PE, the fluorescent dye on the reporter molecule.
- Finally, high-speed digital-signal processors identify each individual microsphere and quantify the result of its bioassay based on fluorescent reporter signals.

The capability of adding multiple conjugated beads to each sample results in the ability to obtain multiple results from each sample. Open-architecture xMAP<sup>®</sup> technology enables multiplexing of many types of bioassays reducing time, labor and costs over traditional methods.

## STORAGE CONDITIONS UPON RECEIPT

- Recommended storage for kit components is 2 - 8°C.
- For long-term storage, freeze reconstituted standards and controls at ≤ -20°C. Avoid multiple (>2) freeze/thaw cycles.
- **DO NOT FREEZE Antibody-Immobilized Beads, Detection Antibodies, and Streptavidin-Phycoerythrin.**

## REAGENTS SUPPLIED

**Note: Store all reagents at 2 – 8°C**

REAGENTS SUPPLIED	CATALOG NUMBER	VOLUME	QUANTITY
Rat Adipokine Standard Cocktail	LRA-8081	lyophilized	1 vial
Rat Adipokine Quality Controls 1 and 2	LRA-6081	lyophilized	2 vials
Set of one 96-Well Filter Plate with 2 Sealers	MX-PLATE	-----	1 plate 2 sealers
Assay Buffer	LE-ABGLP	30 mL	1 bottle
10X Wash Buffer Note: Contains 0.05% Proclin	L-WB	30 mL	1 bottle
Rat Adipocyte Detection Antibodies	LRAC-1082	3.2 mL	1 bottle
Streptavidin-Phycoerythrin	L-SAPE	3.2 mL	1 bottle
Mixing Bottle	-----	-----	1 bottle

Included Rat Adipocyte Antibody-Immobilized Beads are dependent on customizable selection of analytes within the panel.

Rat Adipocyte Antibody-Immobilized Beads:

Bead/Analyte Name	Luminex Bead Region	Customizable 8 Analytes (20X concentration, 200 µL)	
		Available	Cat. #
Anti - MCP-1 Bead	13	✓	RMCP-1
Anti - Leptin Bead	16	✓	RME-LPTN
Anti - IL-6 Bead	38	✓	RIL-6
Anti - Adiponectin Bead	51	✓	HA-ADPN
Anti - TNFα Bead	77	✓	RTNF-A
Anti – PAI-1 (Active) Bead*	86	✓	RACTV-PAI1
Anti - PAI-1 (Total) Bead*	89	✓	RA-PAI1

\*Note that Active PAI-1 and Total PAI-1 cannot be run together in the same assay.

## **MATERIALS REQUIRED BUT NOT PROVIDED**

### Reagents

1. Luminex Sheath Fluid (Luminex Catalogue #40-50000)

### Instrumentation / Materials

1. Adjustable Pipettes with Tips capable of delivering 5  $\mu$ L to 1000  $\mu$ L
2. Multichannel Pipettes capable of delivering 5  $\mu$ L to 50  $\mu$ L or 25  $\mu$ L to 200  $\mu$ L
3. Reagent Reservoirs
4. Polypropylene Microfuge Tubes
5. Rubber Bands
6. Absorbent Pads
7. Laboratory Vortex Mixer
8. Sonicator (Branson Ultrasonic Cleaner Model #B200 or equivalent)
9. Titer Plate Shaker (Lab-Line Instruments Model #4625 or equivalent)
10. Vacuum Filtration Unit (EMD Millipore Vacuum Manifold Catalog #MSVMHTS00 or equivalent with EMD Millipore Vacuum Pump Catalog #WP6111560 or equivalent)
11. Luminex 100™ IS, 200™, or HTS by Luminex Corporation
12. Plate Stand (EMD Millipore Catalog # MX-STAND)

## **SAFETY PRECAUTIONS**

- All blood components and biological materials should be handled as potentially hazardous. Follow universal precautions as established by the Centers for Disease Control and Prevention and by the Occupational Safety and Health Administration when handling and disposing of infectious agents.
- Sodium Azide or Proclin has been added to some reagents as a preservative. Although the concentrations are low, Sodium Azide and Proclin may react with lead and copper plumbing to form highly explosive metal azides. Dispose of unused contents and waste in accordance with international, federal, state, and local regulations.

## **TECHNICAL GUIDELINES**

To obtain reliable and reproducible results, the operator should carefully read this entire manual and fully understand all aspects of each assay step before running the assay. The following notes should be reviewed and understood before the assay is set up.

- FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- Do not use beyond the expiration date on the label.
- Do not mix or substitute reagents with those from other lots or sources.
- The Antibody-Immobilized Beads are light sensitive and must be protected from light at all times. Cover the assay plate containing beads with opaque plate lid or aluminum foil during all incubation steps.
- It is important to allow all reagents to warm to room temperature (20-25°C) before use in the assay.

## TECHNICAL GUIDELINES (continued)

- The bottom of the Microtiter Filter Plate should not be in direct contact with any surface during assay set-up or incubation times. The plate can be set on a plate stand or on the non-flat side of the plate cover or any other plate holder to raise the plate from the surface. A plate stand can be purchased separately from EMD Millipore (EMD Millipore Catalog #MX-STAND).
- Incomplete washing can adversely affect the assay outcome. All washing must be performed with the Wash Buffer provided.
- After the wash steps, keep the bottom of the Microtiter Filter Plate clean by blotting on paper towels or absorbent pads to prevent any leakage due to capillary action.
- Keep the vacuum suction on the plate as low as possible. It is recommended to have a vacuum setting that will remove 200  $\mu$ L of buffer in  $\geq 5$  seconds (equivalent to  $< 100$  mmHg).
- The Standards prepared by serial dilution must be used within 1 hour of preparation. Discard any unused standards except the standard stock which may be stored at  $\leq -20^{\circ}\text{C}$  for 1 month and at  $\leq -80^{\circ}\text{C}$  for greater than one month.
- If samples fall outside the dynamic range of the assay, further dilute the samples with the appropriate diluent and repeat the assay.
- Any unused mixed Antibody-Immobilized Beads may be stored in the Mixing Bottle at  $2-8^{\circ}\text{C}$  for up to one month.
- During the preparation of the standard curve, make certain to mix the higher concentration well before making the next dilution. Use a new tip with each dilution.
- The plate should be read immediately after the assay is finished. If, however, the plate cannot be read immediately, seal the plate, cover with aluminum foil or an opaque lid, and store the plate at  $2-8^{\circ}\text{C}$  for up to 24 hours. Prior to reading, agitate the plate on the plate shaker at room temperature for 10 minutes. Delay in reading a plate may result in decreased sensitivity for some adipocytes.
- The titer plate shaker should be set at a speed to provide maximum orbital mixing without splashing of liquid outside the wells. For the recommended plate shaker, this would be a setting of 5-7 which is approximately 500-800 rpm.
- Ensure that the needle probe is clean. This may be achieved by sonication and/or alcohol flushes. Adjust probe height according to the protocols recommended by Luminex to the kit filter plate using 3 alignment discs prior to reading an assay.
- For cell culture supernatants or tissue extraction, use the culture or extraction medium as the matrix solution in background, standard curve and control wells. If samples are diluted in Assay Buffer, use the Assay Buffer as matrix.
- For cell/tissue homogenate, the final cell or tissue homogenate should be prepared in a buffer that has a neutral pH, contains minimal detergents or strong denaturing detergents, and has an ionic strength close to physiological concentration. Avoid debris, lipids, and cell/tissue aggregates. Centrifuge samples before use.
- Vortex all reagents well before adding to plate.

## SAMPLE COLLECTION AND STORAGE

### Preparation of Tissue/Cell Extract or Supernatant:

- Centrifuge the sample to remove debris and assay immediately or aliquot and store samples at  $\leq -20^{\circ}\text{C}$ .
- Avoid multiple ( $>2$ ) freeze/thaw cycles.
- Appropriate sample dilution with adipokine-free culture media may be required.

### NOTE:

- A maximum of 25  $\mu\text{L}$  per well tissue/cell extract or supernatant samples collected from adipose/adipocytes can be used.
- All samples must be stored in polypropylene tubes. **DO NOT STORE SAMPLES IN GLASS.**

## PREPARATION OF REAGENTS FOR IMMUNOASSAY

### A. Preparation of Antibody-Immobilized Beads

**Prepare the bead mixture right before setting up the assay. Early preparation of the bead mixture ( $>1$  day) should be avoided.** For individual vials of beads, sonicate each antibody-bead vial for 30 seconds; vortex for 1 minute. Add 150  $\mu\text{L}$  from each antibody bead vial to the Mixing Bottle and bring final volume to 3.0 mL with Assay Buffer. Vortex the mixed beads well. **Unused portions of the bead mixture should be discarded after use.** Left-over individual beads can be stored at  $2-8^{\circ}\text{C}$  for up to one month.

Example 1: When using 4 adipocyte antibody-immobilized beads, add 150  $\mu\text{L}$  from each of the 4 bead sets to the Mixing Bottle. Then add 2.40 mL Assay Buffer.

Example 2: When using 8 adipocyte antibody-immobilized beads, add 150  $\mu\text{L}$  from each of the 8 bead sets to the Mixing Bottle. Then add 1.80 mL Assay Buffer.

**Note that Active PAI-1 and Total PAI-1 cannot be run together in the same assay.**

### B. Preparation of Quality Controls

Before use, reconstitute Quality Control 1 and Quality Control 2 with 250  $\mu\text{L}$  deionized water. Invert the vial several times to mix and vortex. Allow the vial to sit for 5-10 minutes and then vortex again. Unused portion may be stored at  $\leq -20^{\circ}\text{C}$  for up to one month.

### C. Preparation of Wash Buffer

Bring the 10X Wash Buffer to room temperature and mix to bring all salts into solution. Dilute 30 mL of 10X Wash Buffer with 270 mL deionized water. Store unused portion at  $2-8^{\circ}\text{C}$  for up to one month.

## PREPARATION OF REAGENTS FOR IMMUNOASSAY (continued)

### D. Preparation of Rat Adipocyte Standard

1.) Prior to use, reconstitute the Rat Adipocyte Standard with 250  $\mu$ L deionized water to give the highest stock concentration of standard. Invert the vial several times to mix. Vortex the vial for 10 seconds. Allow the vial to sit for 5-10 minutes and then vortex again. The unused portions of this standard may be stored at  $\leq -20^{\circ}\text{C}$  for up to one month.

### 2.) Preparation of Working Standards


Label six polypropylene microfuge tubes Standard 1, Standard 2, Standard 3, Standard 4, Standard 5 and Standard 6. Add 150  $\mu$ L of Assay Buffer to each of the six tubes. Prepare 4 time serial dilutions by adding 50  $\mu$ L of the Standard 7 to the Standard 6 tube, mix well and transfer 50  $\mu$ L of the Standard 6 tube to the Standard 5 tube, mix well and transfer 50  $\mu$ L of the Standard 5 tube to the Standard 4 tube, mix well and transfer 50  $\mu$ L of the Standard 4 tube to the Standard 3 tube, mix well and transfer 50  $\mu$ L of the Standard 3 tube to the Standard 2 tube, mix well and transfer 50  $\mu$ L of the Standard 2 tube to the Standard 1 tube and mix well. The 0 standard (Background) will be the Assay Buffer.

Standard Dilution	Volume of Deionized Water to Add	Volume of Standard to Add
Original Standard 7	250 $\mu$ L water	0

Standard Dilution	Volume of Assay Buffer to Add	Volume of Standard to Add
Standard 6	150 $\mu$ L	50 $\mu$ L of Standard 7
Standard 5	150 $\mu$ L	50 $\mu$ L of Standard 6
Standard 4	150 $\mu$ L	50 $\mu$ L of Standard 5
Standard 3	150 $\mu$ L	50 $\mu$ L of Standard 4
Standard 2	150 $\mu$ L	50 $\mu$ L of Standard 3
Standard 1 (optional)	150 $\mu$ L	50 $\mu$ L of Standard 2

## PREPARATION OF REAGENTS FOR IMMUNOASSAY (continued)

### Standards Preparation

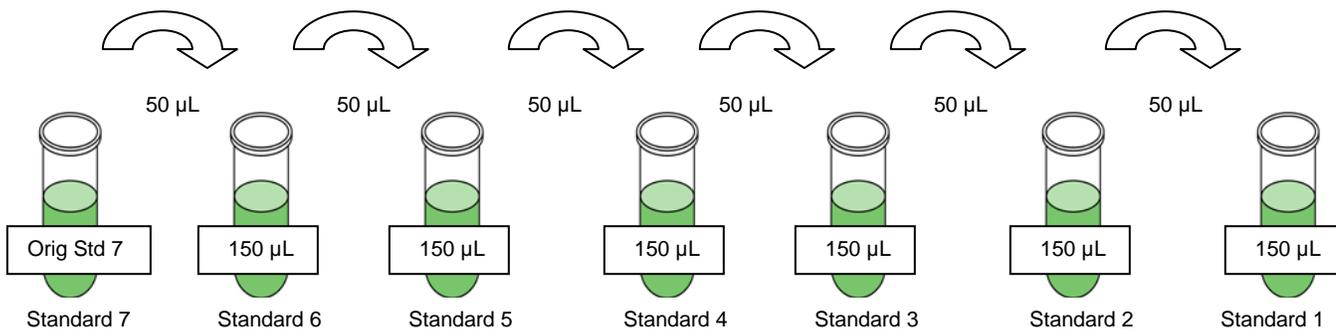

The serial dilutions result in the following concentrations of standards.

Standard Dilution	TNF- $\alpha$ , MCP-1 Standard (pg/ml)	Leptin, IL-6, PAI-1 Standard (pg/ml)	Adiponectin (pg/ml)
Standard 7	20,000	50,000	500,000
Standard 6	5,000	12,500	125,000
Standard 5	1,250	3,125	31,250
Standard 4	312.5	781.25	7,813
Standard 3	78.125	195.31	1,953
Standard 2	19.53	48.83	488
Standard 1	4.88	12.21	122

## IMMUNOASSAY PROCEDURE

- Prior to beginning this assay, it is imperative to read this protocol completely and to thoroughly understand the Technical Guidelines.
  - Allow all reagents to warm to room temperature (20-25°C) before use in the assay.
  - Diagram the placement of Standards [0 (Background), 1, 2, 3, 4, 5, 6 and 7], Controls 1 and 2, and Samples on Well Map Worksheet in a vertical configuration. (Note: Most instruments will only read the 96-well plate vertically by default.) It is recommended to run the assay in duplicate.
  - Set the filter plate on a plate holder at all times during reagent dispensing and incubation steps so that the bottom of the plate does not touch any surface.
1. Prewet the filter plate by pipetting 200  $\mu$ L of Assay Buffer into each well of the Microtiter Filter Plate. Seal and mix on a plate shaker for 10 minutes at room temperature (20-25°C).
  2. Remove Assay Buffer by vacuum. (**NOTE: DO NOT INVERT PLATE.**) Blot excess Assay Buffer from the bottom of the plate with an absorbent pad or paper towels.
  3. Add 25  $\mu$ L of each Standard or Control into the appropriate wells. Assay Buffer should be used for the 0 pg/mL standard (Background).
  4. Add 25  $\mu$ L of Assay Buffer to the sample wells.
  5. Add 25  $\mu$ L of adipokine free culture media to the background, standards, and control wells.
  6. Add 25  $\mu$ L of Sample into the appropriate wells.
  7. Vortex Mixing Bottle and add 25  $\mu$ L of the mixed beads to each well. (Note: During addition of beads, shake bead bottle intermittently to avoid settling.)
  8. Seal the plate with a plate sealer, cover it with the lid. Wrap a rubber band around the plate holder, plate and lid and incubate with agitation on a plate shaker overnight (16-18hours) at 4°C
  9. Gently remove fluid by vacuum. (**NOTE: DO NOT INVERT PLATE.**)
  10. Wash plate 3 times with 200  $\mu$ L/well of Wash Buffer, removing Wash Buffer by vacuum filtration between each wash. Blot excess Wash Buffer from the bottom the plate by with an absorbent pad or paper towels.

Add 200  $\mu$ L Assay Buffer per well



Shake 10 min, RT

Vacuum

- Add 25  $\mu$ L Standard or Control to appropriate wells
- Add 25  $\mu$ L Assay Buffer to background and sample wells
- Add 25  $\mu$ L adipokine free culture media to background, standards and control wells
- Add 25  $\mu$ L Samples to sample wells
- Add 25  $\mu$ L Beads to each well



Incubate overnight at 4°C with shaking

Vacuum and wash 3X with 200  $\mu$ L Wash Buffer

11. Add 50  $\mu$ L of Detection Antibodies into each well. (Note: Allow the Detection Antibodies to warm to room temperature prior to addition.)
12. Seal, cover with lid, and incubate with agitation on a plate shaker for 1 hour at room temperature (20-25°C) **DO NOT VACUUM AFTER INCUBATION.**
13. Add 50  $\mu$ L Streptavidin-Phycoerythrin to each well containing the 50  $\mu$ L of Detection Antibodies.
14. Seal, cover with lid and incubate with agitation on a plate shaker for 30 minutes at room temperature (20-25°C).
15. Gently remove all contents by vacuum. **(NOTE: DO NOT INVERT PLATE.)**
16. Wash plate 3 times with 200  $\mu$ L/well Wash Buffer, removing Wash Buffer by vacuum filtration between each wash. Wipe any excess buffer on the bottom of the plate with a tissue.
17. Add 100  $\mu$ L of Sheath Fluid to all wells. Resuspend the beads on a plate shaker for 5 minutes.
18. Run plate on Luminex 100™ IS, 200™, or HTS.
19. Save and analyze the Median Fluorescent Intensity (MFI) data using a weighted 5-parameter logistic or spline curve-fitting method for calculating adipocyte concentrations in samples.

Add 50  $\mu$ L Detection Antibody per well



Incubate 1 hour at RT

Do Not Vacuum

Add 50  $\mu$ L Streptavidin-Phycoerythrin per well



Incubate for 30 minutes at RT

Vacuum and wash 3X with 200  $\mu$ L Wash Buffer

Add 100  $\mu$ L Sheath Fluid per well

Read on Luminex (50  $\mu$ L, 50 beads per bead set)

## EQUIPMENT SETTINGS

These specifications are for the Luminex 100™ IS v.1.7 or Luminex 100™ IS v2.1/2.2, Luminex 200™ v2.3, xPONENT, and Luminex HTS. Luminex instruments with other software (e.g. MasterPlex, StarStation, LiquiChip, Bio-Plex, LABScan100) would need to follow instrument instructions for gate settings and additional specifications from the vendors.

Events:	50, per bead	
Sample Size:	50 µL	
Gate Settings	8,000 to 15,000	
Reporter Gain:	Default (Low PMT)	
Time Out:	60 seconds	
Bead Set:	8-Plex Beads	
	MCP-1	13
	Leptin	16
	IL-6	38
	Adiponectin	51
	TNFα	77
	PAI-1 active	86
	PAI-1 total	89

## QUALITY CONTROLS

The ranges for each analyte in Quality Control 1 and 2 are provided on the card insert or can be located at the EMD MILLIPORE website

[www.millipore.com/techlibrary/index.do](http://www.millipore.com/techlibrary/index.do) using the catalog number as the keyword.

## ASSAY CHARACTERISTICS

### Cross-Reactivity

There was no or negligible cross-reactivity between the antibodies and any of the other analytes in this panel.

### Assay Sensitivities (minimum detectable concentrations, pg/mL)

MinDC: Minimum Detectable Concentration is calculated by the StatLIA® Immunoassay Analysis Software from Brendan Technologies. It measures the true limits of detection for an assay by mathematically determining what the empirical MinDC would be if an infinite number of standard concentrations were run for the assay under the same conditions.

MCP-1	10.5
Leptin	2.5
IL-6	3.4
Adiponectin	32.7
TNF $\alpha$	1.1
PAI-1 (active)	14.7
PAI-1 (total)	12.5

### Precision

Intra-assay precision is generated from the mean of the %CV's from 8 reportable results across two different concentration of adipocytes in one experiment. Inter-assay precision is generated from the mean of the %CV's from two reportable results each for two different concentrations of adipocytes across 4 different experiments.

Intra-assay variation (%CV)	<8%
Inter-assay variation (%CV)	<20%

### Accuracy

Spike Recovery: The data represent mean percent recovery of 7 levels of spiked standards in serum matrices in 8 independent experiments.

MCP-1	99%
Leptin	97%
IL-6	89%
Adiponectin	99%
TNF- $\alpha$	99%
PAI-1 (active)	92%
PAI-1 (total)	117%

## TROUBLESHOOTING GUIDE

Problem	Probable Cause	Solution
Filter plate will not vacuum	Vacuum pressure is insufficient	Increase vacuum pressure such that 0.2mL buffer can be suctioned in 3-5 seconds.
	Samples have insoluble particles	Centrifuge samples just prior to assay set-up and use supernatant.  If high lipid concentration, after centrifugation, remove lipid layer and use supernatant.
	Sample too viscous	May need to dilute sample.
Insufficient bead count	Vacuum pressure too high	Adjust vacuum pressure such that 0.2mL buffer can be suctioned in 3-5 seconds.
	Bead mix prepared incorrectly	Sonicate bead vials and vortex just prior to adding to bead mix bottle according to protocol. Agitate bead mix intermittently in reservoir while pipetting into the plate.
	Samples cause interference due to particulate matter or viscosity	See above. Also sample probe may need to be cleaned with alcohol flush, backflush and washes; or, if needed, probe should be removed and sonicated.
	Probe height not adjusted correctly	Adjust probe to 3 alignment discs in well H6.
Plate leaked	Vacuum pressure too high	Adjust vacuum pressure such that 0.2mL buffer can be suctioned in 3-5 seconds. May need to transfer contents to a new (prewetted) plate and continue.
	Plate set directly on table or absorbent towels during incubations or reagent additions	Set plate on plate stand or raised edge so bottom of filter is not touching any surface.
	Insufficient blotting of filter plate bottom causing wicking	Blot the bottom of the filter plate well with absorbent towels after each wash step.
	Pipette touching plate filter during additions	Pipette to the side of well.
	Probe height not adjusted correctly	Adjust probe to 3 alignment discs in well H6.
Background is too high	Background wells were contaminated	Avoid cross-well contamination by using sealer appropriately and by pipeting with multichannel pipets without touching reagent in plate.
	Matrix used has endogenous analyte or interference	Check matrix ingredients for crossreacting components (e.g. interleukin modified tissue culture medium).
	Insufficient washes	Increase number of washes.

<p>Beads not in region or gate</p>	<p>Luminex not calibrated correctly or recently</p> <p>Gate settings not adjusted correctly</p> <p>Wrong bead regions in protocol template</p> <p>Incorrect sample type used</p> <p>Instrument not washed or primed</p> <p>Beads were exposed to light</p>	<p>Calibrate Luminex based on instrument manufacturer's instructions at least once a week or if temperature has changed by &gt;3°C.</p> <p>Some Luminex instruments (e.g. Bio-Plex) require different gate settings than those described in the kit protocol. Use instrument default settings.</p> <p>Check kit protocol for correct bead regions or analyte selection.</p> <p>Samples containing organic solvents or if highly viscous should be diluted or dialyzed as required.</p> <p>Prime the Luminex 4 times to eliminate air bubbles. Wash 4 times with sheath fluid or water if there is any remnant alcohol or sanitizing liquid.</p> <p>Keep plate and bead mix covered with dark lid or aluminum foil during all incubation steps.</p>
<p>Signal for whole plate is same as background</p>	<p>Incorrect or no Detection Antibody was added</p> <p>Streptavidin-Phycoerythrin was not added</p>	<p>Add appropriate Detection Antibody and continue.</p> <p>Add Streptavidin-Phycoerythrin according to protocol. If Detection Antibody has already been vacuumed out, sensitivity may be low.</p>
<p>Low signal for standard curve</p>	<p>Detection Antibody may have been vacuumed out prior to adding Streptavidin Phycoerythrin</p> <p>Incubations done at incorrect temperatures, timings or agitation</p>	<p>May need to repeat assay if desired sensitivity not achieved.</p> <p>Assay conditions need to be checked.</p>
<p>Signals too high, standard curves are saturated</p>	<p>Calibration target value set too high</p> <p>Plate incubation was too long with standard curve and samples</p>	<p>With some Luminex instruments (e.g. Bio-Plex) default target setting for RP1 calibrator is set at High PMT. Use low target value for calibration and reanalyze plate.</p> <p>Use shorter incubation time.</p>
<p>Sample readings are out of range</p>	<p>Samples contain no or below detectable levels of analyte</p> <p>Samples contain analyte concentrations higher than highest standard point</p> <p>Standard curve was saturated at higher end of curve</p>	<p>If below detectable levels, it may be possible to use higher sample volume. Check with tech support for appropriate protocol modifications.</p> <p>Samples may require dilution and reanalysis for that particular analyte.</p> <p>See above.</p>

High variation in samples and/or standards	Multichannel pipet may not be calibrated	Calibrate pipets.
	Plate washing was not uniform	Confirm all reagents are vacuumed out completely in all wash steps.
	Samples may have high particulate matter or other interfering substances	See above.
	Plate agitation was insufficient	Plate should be agitated during all incubation steps using a vertical plate shaker at a speed where beads are in constant motion without splashing.
	Cross-well contamination	Check when reusing plate sealer that no reagent has touched sealer.  Care should be taken when using same pipet tips that are used for reagent additions and that pipet tip does not touch reagent in plate.

## REPLACEMENT REAGENTS

## Catalog #

Rat Adipokine Standard  
Rat Adipokine Quality Controls 1 & 2  
Rat Adipocyte Detection Antibodies  
Streptavidin-Phycoerythrin  
Assay Buffer  
Set of two 96-Well Filter Plates with Sealers  
10X Wash Buffer  
Antibody-Immobilized Beads

LRA-8081  
LRA-6081  
LRAC-1082  
L-SAPE  
LE-ABGLP  
MX-PLATE  
L-WB

<u>Analyte</u>	<u>Bead #</u>	<u>Cat. #</u>
MCP-1	13	RMCP-1
Leptin	16	RME-LPTN
IL-6	38	RIL-6
Adiponectin	51	HA-ADPN
TNF $\alpha$	77	RTNF-A
PAI-1 (active)	86	RACTV-PAI1
PAI-1 (total)	89	RA-PAI1

## ORDERING INFORMATION

### To place an order:

To assure the clarity of your custom kit order, please FAX the following information to our customer service department:

Include:

- Your name, telephone and/or fax number
- Customer account number
- Shipping and billing address
- Purchase order number
- Catalog number and description of product
- Quantity of kits
- Selection of MILLIPLEX<sup>®</sup> Analytes

FAX: (636) 441-8050

Toll-Free US: (866) 441-8400  
(636) 441-8400

Mail Orders: EMD Millipore Corporation  
6 Research Park Drive  
St. Charles, Missouri 63304 U.S.A.

### For International Customers:

To best serve our international customers in placing an order or obtaining additional information about MILLIPLEX<sup>®</sup> MAP products, please contact your multiplex specialist or sales representative or email our European Customer Service at [customerserviceEU@merckgroup.com](mailto:customerserviceEU@merckgroup.com).

### Conditions of Sale

For Research Use Only. Not for Use in Diagnostic Procedures.

### Material Safety Data Sheets (MSDS)

Material Safety Data Sheets for EMD Millipore products may be ordered by fax or phone or through our website at [www.emdmillipore.com/techlibrary/index.do](http://www.emdmillipore.com/techlibrary/index.do).

### Technical Services

<http://www.emdmillipore.com/techservices>

To contact by phone

For North America: Toll-Free US: 1-(800) 221-1975 or 1-(781) 533-8045

Outside North America, contact your local office

<http://www.emdmillipore.com/offices>

### WELL MAP

	1	2	3	4	5	6	7	8	9	10	11	12
A	0 Standard (Background)	Standard 4	QC-1 Control	Etc.								
B	0 Standard (Background)	Standard 4	QC-1 Control									
C	Standard 1	Standard 5	QC-2 Control									
D	Standard 1	Standard 5	QC-2 Control									
E	Standard 2	Standard 6	Sample 1									
F	Standard 2	Standard 6	Sample 1									
G	Standard 3	Standard 7	Sample 2									
H	Standard 3	Standard 7	Sample 2									