

## Certified Reference Material Reference material certificate

### Elemental Impurities according to ICH Q3D oral, Standard 2

**TraceCERT®**  
Traceable Certified Reference Materials

**Product no.:** 73108  
**Lot no.:** BCCM5999  
**Description of CRM:** High-purity starting materials in 10% HCl (prepared with HCl suitable for trace analysis and high-purity water, 18.2 MΩ·cm, 0.22 µm filtered).  
**Expiry date:** NOV 2028 (unopened bottle in aluminized bag)  
**Storage:** Store at 5°C-25°C  
**Density at 20°C:** 1048.5 kg m<sup>-3</sup> ± 0.5 kg m<sup>-3</sup> (k = 2)

| Constituent    | Certified values at 20°C and expanded uncertainties, $U = k \cdot u$ ( $k = 2$ ) <sup>[1][2]</sup> |   |  |  |
|----------------|--|---|--|--|
| Gold (Au)      | 95.4 mg kg <sup>-1</sup> ± 0.5 mg kg <sup>-1</sup>   | 100.0 mg L <sup>-1</sup> ± 0.6 mg L <sup>-1</sup> |  |  |
| Iridium (Ir)   | 95.4 mg kg <sup>-1</sup> ± 1.0 mg kg <sup>-1</sup>   | 100.0 mg L <sup>-1</sup> ± 1.0 mg L <sup>-1</sup> |  |  |
| Osmium (Os)    | 95.4 mg kg <sup>-1</sup> ± 0.7 mg kg <sup>-1</sup>   | 100.0 mg L <sup>-1</sup> ± 0.8 mg L <sup>-1</sup> |  |  |
| Palladium (Pd) | 95.4 mg kg <sup>-1</sup> ± 0.3 mg kg <sup>-1</sup>   | 100.0 mg L <sup>-1</sup> ± 0.4 mg L <sup>-1</sup> |  |  |
| Platinum (Pt)  | 95.4 mg kg <sup>-1</sup> ± 0.5 mg kg <sup>-1</sup>   | 100.0 mg L <sup>-1</sup> ± 0.6 mg L <sup>-1</sup> |  |  |
| Rhodium (Rh)   | 95.4 mg kg <sup>-1</sup> ± 0.5 mg kg <sup>-1</sup>   | 100.0 mg L <sup>-1</sup> ± 0.5 mg L <sup>-1</sup> |  |  |
| Ruthenium (Ru) | 95.4 mg kg <sup>-1</sup> ± 0.9 mg kg <sup>-1</sup>   | 100.0 mg L <sup>-1</sup> ± 1.0 mg L <sup>-1</sup> |  |  |

**Metrological traceability:** Certified values are traceable to the International System of units (SI) through a metrologically valid weighing process. Details see "Details on metrological traceability".<sup>[3]</sup>

**Measurement method:** The certified value is determined by high-precision weighing of thoroughly characterized starting materials and verified by measurement against NIST SRMs or similar CRMs in accordance with ISO/IEC 17025.<sup>[4]</sup>

**Intended use:** This reference material is intended for use as a calibration standard in elemental analysis.

**Instructions for handling and correct use:** The bottle's temperature must be 20°C. Shake well before every use. If storage of a partially used bottle is necessary (at the user's risk), the cap should be tightly sealed and the bottle should be stored at reduced temperature (e.g. refrigerator) to minimize transpiration rate.

**Health and safety information:** Please refer to the Safety Data Sheet for detailed information about the nature of any hazard and appropriate precautions to be taken.

**Packaging:** 100 mL FEP bottle, sealed with an aluminized bag

**Accreditation:** Sigma-Aldrich Production GmbH is accredited by the Swiss Accreditation Service SAS as reference material producer under no. SRMS 0001 in accordance with international standard ISO 17034<sup>[5]</sup>

**Certificate issue date:** 17 APR 2025



ISO 17034  
SRMS 0001

Dr. P. Zell – Approving Officer  
(Quality Assurance)



### Certification process details:

To guarantee top reliability of the values for this *TraceCERT*<sup>®</sup> certified reference material, two independent procedures were followed. The values have to agree in the range of their uncertainties, but the value from the gravimetric preparation has been chosen as certified value [3]:

1. Gravimetric preparation using pure materials is a practical realization of concentration units, through conversion of mass to amount of substance [3]. If the purity of the materials is demonstrated and if contamination and loss of material is strictly prevented this approach allows highest accuracy and small uncertainties. The certified value of this *TraceCERT*<sup>®</sup> reference material is based on this approach and directly traceable to the SI unit kilogram. Therefore comprehensively characterized materials of high purity are used. All balances are calibrated annually by an ISO/IEC 17025 accredited laboratory and certified according to DKD guidelines. Calibration is checked daily with OIML Class E2 or F2 weights.
2. The starting material is measured against a certified reference material (i.e. NIST or BAM) followed by gravimetric preparation using balances calibrated with SI-traceable weights. Consequently the value calculated by this unbroken chain of comparisons is traceable to the reference to which the starting material is compared.

### Details on metrological traceability:

Only internationally accepted reference materials e.g. from NIST (USA) or BAM (Germany) have been carefully selected to provide the basis for traceability to the SI unit mole. When no such reference is available, an elemental metal or an adequate salt of highest available purity is used to confirm traceability to this pure material (and therefore to the SI unit kg).

To underpin the certified gravimetric value all traceability measurements are performed with the most accurate and precise analytical technique available. Therefore titrimetry measurement series are applied whenever possible (corrected for trace impurities). When no titrimetric technique is available, the traceability measurements are performed with another analytical technique, e.g. ICP-OES or AAS.

Reference and applied technique used for traceability measurements of the starting material:

| Constituent    | Starting material                                 | Reference                      | Method  |
|----------------|---|--------------------------------|---------|
| Gold (Au)      | H(AuCl <sub>4</sub> )                             | NIST SRM 3121                  | ICP-OES |
| Iridium (Ir)   | IrCl <sub>3</sub> x 3 H <sub>2</sub> O            | Traceable to in-house standard | ICP-OES |
| Osmium (Os)    | (NH <sub>4</sub> ) <sub>2</sub> OsCl <sub>6</sub> | Traceable to in-house standard | ICP-OES |
| Palladium (Pd) | Pd metal  | NIST SRM 3138                  | ICP-OES |
| Platinum (Pt)  | H <sub>2</sub> PtCl <sub>6</sub>                  | NIST SRM 3140                  | ICP-OES |
| Rhodium (Rh)   | RhCl <sub>3</sub>                                 | NIST SRM 3144                  | ICP-OES |
| Ruthenium (Ru) | RuCl <sub>3</sub>                                 | Traceable to in-house standard | ICP-OES |

### Details on starting materials:

For high purity materials ( $P > 99.9\%$ ) the most appropriate way of purity determination is to quantify the impurities ( $w_i$ ) and to subtract the sum from 100%. Impurities below the detection limit are considered with a contribution of half of the detection limit ( $DL_j$ ).

$$P = 100\% - \sum_i w_i - \sum_j \left( \frac{DL_j}{2} \right)$$

Water containing materials were dried to absolute dryness by individual drying conditions (up to 600°C). When drying is impossible due to decomposition water was determined by high-precision KF-titration.

### Homogeneity assessment:

Due to the production process, a homogeneous solution derives. Nevertheless a small homogeneity contribution is included into the calculation of content uncertainty of this CRM.

### Stability assessment:

The storage behavior of standard solutions is an important factor regarding the certified value. An ideal container for a CRM solution is chemically inert, completely tight and does neither adsorb the analyte nor leach trace impurities into the solution. Therefore the most important storing effects were investigated by in-depth studies.

Every plastic bottle (e.g. PP, HDPE, FEP, PFA) is untight for a certain range due to the transpiration of solvent through the wall which leads to an increase of the analyte concentration over time. To avoid significant loss of solvent through transpiration the bottle is delivered in aluminum coated bags. When the CRM is in the sealed Al-bag the change in concentration at 23°C is less than 0.02% per year. After the bottle has been removed from the bag, transpiration will occur at an accelerated rate depending on the temperature. We highly recommend not to open the bag until the solution is needed. Once the bottle is opened the solution should be stored at reduced temperature (4°C) to minimize the transpiration rate.

### Density Measurement:


The density measurement is carried out in accordance with ISO/IEC 17025<sup>[4]</sup> and ISO 15212-1<sup>[6]</sup> using the digital density meter DMA 4500M from Anton Paar with an oscillating U-tube installed. The measurement uncertainty is calculated according to Eurachem/CITAC Guide and reported as combined expanded uncertainty at the 95% confidence level, using a coverage factor of  $k = 2$ .

### Uncertainty evaluation:

The uncertainty contributions are illustrated by the following cause-effect diagram<sup>[7]</sup>:

Typical relative contributions are:

|                |          |
|----------------|----------|
| $u(m_{SM})$    | < 0.01 % |
| $u(m_{Batch})$ | < 0.01 % |
| $u(P_{SM})$    | < 0.05 % |
| $u_{hom}$      | < 0.03 % |
| $u_{stab}$     | < 0.09 % |
| $u(\rho)$      | < 0.05 % |

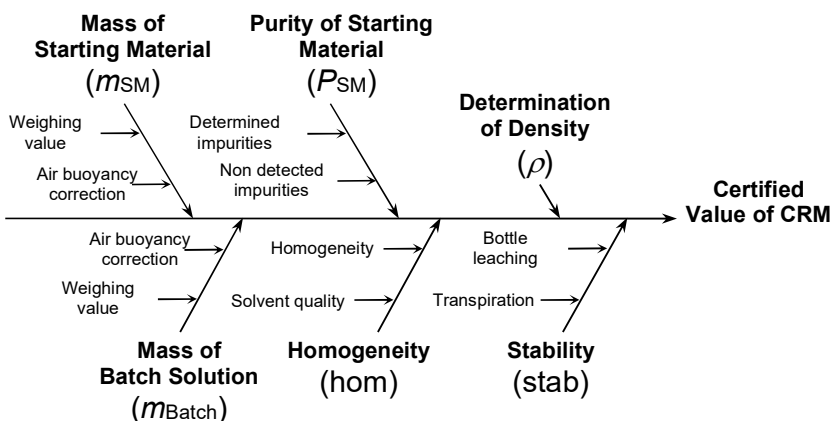

The combined standard uncertainty is calculated by combination of the standard uncertainties of the input estimates according to Eurachem/CITAC Guide "Quantifying Uncertainty in Analytical Measurement" and ISO 17034.<sup>[2][5]</sup>

Expanded uncertainty is then calculated to a confidence level of 95%, typically by multiplying with a coverage factor of  $k=2$ .

### References:

- [1] ISO Guide 35:2017, "Reference materials - Guidance for characterization and assessment of homogeneity and stability"
- [2] Eurachem/CITAC Guide, 3<sup>rd</sup> Ed. (2012), "Quantifying uncertainty in analytical measurement"
- [3] Eurachem/CITAC Guide, 2<sup>nd</sup> Ed. (2019), "Metrological Traceability in chemical measurement"
- [4] The accredited testing laboratory STS 0490 performs the measurements and weighing steps for the certification of this CRM under ISO/IEC 17025:2017, "General requirements for the competence of testing and calibration laboratories"
- [5] ISO 17034:2016, "General requirements for the competence of reference material producers"
- [6] DIN EN ISO 15212-1:1998, Oscillation-type density meters - Part 1: Laboratory instruments
- [7] Reichmuth, A., Wunderli, S., Weber, M., Meyer, V. R. (2004), "The uncertainty of weighing data obtained with electronic analytical balances", *Microchimica Acta* 148: 133-141.

**Certificate of analysis revision history:**

| Certificate version | Certificate issue date | Reason for version |
|---------------------|------------------------|--------------------|
| 01                  | 17 APR 2025            | Initial version    |

The most recent version of the Certificate is available online ([www.sigmaaldrich.com](http://www.sigmaaldrich.com)).

**Disclaimer:**

The purchaser must determine the suitability of this product for its particular use. Sigma-Aldrich Production GmbH makes no warranty of any kind, express or implied, other than its products meet all quality control standards set by Sigma-Aldrich Production GmbH. We do not guarantee that the product can be used for a special application.

The vibrant M, Supelco, TraceCERT and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources. © 2018 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.



The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the US and Canada.