

# UVRAG Hu-Cy5 SmartFlare™ RNA Detection Probe

Cat. # SF-753

pack size: 50µL (250 rxns)

Store at 2-8°C,  
after reconstitution store at 23-27°C  
DO NOT FREEZE

FOR RESEARCH USE ONLY  
NOT FOR USE IN DIAGNOSTIC PROCEDURES  
NOT FOR HUMAN OR ANIMAL CONSUMPTION



## Product Data Sheet

page 1 of 1

Validated Accession NM\_003369.3


Species Hu

Gene Aliases p63, DHTX, VPS38

### Confirmation of UVRAG SmartFlare Performance:

UVRAG SmartFlare has been tested in a buffer system to detect the release of the fluorophore in the presence of a complementary base pair sequence for each lot to confirm target specificity.

UVRAG SmartFlare has also been tested in a cell model system and demonstrated increased fluorescence in cells expressing the target compared to a scrambled negative control SmartFlare (Figure 1). For additional accession numbers predicted to react with this SmartFlare Probe, please visit <http://www.millipore.com/catalogue/item/SF-753>

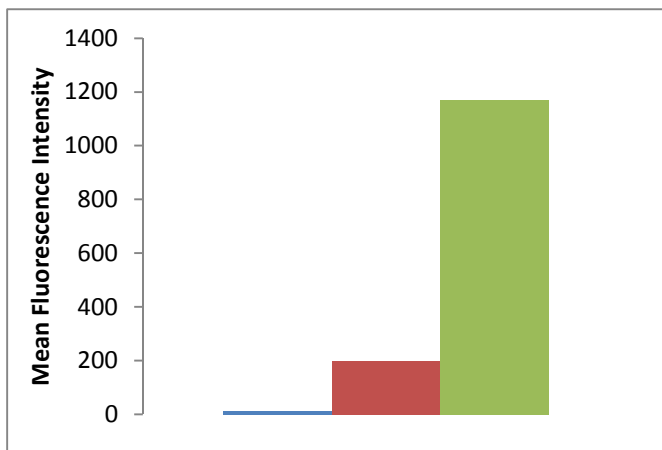

**Figure 1:** UVRAG Mean Fluorescence Intensity (green) measured by flow cytometry in living A431 cells demonstrated a significant increase over unflared cells (blue) as well as scramble control (red). Data shown in graph is representative.

### Storage and Handling:

Material has been 0.22µm filtered. Stable for 5 years at 2-8°C degrees in lyophilized format ONLY. Room temperature is required for reconstituted product.

**Warning-after reconstitution product is sensitive to cold and hot temperatures, a stable room temperature of 23-27°C is required.**

### Handling Recommendations:

Reconstitute with sterile nuclease free water in a drop wise fashion and tap tube repeatedly to fully dissolve lyophilized material. Vortex for 5-10 sec.

Upon reconstitution, store at room temperature for up to 1 year protected from light. Product must be handled with gloves as product can be absorbed through the skin.

### Recommended Cell Testing Protocol:

(example: 30,000 cells in a 200µL media volume within each well of a 96 well plate)

- Reconstitute reagent in 50µL of sterile nuclease free water.
- Create a working solution based on your experiment by diluting 1:20 in sterile PBS.
- Add 4µL directly to cells (at approx 80% confluency).
- Allow to incubate overnight for 16 hrs.
- Detect using fluorescence detection platform of choice.

**SPECIES LEGEND:** Hu Human Ms Mouse Rt Rat Rb Rabbit WR Most Common Vertebrates ( ) Predicted Reactivity  
Please visit [www.millipore.com/SmartFlare](http://www.millipore.com/SmartFlare) for additional product information.

EMD Millipore Corporation, 28820 Single Oak Drive, Temecula, CA 92590, USA 1-800-437-7500

Technical Support: T: 1-800-MILLIPORE (1-800-645-5476) • F: 1-800-437-7502

**FOR RESEARCH USE ONLY.** Not for use in diagnostic procedures. Not for human or animal consumption. Purchase of this Product does not include any right to resell or transfer, either as a stand-alone product or as a component of another product. Any use of this Product for purposes other than research is strictly prohibited.

EMD Millipore, the M mark and SmartFlare are registered trademarks of Merck KGaA, Darmstadt, Germany.

© 2012 EMD Millipore Corporation, Billerica MA USA. All rights reserved.



We Buy 100% Certified  
Renewable Energy

Rev A/2013/03/12/SF-753DSCA/AK

