


**RABBIT ANTI-SV2B
POLYCLONAL ANTIBODY**

CATALOG NUMBER: AB15225
LOT NUMBER:
QUANTITY: 100 μ L
SPECIFICITY: SV2B. By Western blot the antibody recognizes a band at ~78 kDa in rat brain extract.
IMMUNOGEN: Synthetic peptide from rat.
APPLICATIONS: Western blot: 1:5,000-1:10,000 on extracts from rat brain.

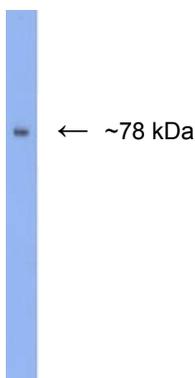

Optimal working dilutions must be determined by the end user.

SPECIES REACTIVITY: Rat. Other species have not been tested.
FORMAT: Rabbit serum.
PRESENTATION: Liquid.
STORAGE/HANDLING: Maintain at -20°C in undiluted aliquots for up to 12 months after date of receipt. Avoid repeated freeze/thaw cycles.

Important Note: *During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 μ L or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the container's cap.*

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.

For research use only; not for use as a diagnostic.

© 2006: CHEMICON International, Inc. - By CHEMICON International, Inc. All rights reserved. No part of these works may be reproduced in any form without permissions in writing.