

RABBIT ANTI-RabGDI alpha AFFNITY PURIFIED POLYCLONAL ANTIBODY

CATALOG NUMBER: AB15850

LOT NUMBER:

QUANTITY: $100 \mu g$


CONCENTRATION: 1 mg/mL

SPECIFICITY: RabGDI alpha (Rab GDP dissociation inhibitor alpha). By Western blot the antibody reacts

with a band at ~ 51-55 kDa on rat brain lysate.

IMMUNOGEN: Synthetic peptide from the C-terminus of rat RabGDI alpha.

APPLICATIONS: Western blot: 0.01-0.05 μg/mL on rat brain lysate.

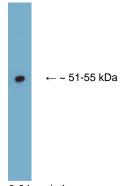

 $0.01~\mu g/mL$

Optimal working dilutions must be determined by the end user.

SPECIES REACTIVITY: Rat. Other species have not been tested. The immunogen sequence has 100% homology

with mouse.

FORMAT: Affinity purified immunoglobulin.

PRESENTATION: Liquid in PBS with 0.1% sodium azide.

STORAGE/HANDLING: Maintain at 2-8°C in undiluted aliquots for up to 6 months after date of receipt.

Important Note: During shipment, small volumes of product will occasionally become entrapped in the seal of the

product vial. For products with volumes of 200 μ L or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the

container's cap.

For research use only; not for use as a diagnostic.

© 2007: Millipore Corporation. All rights reserved. No part of these works may be reproduced in any form without permissions in writing.





For research use only; not for use as a diagnostic.

 $@\ 2007{:}\ Millipore\ Corporation.\ All\ rights\ reserved.\ No\ part\ of\ these\ works\ may\ be\ reproduced\ in\ any\ form\ without\ permissions\ in\ writing.$