


Stick with PureProteome™ magnetic beads.

More choices for flexible, fast affinity purification.

INTRODUCING TWO NEW PURIFICATION TOOLS:

1. Streptavidin magnetic beads capture biotinylated molecules.

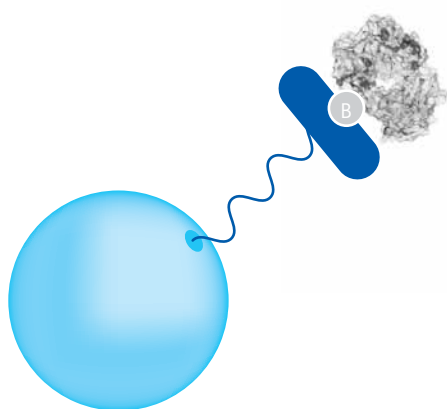


Streptavidin PureProteome™ Magnetic Beads

- SPECIFIC: Low nonspecific binding
- HIGH CAPACITY: Binds >20 µg Biotin-BSA/mg settled beads
- FLEXIBLE: Wide array of commercially available biotinylated targets

FREE Streptavidin bead samples available on request

2. NHS FlexiBind magnetic beads can be customized to capture your protein of interest.

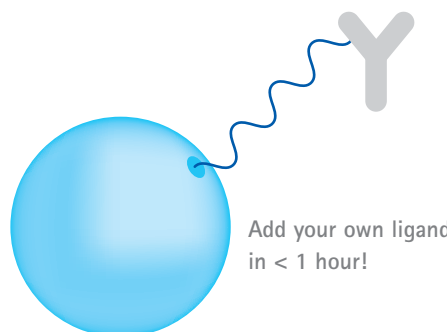

Add your own ligand in < 1 hour!

NHS FlexiBind PureProteome™ Magnetic Beads


- SIMPLE: Customize your own affinity beads for protein and other biomolecule ligands in <60 min
- EASY TO USE: Kit contains beads and all buffers for washing, coupling and quenching as well as Amicon® Ultra centrifugal devices for purification.
- FLEXIBLE: Reacts with any primary amine-containing molecule, including proteins, peptides, nucleic acids, lectins, and amine-containing drug targets.

Fast, simple, flexible affinity purification.

PureProteome™ NHS FlexiBind magnetic beads offer you flexibility in binding your target ligand. The kit contains everything you need, and offers high binding capacities (NHS density > 17 µmoles/mL of settled beads) with high specificity due to covalent linkages. NHS FlexiBind beads are perfect for applications involving targets that do not have affinity for common preconjugated magnetic beads.

For additional technical information on the NHS FlexiBind beads, visit: www.millipore.com/NHS

Customize your own bead in an hour!

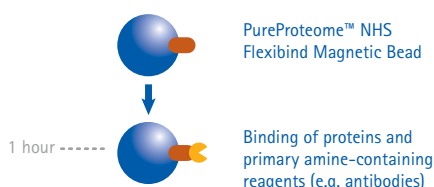

Description	Contents	Catalogue No.
PureProteome™ NHS FlexiBind Magnetic Beads Kit	0.5 mL NHS magnetic beads, 20% slurry 4 Pk Amicon® Ultra-0.5 Centrifugal Filters, 30,000 MWCO 5 mL Equilibration Buffer 25 mL Coupling Et Wash Buffer 5 mL Quenching Buffer	LSKMAGN01
PureProteome™ NHS FlexiBind Magnetic Beads	2 mL NHS magnetic beads, 20% slurry	LSKMAGN04
PureProteome™ Streptavidin Magnetic Beads	2 mL Streptavidin-conjugated magnetic beads, 5% slurry	LSKMAGT02
PureProteome™ Streptavidin Magnetic Beads	10 mL Streptavidin-conjugated magnetic beads, 5% slurry	LSKMAGT10

To Place an Order or Receive Technical Assistance

In Europe, please call Customer Service:

France: 0825 045 645 Italy: 848 845 645
Germany: 01805 045 645 UK: 0870 900 46 45
Spain: 901 516 645 Option 1 Switzerland: 0848 645 645

For other countries across Europe, please call: +44 (0) 115 943 0840,
or visit: www.merckmillipore.com/offices

For Technical Service, please visit: www.merckmillipore.com/techservice



www.merckmillipore.com

Get Connected!

Join Merck Millipore Bioscience on your favorite social media outlet for the latest updates, news, products, innovations, and contests!

 facebook.com/MerckMilliporeBioscience

 twitter.com/Merck4Bio