

3050 Spruce Street, St. Louis, MO 63103 USA
Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757
email: techservice@sial.com sigma-aldrich.com

# **Product Information**

PHKG1, active, GST tagged, human PRECISIO® Kinase recombinant, expressed in *Sf*9 cells

Catalog Number **SRP5061** Storage Temperature –70°C

Synonyms: PHKG

# **Product Description**

PHKG1 or phosphorylase kinase gamma 1 is a member of the Ser/Thr protein kinase family that encodes a protein with one protein kinase domain and two calmodulin-binding domains. Phosphorylase kinase is a crucial glycogenolytic regulatory enzyme. PHKG1 is the catalytic member of a 16-subunit protein kinase complex that contains equimolar ratios of 4 subunit types known as alpha, beta, gamma and delta. Skeletal muscle contains the highest amount of phosphorylase kinase enzymatic activity, although activity is also observed in liver, cardiac muscle, brain, and several other tissues.

Recombinant, full-length, human PHKG1 was expressed by baculovirus in *Sf*9 insect cells using an N-terminal GST tag. The gene accession number is NM\_006213. Recombinant protein stored in 50 mM Tris-HCl, pH 7.5, 150 mM NaCl, 10 mM glutathione, 0.1 mM EDTA, 0.25 mM DTT, 0.1 mM PMSF, and 25% glycerol.

Molecular mass: ~70 kDa

Purity: 70–95% (SDS-PAGE, see Figure 1)

Specific Activity: 25-33 nmole/min/mg (see Figure 2)

# **Precautions and Disclaimer**

This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

### Storage/Stability

The product ships on dry ice and storage at  $-70~^{\circ}$ C is recommended. After opening, aliquot into smaller quantities and store at  $-70~^{\circ}$ C. Avoid repeated handling and multiple freeze/thaw cycles.

Figure 1.

SDS-PAGE Gel of Typical Lot 70–95% (densitometry)


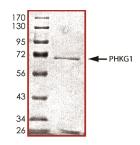


Figure 2.
Specific Activity of Typical Lot 25–33 nmole/min/mg

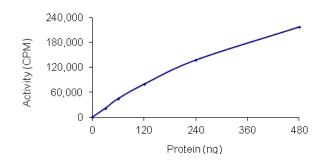

#### **Procedure**

# **Preparation Instructions**

Kinase Assay Buffer – 25 mM MOPS, pH 7.2, 12.5 mM glycerol 2-phosphate, 25 mM MgCl<sub>2</sub>, 5 mM EGTA, and 2 mM EDTA. Just prior to use, add DTT to a final concentration of 0.25 mM.

Kinase Dilution Buffer – Dilute the Kinase Assay Buffer 5-fold with a 50 ng/µl BSA.

Kinase Solution – Dilute the active PHKG1 (0.1  $\mu$ g/ $\mu$ l) with Kinase Dilution Buffer to the desired concentration. Note: The specific activity plot may be used as a guideline (see Figure 2). It is recommended the researcher perform a serial dilution of active PHKG1 kinase for optimal results.

10 mM ATP Stock Solution – Dissolve 55 mg of ATP in 10 ml of Kinase Assay Buffer. Store in 200  $\mu$ l aliquots at –20 °C.

 $\gamma$ -<sup>33</sup>P-ATP Assay Cocktail (250 μM) – Combine 5.75 ml of Kinase Assay Buffer, 150 μl of 10 mM ATP Stock Solution, 100 μl of  $\gamma$ -<sup>33</sup>P-ATP (1 mCi/100 μl). Store in 1 ml aliquots at –20 °C.

Substrate Solution – Dissolve the synthetic peptide substrate in distilled water at a final concentration of 1 mg/ml.

1% phosphoric acid solution – Dilute 10 ml of concentrated phosphoric acid to a final volume of 1 L with water.

# Kinase Assay

This assay involves the use of the <sup>33</sup>P radioisotope. All institutional guidelines regarding the use of radioisotopes should be followed.

- 1. Thaw the active PHKG1, Kinase Assay Buffer, Substrate Solution, and Kinase Dilution Buffer on ice. The  $\gamma$ -<sup>33</sup>P-ATP Assay Cocktail may be thawed at room temperature.
- In a pre-cooled microcentrifuge tube, add the following solutions to a volume of 20 μl:

10  $\mu$ l of Kinase Solution 5  $\mu$ l of Substrate Solution 5  $\mu$ l of cold water (4 °C)

- Set up a blank control as outlined in step 2, substituting 5 μl of cold water (4 °C) for the Substrate Solution.
- 4. Initiate each reaction with the addition of 5  $\mu$ l of the  $\gamma$ - $^{33}$ P-ATP Assay Cocktail, bringing the final reaction volume to 25  $\mu$ l. Incubate the mixture in a water bath at 30 °C for 15 minutes.
- 5. After the 15 minute incubation, stop the reaction by spotting 20  $\mu$ l of the reaction mixture onto an individually precut strip of phosphocellulose P81 paper.

- 6. Air dry the precut P81 strip and sequentially wash in the 1% phosphoric acid solution with constant gentle stirring. It is recommended the strips be washed a total of 3 times of ~10 minutes each.
- 7. Set up a radioactive control to measure the total  $\gamma$ - $^{33}$ P-ATP counts introduced into the reaction. Spot 5  $\mu$ l of the  $\gamma$ - $^{33}$ P-ATP Assay Cocktail on a precut P81 strip. Dry the sample for 2 minutes and read the counts. Do not wash this sample.
- Count the radioactivity on the P81 paper in the presence of scintillation fluid in a scintillation counter.
- 9. Determine the corrected cpm by subtracting the blank control value (see step 3) from each sample and calculate the kinase specific activity

#### Calculations:

1. Specific Radioactivity (SR) of ATP (cpm/nmole)

SR = 
$$\frac{\text{cpm of 5} \, \mu \text{l of } \gamma^{-33}\text{P-ATP Assay Cocktail}}{\text{nmole of ATP}}$$

cpm – value from control (step 7) nmole – 1.25 nmole (5 μl of 250 μM ATP Assay Cocktail)

2. Specific Kinase Activity (SA) (nmole/min/mg)

nmole/min/mg = 
$$\frac{\Delta cpm \times (25/20)}{SR \times E \times T}$$

SR = specific radioactivity of the ATP (cpm/nmole ATP)  $\Delta$ cpm = cpm of the sample – cpm of the blank (step 3) 25 = total reaction volume

20 = spot volume

T = reaction time (minutes)

E = amount of enzyme (mg)

## References

- Burwinkel, B. et al., Muscle glycogenosis with low phosphorylase kinase activity: mutations in PHKA1, PHKG1 or six other candidate genes explain only a minority of cases. Europ. J. Hum. Gen., 11, 516-526 (2003).
- Wehner, M. et al., Human cDNA encoding the muscle isoform of the phosphorylase kinase gamma subunit (PHKG1). Hum. Genet., 96, 616-618 (1995).

PRECISIO is a registered trademark of Sigma-Aldrich Co. LLC.

TD,MAM 11/11-1