3050 Spruce Street, St. Louis, MO 63103 USA
Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757
email: techservice@sial.com sigma-aldrich.com

Product Information

MEKK6, active, GST-tagged, human PRECISIO® Kinase recombinant, expressed in *Sf*9 cells

Catalog Number **SRP5294** Storage Temperature –70 °C

Synonyms: MAP3K6, MAPKKK6, ASK2

Product Description

MEKK6 or MAP3K6 (mitogen-activated protein kinase kinase kinase 6) is a member of the serine/threonine protein kinase family. MEKK6 was identified by its interaction with MAP3K5/ASK, a protein kinase and an activator of c-Jun kinase (MAPK7/JNK) and MAPK14/p38 kinase. MEKK6 cascades relay signals from the cell surface to the nucleus which result in gene transcription. MEKK6 activates the JNK and p38 pathway, and induces apoptosis when overexpressed in stably transfected cells. The MEKK6 transcript is expressed in several tissues, with the highest expression in heart and skeletal muscle.

Recombinant human MEKK6 (620-950) was expressed by baculovirus in *Sf*9 insect cells using an N-terminal GST-tag. The MEKK6 gene accession number is NM_004672. It is supplied in 50 mM Tris-HCl, pH 7.5, 150 mM NaCl, 10 mM glutathione, 0.1 mM EDTA, 0.25 mM DTT, 0.1 mM PMSF, and 25% glycerol.

Molecular mass: ~67 kDa

Precautions and Disclaimer

This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

The product ships on dry ice and storage at $-70~^{\circ}$ C is recommended. After opening, aliquot into smaller quantities and store at $-70~^{\circ}$ C. Avoid repeated handling and multiple freeze/thaw cycles.

Figure 1.

SDS-PAGE Gel of Typical Lot:

≥70% (SDS-PAGE, densitometry)


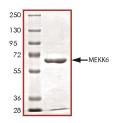


Figure 2.Specific Activity of Typical Lot: 8–12 nmole/min/mg

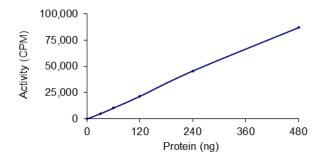

Procedure

Preparation Instructions

Kinase Assay Buffer – 25 mM MOPS, pH 7. 2, 12.5 mM glycerol 2-phosphate, 25 mM MgC1₂, 5 mM EGTA, and 2 mM EDTA. Just prior to use, add DTT to a final concentration of 0.25 mM.

Kinase Dilution Buffer – Dilute the Kinase Assay Buffer 5-fold with a 50 ng/μL BSA solution.

Kinase Solution – Dilute the active MEKK6 (0.1 μ g/ μ L) with Kinase Dilution Buffer to the desired concentration. Note: The specific activity plot may be used as a guideline (see Figure 2). It is recommended the researcher perform a serial dilution of active MEKK6 kinase for optimal results.

10 mM ATP Stock Solution – Dissolve 55 mg of ATP in 10 mL of Kinase Assay Buffer. Store in 200 μ L aliquots at –20 °C.

 $\gamma\text{-}^{33}\text{P-ATP}$ Assay Cocktail (250 $\mu\text{M})$ – Combine 5.75 mL of Kinase Assay Buffer, 150 μL of 10 mM ATP Stock Solution, 100 μL of $\gamma\text{-}^{33}\text{P-ATP}$ (1 mCi/100 μL). Store in 1 mL aliquots at –20 °C.

Substrate Solution – MBP Protein substrate diluted in distilled water to a final concentration of 1 mg/mL.

1% phosphoric acid solution – Dilute 10 mL of concentrated phosphoric acid to a final volume of 1 L with water.

Kinase Assay

This assay involves the use of the ³³P radioisotope. All institutional guidelines regarding the use of radioisotopes should be followed.

- 1. Thaw the active MEKK6, Kinase Assay Buffer, Substrate Solution, and Kinase Dilution Buffer on ice. The γ -³³P-ATP Assay Cocktail may be thawed at room temperature.
- 2. In a pre-cooled microcentrifuge tube, add the following solutions to a volume of 20 μ L:

10 μL of Kinase Solution

5 μL of Substrate Solution

5 μL of cold water (4 °C)

- 3. Set up a blank control as outlined in step 2, substituting 5 μ L of cold water (4 °C) for the Substrate Solution.
- 4. Initiate each reaction with the addition of 5 μ L of the γ - 33 P-ATP Assay Cocktail, bringing the final reaction volume to 25 μ L. Incubate the mixture in a water bath at 30 °C for 15 minutes.
- After the 15 minute incubation, stop the reaction by spotting 20 μL of the reaction mixture onto an individually precut strip of phosphocellulose P81 paper.

- Air dry the precut P81 strip and sequentially wash in the 1% phosphoric acid solution with constant gentle stirring. It is recommended the strips be washed a total of 3 times of ~10 minutes each.
- 7. Set up a radioactive control to measure the total γ - 33 P-ATP counts introduced into the reaction. Spot 5 μ L of the γ - 33 P-ATP Assay Cocktail on a precut P81 strip. Dry the sample for 2 minutes and read the counts. Do not wash this sample.
- 8. Count the radioactivity on the P81 paper in the presence of scintillation fluid in a scintillation counter.
- 9. Determine the corrected cpm by subtracting the blank control value (see step 3) from each sample and calculate the kinase specific activity

Calculations:

1. Specific Radioactivity (SR) of ATP (cpm/nmole)

SR =
$$\frac{\text{cpm of 5} \mu \text{L of } \gamma^{-33} \text{P-ATP Assay Cocktail}}{\text{nmole of ATP}}$$

cpm – value from control (step 7) nmole – 1.25 nmole (5 μ L of 250 μ M ATP Assay Cocktail)

2. Specific Kinase Activity (SA) (nmole/min/mg)

nmole/min/mg =
$$\Delta cpm \times (25/20)$$

SR \times E \times T

SR = specific radioactivity of the ATP (cpm/nmole ATP) Δ cpm = cpm of the sample – cpm of the blank (step 3) 25 = total reaction volume

20 = spot volume

T = reaction time (minutes)

E = amount of enzyme (mg)

References

- Wang, X.S. et al., MAPKKK6, a novel mitogenactivated protein kinase kinase kinase, that associates with MAPKKK5. Biochem. Biophys. Res. Commun., 253, 33-37 (1998).
- Takeda, K. et al., Apoptosis signal-regulating kinase (ASK) 2 functions as a mitogen-activated protein kinase kinase kinase in a heteromeric complex with ASK1. J. Biol. Chem., 282(10), 7522-31 (2007).

PRECISIO is a registered trademark of Sigma-Aldrich Co. LLC.

RC.MAM 12/12-1