

#### sigma-aldrich.com

3050 Spruce Street, St. Louis, MO 63103 USA Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757 email: techservice@sial.com sigma-aldrich.com

# **Product Information**

Membrane Scaffold Protein 1E1D1 recombinant, expressed in *E. coli* buffered aqueous solution

Catalog Number **MSP04** Storage Temperature –20 °C

Synonym: MSP1E1D1

### **Product Description**

Nanodisc technology is a widely applicable approach to render membrane proteins soluble in aqueous solutions in a native-like bilayer environment, where the membrane proteins remain stable and active.

The Nanodisc concept is derived from high density lipoprotein (HDL) particles and their primary protein component, apolipoprotein. The Nanodisc is a non-covalent structure of phospholipid bilayer and membrane scaffold protein (MSP), a genetically engineered protein, which mimics the function of Apolipoprotein A-1 (ApoA-1).<sup>1,2</sup> A soluble Nanodisc assembles as the phospholipid forms a bilayer, which is encircled by two amphipathic MSP molecules covering the hydrophobic alkyl chains of the bilayer. The length of the MSP controls the size of the Nanodisc structure.<sup>3</sup>


Incorporation of a membrane protein target into a Nanodisc is accomplished by first solubilizing the target protein in detergent, and then mixing with cholatesolubilized phospholipid and MSP. Removal of the detergents results in the assembly of the Nanodisc structure around the target protein (see Figure 1).

The resulting structure renders the protein soluble in a model membrane system of defined phospholipid content and size with a native-like bilayer. This provides stability and functional requirements for the protein target. The size of the Nanodisc, determined by the length of the MSP, allows control of the oligomerization state of the membrane protein.

Some examples of membrane-bound proteins incorporated into Nanodiscs are GPCRs,<sup>4</sup> cytochromes P450,<sup>5-7</sup> bacteriorhodopsin,<sup>8</sup> coagulation factors,<sup>9</sup> cholera toxin,<sup>10</sup> and TAR receptor.<sup>11</sup>

#### Figure 1.

Example of a Nanodisc containing a 7-transmembrane protein

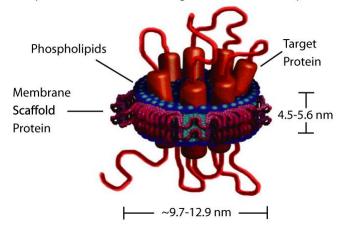

A critical component of Nanodiscs is the membrane scaffold protein (MSP), which forms the encircling amphipathic helical protein belt.<sup>12,13</sup> The sequence similarities of MSPs to Apolipoprotein A-1 are presented in Figure 2.

### Figure 2.

Protein Maps of Apolipoprotein A-1 and MSP1E1D1

			Арс	olipo	protei	n A1					
Globu	ilar Region	H1	H2	H3	H4	H5	H6	H7	H8	H9	H10
TE	4 V protease ognition si		100 - 112 	88 99 ma M	121 ISP04 :	143 = MSF	<sup>16</sup> 1E1	8552 8833 Marca	38 20	09 220	243
Hist. tag	TEV H	0.5 H2	НЗ	H4	H4	H5	Н	6 H	7 Н	8 H	9 H10
1 9	23	34	56 6	7	89 1	11	133	155	177	199	210 233

The first MSP, MSP1, was engineered with its sequence based on the sequence of A-1, but without the globular N-terminal domain of native A-1.<sup>14</sup> The MSP1E1D1 variant of MSP1 deletes the first 11 amino acids in the Helix 1 portion (referred to as "H0.5" in Figure 2) of the original MSP1 sequence <sup>3</sup> (which is known separately as MSP1D1). The MSP1E1D1 protein is supplied as an N-terminal histidine-tagged protein with a TEV protease cleavage site between the His-tag and the protein sequence. The MSP1E1D1 variant also repeats the Helix 4 sequence of the original MSP1 sequence of the original MSP1 sequence between the parent Helix 4 and Helix 5 segments of MSP1D1.<sup>15</sup>

MSP1E1D1 yields Nanodiscs of ~10.5 nm.<sup>16</sup> The number of lipid molecules in an MSP1E1D1 can range from 180–250.<sup>17</sup> The thickness of a Nanodisc is dependent on the type of phospholipid incorporated (typically 4.6–5.6 nm).<sup>3</sup>

This product is supplied as a buffered aqueous solution, at a protein concentration of ~5 mg/mL, and containing 20 mM Tris-HCl, pH 7.4, 100 mM NaCl, 0.5 mM EDTA, 5 mM sodium cholate, and 0.1% sodium azide.

## MSP1E1D1 sequence:<sup>18</sup>

GHHHHHHHDYDIPTTENLYFQGSTFSKLREQLGPVT QEFWDNLEKETEGLRQEMSKDLEEVKAKVQPYLDDF QKKWQEEMELYRQKVEPLRAELQEGARQKLHELQE KLSPLGEEMRDRARAHVDALRTHLAPYLDDFQKKWQ EEMELYRQKVEPLRAELQEGARQKLHELQEKLSPLG EEMRDRARAHVDALRTHLAPYSDELRQRLAARLEAL KENGGARLAEYHAKATEHLSTLSEKAKPALEDLRQGL LPVLESFKVSFLSALEEYTKKLNTQ

Molecular mass: 27,546.9 Da (based on sequence) Molecular mass of histidine-tag cleaved variant: 24,929.1 (based on sequence)

The MSP1E1D1 protein concentration can be determined spectrophotometrically using the following extinction coefficients:

## MSP1E1D1, $\varepsilon_{280} = 29,910 \text{ M}^{-1} \text{ cm}^{-1}$

Histidine-tag cleaved variant,  $\varepsilon_{280} = 26,930 \text{ M}^{-1} \text{ cm}^{-1}$  (Solvent: 20 mM Tris, pH 7.4, with 0.1 M NaCl, 0.5 mM EDTA, and 0.01% NaN<sub>3</sub>)

## **Reagents Required, but Not Provided**

- Sodium Cholate (Catalog No. C6445)
- Amberlite<sup>®</sup> XAD2 (Catalog No. 20275)
- Phospholipid(s)
  - 2-Oleoyl-1-palmitoyl-*sn*-glycero 3-phosphocholine (POPC, Catalog No. P3017)
  - 1,2-Dipalmitoyl-*rac*-glycero-3-phosphocholine (DPPC, Catalog No. P5911)
  - 1,2-Dimyristoyl-*rac*-glycero-3-phosphocholine (DMPC, Catalog No. P7930)
- Triton™ X-100 (optional, Catalog No. X100)
- Superdex<sup>®</sup> Gel Filtration Column (Catalog No. 54801-U) or Superdex 200 (Catalog No. S6782)

## **Precautions and Disclaimer**

This product is for R&D use only, not for drug, household, or other uses. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

## Storage/Stability

The product ships at 2-8 °C. Storage at -20 °C is recommended. When stored at -20 °C, the protein remains functional for at least 1 year.<sup>15</sup>

The MSP1E1D1 Stock Solution can be stored at 2–8 °C for up to several days. Recommended long-term storage is in frozen aliquots at –20 °C.<sup>15</sup>

## Procedure

The successful assembly of Nanodiscs is highly dependent on the properties of the target protein. The choice of the types and ratios of MSPs, phospholipids, and detergents might require trials and titrations to optimize disc assembly. A thorough review of the literature is essential for this aspect of Nanodisc construction.<sup>3,12,13,16</sup>

Successful incorporation of the target membrane protein into Nanodisc most likely will require optimization of the following important parameters:

- Choice and concentration of secondary detergent for preparation of target membrane protein
- Choice and preparation of phospholipid
- Temperature of reconstitution
- Lipid to MSP ratio
- MSP to target protein ratio

#### References

- Brouillette, C.G., et al., Biochim. Biophys. Acta, 1531(1-2), 4-46 (2001).
- 2. Jonas, A., *Methods Enzymol.*, **128**, 553-582 (1986).
- 3. Denisov, I.G., *et al.*, *J. Am. Chem. Soc.*, **126**(11), 3477–3487 (2004).
- Bayburt, T.H., et al., J. Biol. Chem., 282(20), 14875-14881 (2007).
- Denisov, I.G., et al., J. Biol. Chem., 282(10), 7066-7076 (2007).
- Grinkova, Y.V., et al., Biochem. Biophys. Res. Commun., 372(2), 379-382 (2008).
- Denisov, I.G., et al., Biochim Biophys Acta, 1814(1), 223-229 (2011).
- Bayburt, T.H., *et al.*, *Arch. Biochem. Biophys.*, 450(2), 215-222 (2006).
- Shaw, A.W., et al., J. Biol. Chem., 282(9), 6556-6563 (2007).
- 10. Borch, J., *et al.*, *Anal. Chem.*, **80(16)**, 6245-6252 (2008).
- Boldog, T., et al., Proc. Nat. Acad. Sci. USA, 103(31), 11509-11514 (2006).
- Nath, A., et al., Biochemistry, 46(8), 2059-2069 (2007).
- 13. Bayburt, T.H., and Sligar, S.G., *FEBS Lett.*, **584(9)**,1721-1727 (2010).
- 14. Bayburt, T.H., *et al., Nano Letters*, **2(8)**, 853-856 (2002).
- Grinkova, Y.V., et al., Protein Eng. Des. Sel., 23(11), 843-848 (2010).
- Ritchie, T.K., et al., Methods Enzymol., 464, 211-231 (2009).
- 17. Siuda, I., and Tieleman, D.P., *J. Chem. Theory Comput.*, **11(10)**, 4923-4932 (2015).
- 18. Supplier information.

For a further and extensive list of citations and protocols, visit the website of Professor Stephen Sligar, University of Illinois: http://sligarlab.life.uiuc.edu/nanodisc.html Nanodisc technology, and many of its uses, are covered by the following patents held by the University of Illinois:

- 7,691,414: Membrane scaffold proteins
- 7,662,410: Membrane scaffold proteins and embedded membrane proteins
- 7,622,437: Tissue factor compositions and methods
- 7,592,008: Membrane scaffold proteins
- 7,575,763: Membrane scaffold proteins and tethered membrane proteins
- 7,083,958: Membrane scaffold proteins
- 7,048,949: Membrane scaffold proteins

Triton is a trademark of The Dow Chemical Company or an affiliated company of Dow.

Amberlite and XAD are registered trademarks of The Dow Chemical Company or an affiliated company of Dow.

Superdex is a registered trademark of GE Healthcare.

GCY, RBG, TD, RC, MAM 02/17-1

©2017 Sigma-Aldrich Co. LLC. All rights reserved. SIGMA-ALDRICH is a trademark of Sigma-Aldrich Co. LLC, registered in the US and other countries. Sigma brand products are sold through Sigma-Aldrich, Inc. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see product information on the Sigma-Aldrich website at www.sigmaaldrich.com and/or on the reverse side of the invoice or packing slip.