


Product Information

COCAINE HYDROCHLORIDE


Methanol Solution

Product Number	C 1528	Synonyms	3-(Benzyloxy)-8-azabicyclo-[3.2.1]octane-2-carboxylic acid methyl ester hydrochloride; l-Cocaine hydrochloride; Ecgonine methyl ester bezoate hydrochloride; Benzoylmethylecgonine hydrochloride
Lot Number	084K8800	Stability	If stored at 2-8°C in an air tight container, protected from light, decomposition of the product is less than 1% in 24 months. After opening, concentration may change due to loss of solvent.
CAS Number	[53-21-4]	Expiration Date	April 2007
Concentration	1.12 mg Cocaine Hydrochloride (equivalent to 1.00 mg free base) /mL GC grade methanol (conc. verified by UV)	Storage	2-8°C.
Molecular Formula	C ₁₇ H ₂₁ NO ₄ • HCl		
Molecular Weight	339.8		
DEA Class	Exempt preparation of Class II	70 eV Electron Ionization Mass Spectrum	

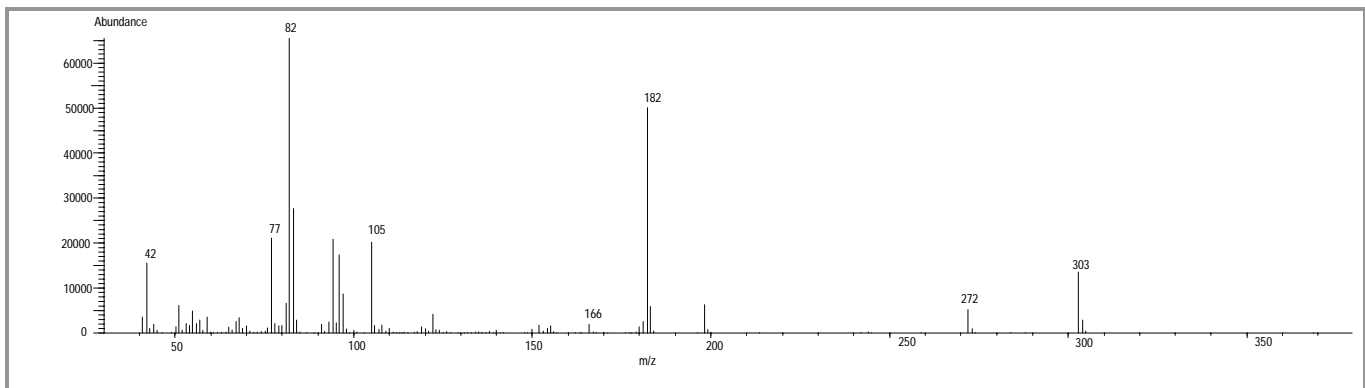

Hewlett-Packard 5890 Series II Gas Chromatograph
 Hewlett-Packard 5971 Mass Selective Detector
 Hewlett-Packard G1034C MS ChemStation (HP-DOS series)

FT-IR

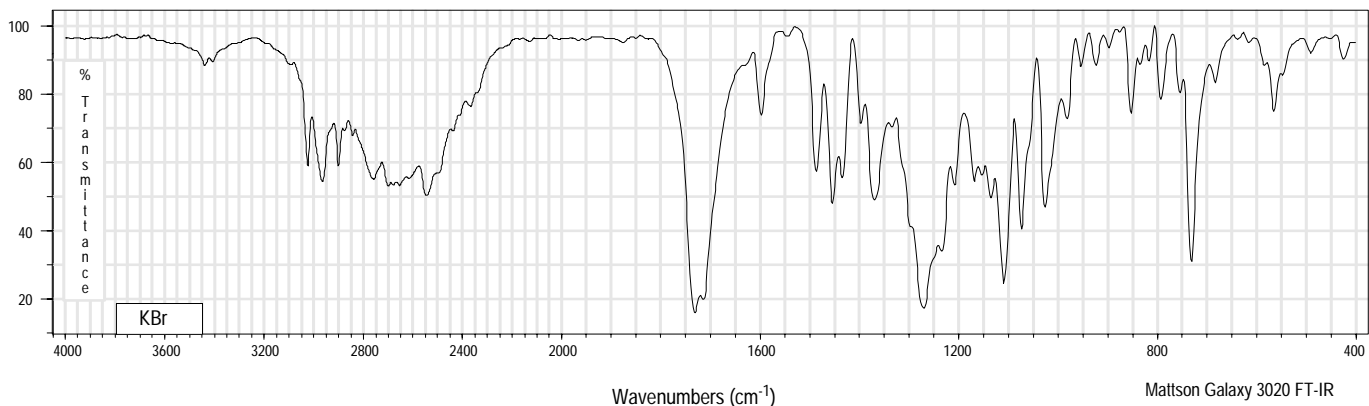


COCAINE HYDROCHLORIDE

Methanol Solution


Product Number C 1528

Lot Number 084K8800

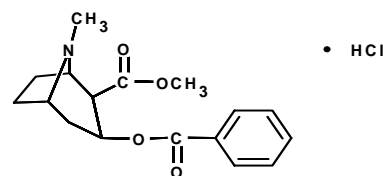

Ultraviolet Spectrum


Peak	233.8 nm
Absorbance	0.915
Solvent	0.2 N H ₂ SO ₄
Dilution	1:49 (v/v)

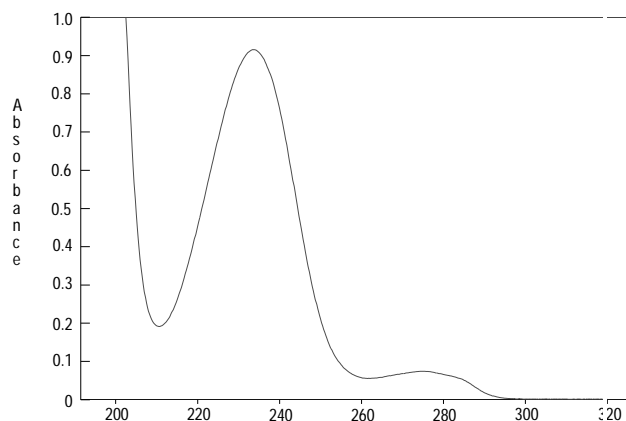

Bio-Tek Uvikon XL

Capillary Gas Chromatography

Column Supelco SPB-1 (Cat. No. 2-4162)
I.D. 0.20 mm **Length** 15 m **Film Thickness** 0.20 μ m
Oven Temperature 50°C (2 min.), then 15°/min. to 300°C, hold 11 min.
Injector Temperature 250°C, Split
Detector Hewlett-Packard 5971
Mass Selective Detector (70 eV, 40-500 m/z)

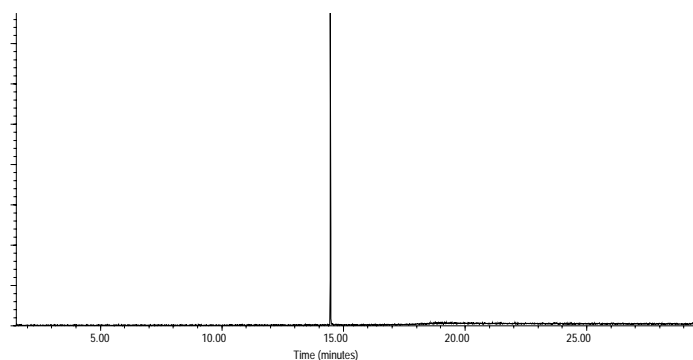

Hewlett-Packard 5890

C 1528/084K8800/08/04/04

This product is not for *in vitro* diagnostic procedures. For R&D use only. Not for drug, household or other uses.
Sigma brand products are sold through Sigma-Aldrich, Inc.
Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications.
Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply.
Please see reverse side of the invoice or packing slip.