


Anti-HDAC7 (KG-17) (H 2662)

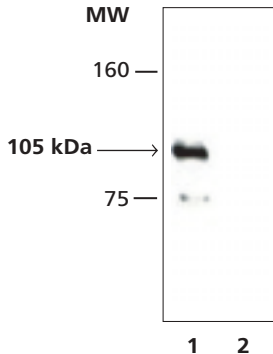

Immunoprecipitation

HDAC7 was immunoprecipitated from whole cell extract of human HEK 293T cells expressing HDAC7-FLAG using Rabbit Anti-HDAC7 (KG-17) (Product No. H 2662). The IP samples were separated on SDS-PAGE and blotted with Rabbit Anti-FLAG™ (F 7425). The antibody was developed with Goat Anti-Rabbit IgG, Peroxidase Conjugate (Product No. A 0545) and chemiluminescent substrate.

Lane 1: 1.8 µg Antibody immunoprecipitates HDAC7 from a RIPA extract of 180 µg of 293T cells expressing HDAC7-FLAG.

Lane 2: negative control (extract without antibody)