

ProductInformation

TESTOSTERONE STEROID

Methanol Solution

Concentration

Product NumberT 5411Synonyms4-Androsten-17β-ol-3-one;17β-Hydroxy-4-

Lot Number 065K8803 androsten-3-one; trans-Testosterone

CAS Number [58-22-0] Stability If stored at 2-8°C in an air tight container,

protected from light, decomposition of the product is less than 1% in 12 months. After 1.00 mg Testosterone/mL HPLC grade

opening, concentration may change due to


1,2- dimethoxyethane (conc. Verified loss of solvent.

Molecular Formula $C_{19}H_{28}O_2$ Expiration Date December 2007


Molecular Weight 288.4 Storage 2-8°C

DEA Class Exempt preparation of Class III

70 eV Electron Ionization Mass Spectrum

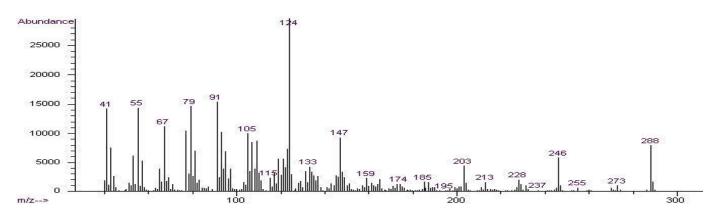

Hewlett-Packard 5890 Series II Gas Chromatograph Hewlett-Packard 5971 Mass Selective Detector Hewlett-Packard G1034C MS ChemStation (HP-DOS series)

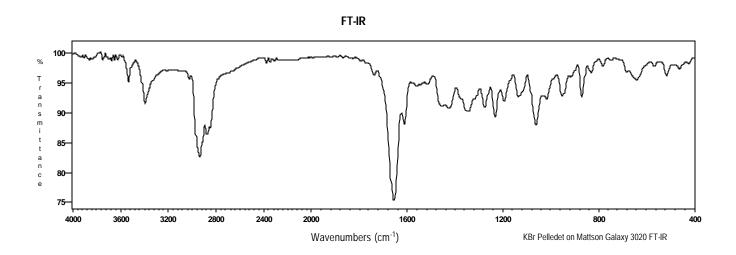


TESTOSTERONE STEROID

Methanol Solution **Product Number** T 5411 **Lot Number** 065K8803

CH₃ OH

Ultraviolet Spectrum

Peak	241.0 nm
Absorbance	0.569
Solvent	Methanol
Dilution	1:99 (v/v)

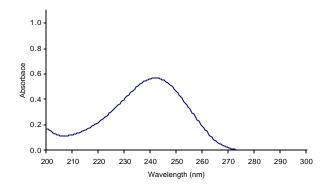

Bio-Tek Uvikon XL


Capillary Gas Chromatography

 $\begin{tabular}{ll} \textbf{Column} & Supelco SPB-1 (Cat. No. 2-4162) \\ \textbf{I.D.} & 0.20 \ mm & \textbf{Length} & 15 \ m \end{tabular} \begin{tabular}{ll} \textbf{Film Thickness} & 0.20 \ \mu m \\ \textbf{Oven Temperature} & 150 \end{tabular} \begin{tabular}{ll} \textbf{Cat. No. 2-4162)} \\ \textbf{No. 2-4162)} \\$

Injector Temperature 250°C, Split

Detector Hewlett-Packard 5971

Mass Selective Detector (70 eV, 40-500 m/z)

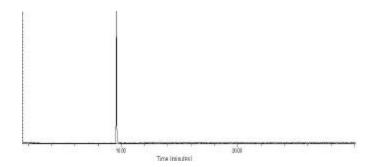

Hewlett-Packard 5890

T 5411/065k8803/07/13/05

This product is not for *in vitro* diagnostic procedures. For R&D use only. Not for drug, household or other uses. Sigma brand products are sold through Sigma-Aldrich, Inc.

Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply.

Please see reverse side of the invoice or packing slip.