




## Product Information

### SIGMA QUALITY CONTROL TEST PROCEDURE

#### Enzymatic Assay of TRYPSIN<sup>1</sup> (EC 3.4.21.4)

##### PRINCIPLE:

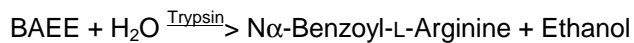

Abbreviation used:

BAEE = N $\alpha$ -Benzoyl-L-Arginine Ethyl Ester

**CONDITIONS:** T = 25°C, pH = 7.6, A<sub>253nm</sub>, Light path = 1 cm

**METHOD:** Continuous Spectrophotometric Rate Determination

##### REAGENTS:

- A. 67 mM Sodium Phosphate Buffer, pH 7.6 at 25°C  
(Prepare 100 ml in deionized water using Sodium Phosphate, Monobasic, Anhydrous, Sigma Prod. No. S-0751. Adjust to pH 7.6 at 25°C with 1 M NaOH.)
- B. 0.25 mM N $\alpha$ -Benzoyl-L-Arginine Ethyl Ester Solution (BAEE)  
(Prepare 50 ml in Reagent A using N $\alpha$ -Benzoyl-L-Arginine Ethyl Ester, Hydrochloride, Sigma Prod. No. B-4500.)
- C. 1 mM Hydrochloric Acid Solution (HCl)  
(Prepare 50 ml in deionized water using concentrated Hydrochloric Acid, Sigma Prod. No. H-7020.)
- D. Trypsin Enzyme Solution  
(Immediately before use, prepare a solution containing 500 BAEE units/ml of Trypsin in cold Reagent C.)

## Enzymatic Assay of TRYPsin<sup>1</sup> (EC 3.4.21.4)

### PROCEDURE:

Pipette (in milliliters) the following reagents into suitable quartz cuvettes:

	Test	Blank
Reagent B (BAEE)	3.00	3.00

Equilibrate to 25°C. Monitor the  $A_{253\text{nm}}$  until constant, using a suitably thermostatted spectrophotometer. Then add:

Reagent C (Hcl)	-----	0.20
Reagent D (Enzyme Solution)	0.20	-----

Immediately mix by inversion and record the increase in  $A_{253\text{nm}}$  for approximately 5 minutes. Obtain the  $\Delta A_{253\text{nm}}/\text{minute}$  using the maximum linear rate for both the Test and Blank.

### CALCULATIONS:

$$\text{BAEE Units/ml enzyme} = \frac{(\Delta A_{253\text{nm}}/\text{min Test} - \Delta A_{253\text{nm}}/\text{min Blank})(\text{df})}{(0.001)(0.20)}$$

df = Dilution factor

0.001 = The change in  $A_{253\text{nm}}/\text{minute}$  per unit of Trypsin at pH 7.6 at 25°C in a 3.2 ml reaction mix

0.20 = Volume (in milliliters) of enzyme used

$$\text{Units/mg solid} = \frac{\text{units/ml enzyme}}{\text{mg solid/ml enzyme}}$$

$$\text{Units/mg protein} = \frac{\text{units/ml enzyme}}{\text{mg protein/ml enzyme}}$$

### UNIT DEFINITION:

One BAEE unit will produce a  $\Delta A_{253\text{nm}}$  of 0.001 per minute with BAEE as substrate at pH 7.6 at 25°C in a reaction volume of 3.2 ml.

### FINAL ASSAY CONCENTRATION:

In a 3.2 ml reaction mix, the final concentrations are 63 mM sodium phosphate, 0.23 mM N $\alpha$ -benzoyl-L-arginine ethyl ester, 0.06 mM hydrochloric acid, and 100 units trypsin.

## Enzymatic Assay of TRYPsin<sup>1</sup> (EC 3.4.21.4)

### REFERENCE:

Bergmeyer, H.U., Gawehn, K., and Grassl, M. (1974) in *Methods of Enzymatic Analysis* (Bergmeyer, H.U., ed) Volume I, 2nd ed., 515-516, Academic Press, Inc., New York, NY

### NOTES:

1. This assay procedure is not to be used to assay Sigma Prod. Nos. T-1763, T-4019, T-8386, T-8899, and T-9906.
2. This assay is based on the cited reference.
3. Where Sigma Product or Stock numbers are specified, equivalent reagents may be substituted.

**Sigma warrants that the above procedure information is currently utilized at Sigma and that Sigma products conform to the information in Sigma publications. Purchaser must determine the suitability of the information and products for its particular use. Upon purchase of Sigma products, see reverse side of invoice or packing slip for additional terms and conditions of sale.**