3050 Spruce Street, St. Louis, MO 63103 USA
Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757
email: techservice@sial.com sigma-aldrich.com

Product Information

Psilocin solution 100 μg/mL±5% in methanol drug standard

Catalog Number **P4054**Lot Number 039K8720
Storage Temperature –20 °C

CAS RN 520-53-6

Synonyms: 3-[2-(Dimethylamino)ethyl]-1-H-indol-4-ol; 4-Hydroxy-N,N-dimethytryptamine

Product Description

Molecular Formula: C₁₂H₁₆N₂O Molecular Weight: 204.27

DEA Class: Exempt preparation of Class I

Concentration: 104 μg of psilocin/mL of GC grade methanol, concentration verified by UV.

Expiration Date: April 2013

Precautions and Disclaimer


This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

Store the product at -20 °C. If stored at -20 °C in an airtight container, protected from light, decomposition of the product is less than 1% in 36 months. After opening, concentration may change due to loss of solvent.


CMW, MAM 05/09-1

70 eV Electron Ionization Mass Spectrum

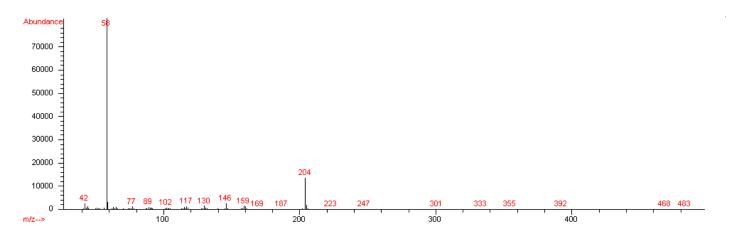

Hewlett-Packard 5890 Series II Gas Chromatograph Hewlett-Packard 5971 Mass Selective Detector Hewlett-Packard G1034C MS ChemStation (HP-DOS series)



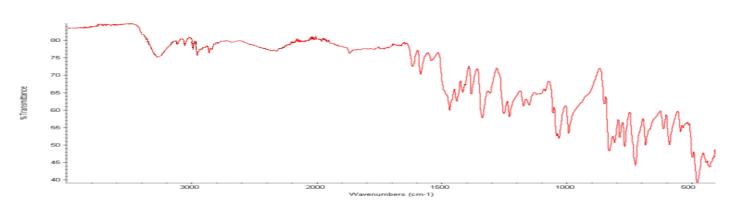


Nicolet 380 FT-IR Thermo Smart Orbit


Ultraviolet Spectrum

Peak 266.5 nm Absorbance 0.4409 Solvent Methanol Dilution 1:9 (v/v)

Capillary Gas Chromatography

Column Supelco SPB-1 (Catalog Number 24162)
I.D. 0.20 mm Length 15 m Film Thickness 0.20 μm
Oven Temperature 50 °C (2 minutes), then
15 °C/minute to 300 °C, hold 11 minutes
Injector Temperature 250 °C, Split
Detector Hewlett-Packard 5971
Mass Selective Detector (70 eV, 40-500 m/z)

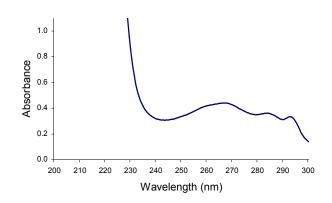

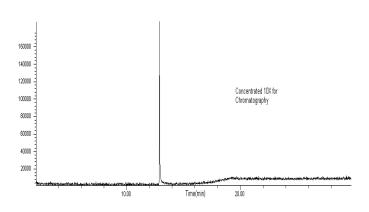

Bio-Tek Uvikon XL

Hewlett-Packard 5890

Sigma brand products are sold through Sigma-Aldrich, Inc.
Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser