


Product Information

Heroin hydrochloride solution
100 µg/mL±5% in methanol
drug standard

Catalog Number **H5144**
Lot Number 050M8705V
Storage Temperature -20 °C

CAS RN 1502-95-0
Synonyms: Diacetylmorphine Hydrochloride,
Acetomorphine Hydrochloride, Diamorphine
Hydrochloride

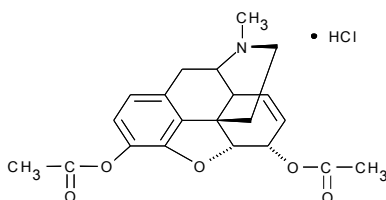

Product Description

Molecular Formula: C₂₁H₂₃NO₅ · HCl
Molecular Weight: 405.87

DEA Class: Exempt preparation of Class I

Concentration: 0.114 mg of heroin hydrochloride
(equivalent to 0.104 mg of free base)/mL of GC grade
methanol (concentration verified by UV)

Expiration Date: August 2013

Precautions and Disclaimer


This product is for R&D use only, not for drug,
household, or other uses. Please consult the Material
Safety Data Sheet for information regarding hazards
and safe handling practices.

Storage/Stability

Store the product at -20 °C. If stored at -20 °C in an
airtight container, protected from light, decomposition of
the product is less than 1% in 12 months. After
opening, concentration may change due to loss of
solvent.

CMW,MAM 12/12-1

70 eV Electron Ionization Mass Spectrum

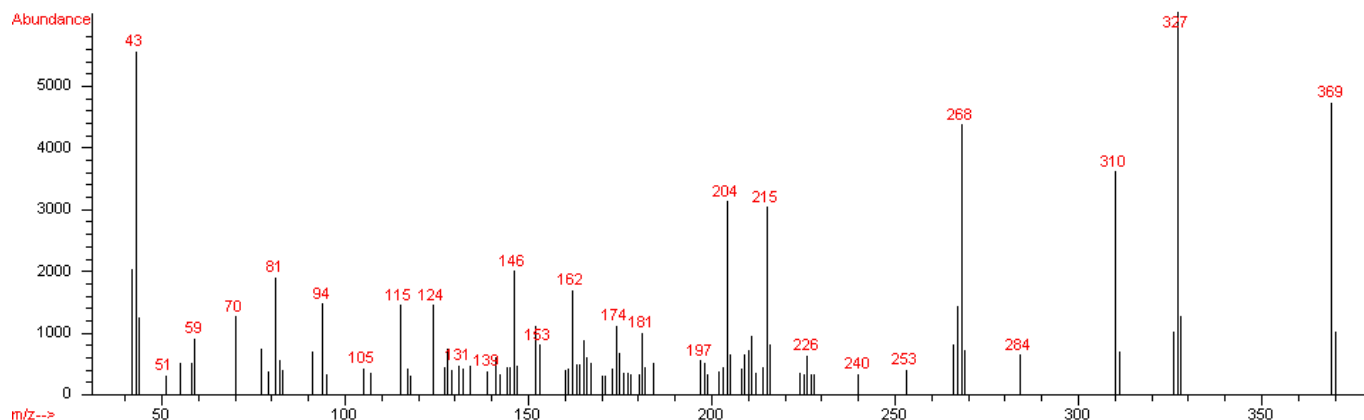

Hewlett-Packard 5890 Series II Gas Chromatograph
Hewlett-Packard 5971 Mass Selective Detector
Hewlett-Packard G1034C MS ChemStation (HP-DOS series)


FT-IR



Nicolet 380 FT-IR Smart Orbit(ATR)

Ultraviolet Spectrum


Peak	279.5 nm
Absorbance	0.5966
Dilution	Neat

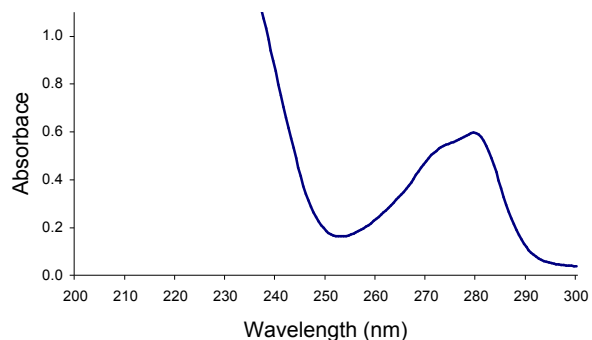

Beckman-Coulter DU 800

Capillary Gas Chromatography

Column Supelco SPB-1 (Catalog Number 24162)
I.D. 0.20 mm **Length** 15 m **Film Thickness** 0.20 μ m
Oven Temperature 50 °C (2 minutes), then
15 °C/minute to 300 °C, hold 11 minutes
Injector Temperature 250 °C, Split
Detector Hewlett-Packard 5971
Mass Selective Detector (70 eV, 40-500 m/z)

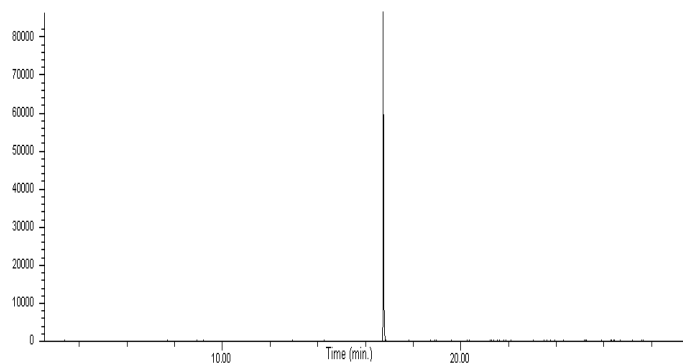

Hewlett-Packard 5890

Sigma brand products are sold through Sigma-Aldrich, Inc.

Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip.