


Anti-VDAC/Porin (V 2139)

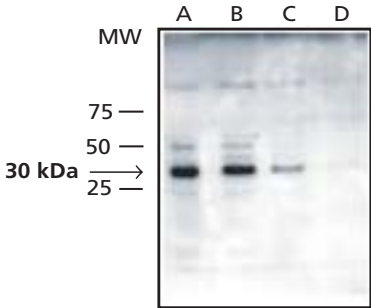

Immunoblot

Whole cell extract of rat liver mitochondria was separated on SDS-PAGE and blotted with Rabbit Anti-VDAC1/Porin (Product No. V 2139).

The antibody was developed with Goat Anti-Rabbit IgG Peroxidase (Product No. A 0545) and chemiluminescent substrate.

Lane A: Antibody dilution 1:2,000

Lane B: Antibody dilution 1:3,000

Lane C: Antibody dilution 1:4,000

Lane D: Antibody dilution 1:2,000
+ 20 µg/ml VDAC1/Porin peptide