

3050 Spruce Street, St. Louis, MO 63103 USA Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757 email: techservice@sial.com sigma-aldrich.com

Product Information

89740 Atto 390 maleimide

Application

Atto 390 is a novel fluorescent label with a coumarin structure. The dye is intended for application in the area of life science, e.g. labeling of DNA, RNA or proteins. Characteristic features of the label are high fluorescence quantum yield, large Stokes-shift, good photostability and low molecular weight.

Maleimides are well suited for coupling to thiol groups. This is similar to iodacetamides, but maleimides do react more thiol selective. They do not show significant reaction with histidine or methionine. Hydrolysis of maleimides to a mixture of isomeric nonreactive maleamic acids can compete significantly with thiol modification, particularly above pH 8. Maleimides may be used for labelling of amines, which usually requires a higher pH than reaction of maleimides with thiols.


After coupling to a substrate the label carries a net electrical charge of -1.

Product Description

$\begin{array}{l} \text{MW} \\ \lambda_{abs} \\ \epsilon_{max} \\ \lambda_{fl} \\ \eta_{fl} \\ \tau_{fl} \\ \text{CF}_{260} \\ \end{array}$	466 g/mol 390 nm 2.4 x 10 ⁴ M ⁻¹ cm ⁻¹ 479 nm 90 % 5.0 ns 0.52
CF ₂₆₀	0.52
CF ₂₈₀	0.08



Optical data of the carboxy derivative (in water)

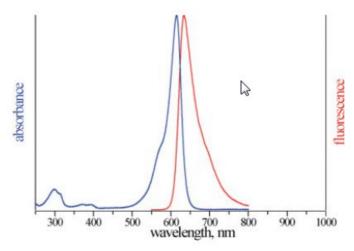

Storage: Store at -20 ℃ and protected from light.

General procedure for labelling proteins with maleimides

- 1) Dissolve the protein at $50-100~\mu\text{M}$ in a suitable buffer at pH 7.0-7.5 at room temperature. Common buffers include 10-100~mM phosphate, Tris, HEPES. Under those conditions, the protein thiol groups are sufficiently nucleophilic so that they react almost exclusively with the reagent. Other protein amines mostly remain protonated and relatively unreactive.
- 2) Reduce disulfide bonds in the protein. A 10-fold molar excess of a reducing agent such as DTT (43817) or TCEP (93284) is usually sufficient. If DTT is used, then dialysis is required to remove the excess DTT and TCEP prior to introducing the reactive dye.
- 3) As thiols can be oxidized to disulfides, It may be advisable to carry out thiol modifications in an oxygen-free environment. This is particularly important if the protein has been treated with a reagent such as dithiothreitol prior to thiol modification. In this case, all buffers should be deoxygenated and the reactions carried out under an inert atmosphere to prevent reformation of disulfides.
- **4)** Prepare a 10–20 mM stock solution of the reactive dye in a suitable solvent immediately prior to use (DMSO is the most common choice). Protect all stock solutions from light as much as possible by wrapping containers in aluminum foil.
- **5)** Add sufficient protein-modification reagent from a stock solution to achieve an 10–20 molar excess compared to protein. Add the reagent dropwise to the protein solution as it is stirring.
- 6) Let the reaction proceed for 2 hours at room temperature or overnight at 4 ℃. In both cases reaction should take place in the dark.
- 7) Upon completion of the reaction with the protein, an excess soluble low molecular weight thiol (e.g. glutathione, mercaptoethanol) can be added to consume excess thiol-reactive reagent, thus ensuring that no reactive species are present during the purification step.
- 8) Separate the conjugate on a gel filtration column, such as a Sephadex G-25 column or equivalent matrix, or by extensive dialysis at 4°C in an appropriate buffer.

Precautions and Disclaimer

This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

©2013 Sigma-Aldrich Co. LLC. All rights reserved. SIGMA-ALDRICH is a trademark of Sigma-Aldrich Co. LLC, registered in the US and other countries. Sigma brand products are sold through Sigma-Aldrich, Inc. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see product information on the Sigma-Aldrich website at www.sigmaaldrich.com and/or on the reverse side of the invoice or packing slip.