





## Product Information

### SIGMA QUALITY CONTROL TEST PROCEDURE

#### Enzymatic Assay of CELLULASE (EC 3.2.1.4)

##### PRINCIPLE:

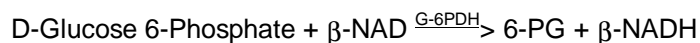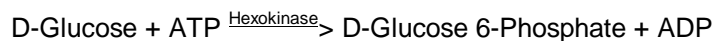

##### Abbreviations used:

ATP = Adenosine 5'-Triphosphate

ADP = Adenosine 5'-Diphosphate

$\beta$ -NAD =  $\beta$ -Nicotinamide Adenine Dinucleotide, Oxidized Form

$\beta$ -NADH =  $\beta$ -Nicotinamide Adenine Dinucleotide, Reduced Form

G-6PDH = Glucose 6-Phosphate Dehydrogenase

G-PG = 6-Phospho-D-Gluconate

**CONDITIONS:** T = 37°C, pH = 5.0, A<sub>340nm</sub>, Light path = 1 cm

**METHOD:** Spectrophotometric Stop Rate Determination

##### REAGENTS:

- A. 50 mM Sodium Acetate Buffer, pH 5.0 at 37°C  
(Prepare 200 ml in deionized water using Sodium Acetate, Trihydrate, Sigma Prod. No. S-8625. Adjust to pH 5.0 at 37°C with 1 M HCl.)
- B. 5% (w/v) Sigmacell Solution (Sigmacell)  
(Prepare 100 ml in Reagent A using Cellulose (Sigmacell), Type 20, Sigma Prod. No. S-3504. Mix and heat gently to make a uniform suspension.)
- C. Cellulase Enzyme Solution (Cellulase)  
(Immediately before use, prepare a solution containing 2 - 6 units/ml of Cellulase in cold deionized water.)
- D. Glucose (HK) Determination Vial (16-10)  
(Use Glucose (HK) 10, Sigma Stock No. 16-10. Dissolve the contents in 10 ml of deionized water.)

**Enzymatic Assay of CELLULASE  
(EC 3.2.1.4)**

**PROCEDURE:**

Step 1:

Pipette (in milliliters) the following reagents into suitable tubes:

	<u>Test</u>	<u>Blank</u>
Reagent B (Sigmacell)	4.00	4.00

Equilibrate to 37°C. Then add:

Reagent C (Cellulase)	1.00	-----
Deionized Water	-----	1.00

Immediately mix by swirling and incubate at 37°C for exactly 120 minutes with moderate shaking.

Immediately transfer into an ice bath. Allow to stand until the suspension is settled. Centrifuge for 2 minutes to clarify. Use the supernatant in Step 2.

Step 2:

Pipette (in milliliters) the following reagents into suitable cuvettes:

	<u>Test</u>	<u>Blank</u>
Reagent D (16-10)	3.00	3.00

Equilibrate to 25°C. Monitor the  $A_{340nm}$  until constant, using a suitably thermostatted spectrophotometer. Record the initial  $A_{340nm}$  for both the Test and Blank. Then add:

Test Supernatant (Step 1)	0.10	-----
Blank Supernatant (Step 1)	-----	0.10

Immediately mix by inversion and record the increase in  $A_{340nm}$  until complete (for approximately 5 minutes). Obtain the final  $A_{340nm}$  for both the Test and Blank.

**CALCULATIONS:**

$$\Delta A_{340nm} \text{ Test} = A_{340nm} \text{ Test Final} - A_{340nm} \text{ Test Initial}$$

$$\Delta A_{340nm} \text{ Blank} = A_{340nm} \text{ Blank Final} - A_{340nm} \text{ Blank Initial}$$

**Enzymatic Assay of CELLULASE  
(EC 3.2.1.4)**

**CALCULATIONS:** (continued)

$$\text{Units/ml enzyme} = \frac{(\Delta A_{340\text{nm}} \text{ Test} - \Delta A_{340\text{nm}} \text{ Blank})(3.1)(5)(\text{df})}{(6.22)(2)(1)(0.1)}$$

3.1 = Final volume (in milliliters) of Step 2

5 = Total volume (in milliliters) of reaction mix (Step 1)

df = Dilution factor

6.22 = Millimolar extinction coefficient of  $\beta$ -NADH at 340nm

2 = Conversion factor from 2 hours to 1 hour as per the Unit Definition

1 = Volume (in milliliter) of cellulase used in Step 1

0.1 = Volume (in milliliter) from Step 1 used in Step 2

$$\text{Units/mg solid} = \frac{\text{units/ml enzyme}}{\text{mg solid/ml enzyme}}$$

**UNIT DEFINITION:**

One unit will liberate 1.0  $\mu$ mole of glucose from cellulose in one hour at pH 5.0 at 37°C (2 hour incubation time).

**FINAL ASSAY CONCENTRATION:**

In a 5.00 ml reaction mix, the final concentrations are 40 mM sodium acetate, 4% (w/v) Sigmacell and 2 - 6 units of cellulase.

**REFERENCE:**

Worthington, C.E. (1988) *Worthington Enzyme Manual*, pp. 76-79, Worthington Biochemical Corporation, Freehold, NJ

**NOTES:**

1. This assay is based on the cited reference.
2. Where Sigma Product or Stock numbers are specified, equivalent reagents may be substituted.

**Sigma warrants that the above procedure information is currently utilized at Sigma and that Sigma products conform to the information in Sigma publications. Purchaser must determine the suitability of the information and products for its particular use. Upon purchase of Sigma products, see reverse side of invoice or packing slip for additional terms and conditions of sale.**