

Product Information

Cannabidiol solution

1.0 mg/mL±5% in methanol
drug standard


Catalog Number **C6395**

Lot Number 090M8725

Storage Temperature 2–8 °C

CAS RN 13956-29-1

Synonyms: CBD, 2-[3-Methyl-6-(1-methylethenyl)-
2-cyclohexen-1yl]-5-pentyl-1,3-benzenediol

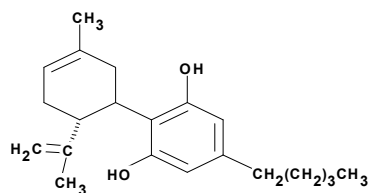

Product Description

Molecular Formula: C₂₁H₃₀O₂

Molecular Weight: 314.46

DEA Class: Exempt preparation of Class I

Concentration: 1.01 mg of cannabidiol/mL of GC grade methanol, concentration verified by UV.

Expiration Date: September 2013

Precautions and Disclaimer


This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

Store the product at 2–8 °C. If stored at 2–8 °C in an airtight container, protected from light, decomposition of the product is less than 1% in 48 months. After opening, concentration may change due to loss of solvent.


CMW,MAM 09/10-1

70 eV Electron Ionization Mass Spectrum

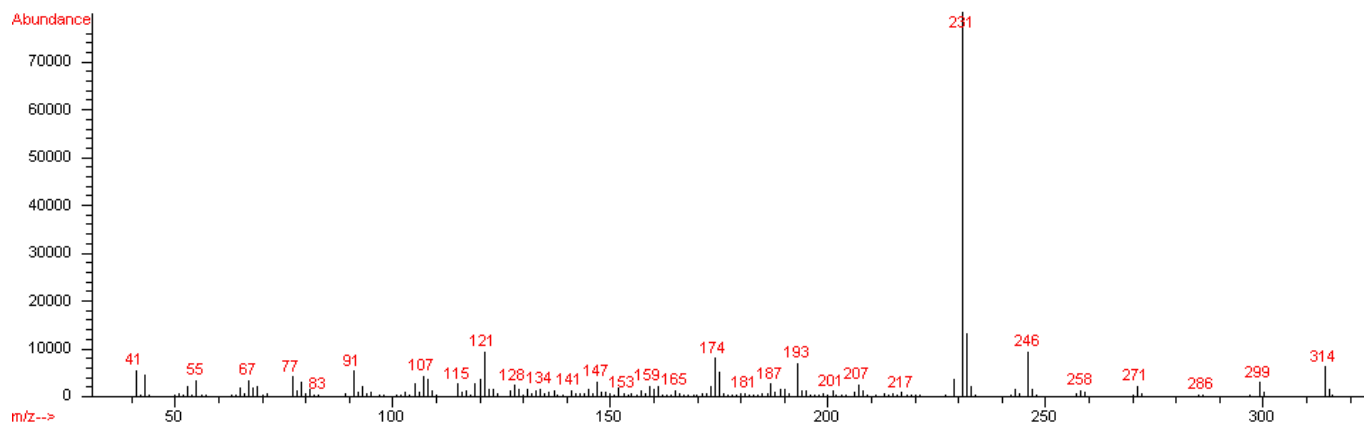

Hewlett-Packard 5890 Series II Gas Chromatograph
Hewlett-Packard 5971 Mass Selective Detector
Hewlett-Packard G1034C MS ChemStation (HP-DOS series)


FT-IR

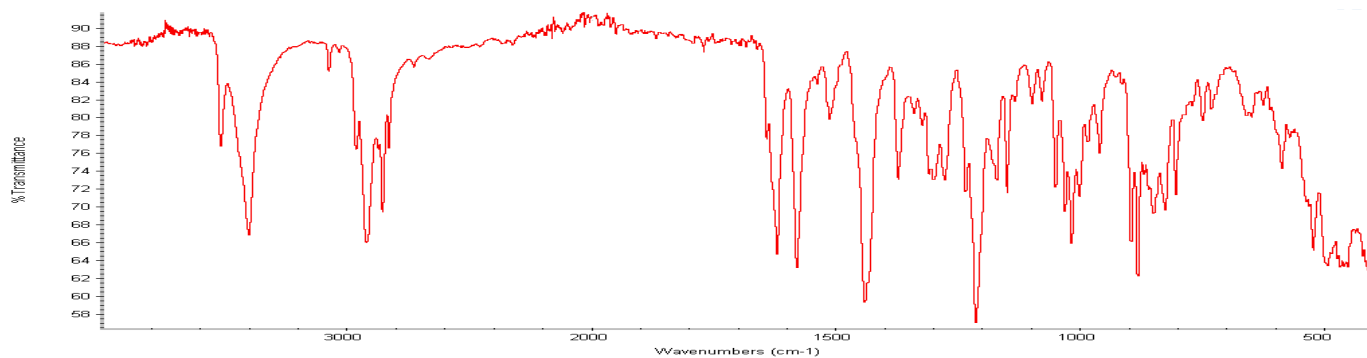

Nicolet 380 FT-IR Thermo Smart Orbit

Ultraviolet Spectrum


Peak	273.5 nm
Absorbance	0.4012
Solvent	GC grade methanol
Dilution	1:9 (v/v)

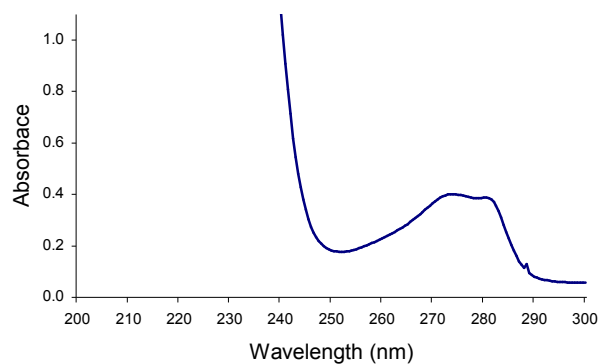

Beckman-Coulter DU 800

Capillary Gas Chromatography

Column Supelco SPB-1 (Catalog Number 24162)
I.D. 0.20 mm **Length** 15 m **Film Thickness** 0.20 μ m
Oven Temperature 150 $^{\circ}$ C (2 minutes), then
15 $^{\circ}$ C/minute to 300 $^{\circ}$ C, hold 11 minutes
Injector Temperature 250 $^{\circ}$ C, Split
Detector Hewlett-Packard 5971
Mass Selective Detector (70 eV, 40-500 m/z)

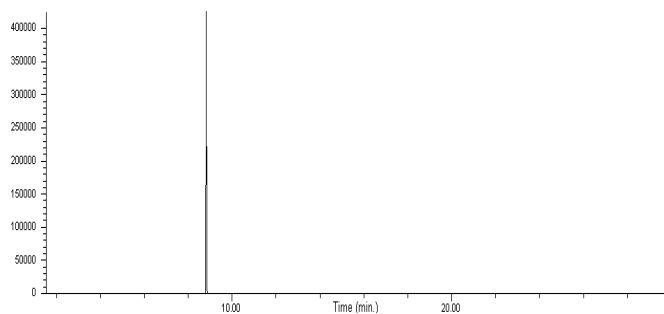

Hewlett-Packard 5890

Sigma brand products are sold through Sigma-Aldrich, Inc.

Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip.