


PUREGRIP® Carboys www.foxxlifesciences.com

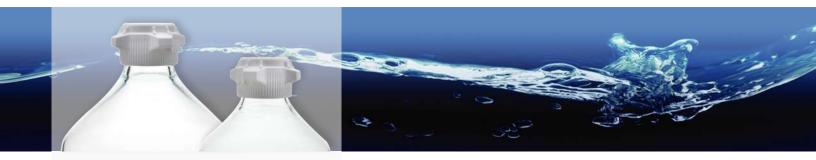

- World's first glass Carboy with 83B and 120mm caps
- · Sizes range from 5 L 50 L
- Large, visible writing surface on bottle for additional labeling; white graduations for identification of approximate volume; retrace code for identification
- Borosil® 3.3 Borosilicate Glass is mechanically strong, chemically resistant, and withstands thermal shock
- Bottle, VersaCap®, and Pouring Ring are autoclavable and comply with USP-VI standards

FOXX
PUREGRIP™
GLASS CARBOY
10000 mL

Standard Company

GLASS CARBOY
10000 mL



PUREGRIP® is a registered trademark of Foxx Life Sciences.

The PUREGRIP® Glass Carboy encompasses the latest offering in Foxx's collaboration with Borosil®. Combining the quality of Borosil® glassware with Foxx's innovative, user-focused designs, PUREGRIP® Carboys are the first ever 83mm Glass Carboy on the market. This product line features the largest writing surfaces and visible graduations on a glass carboy, while including our patented VersaCap® made from Class VI, Animal Free, Polypropylene.

A platinum-cured silicone O-ring ensures our leak-proof seal while the large handles with finger grips allow for easy tightening and opening. Our Open VersaCap® technology, sold separately, is also compatible with a variety of adapter inserts, enabling users to quickly change top connections to match their application needs. This product line is autoclavable and sterilizable.

Frosted surface on top of cap for convenient labeling

Patented VersaCap® is leak-proof, has replaceable platinum-cured silicone gaskets, and provides the ultimate in top-connection flexibility

3 USP Class VI, Animal-Free Polypropylene VersaCap® and Platinum Cured Silicone Gasket for leakproof, drip-free experience

4 Sloped shoulders for smooth, natural pouring motion

5 Largest enamel in the industry for additional labeling; white graduations for identification of approximate volume; retrace code for identification

6 Easy to read graduation marks





PUREGRIP® Carboys www.foxxlifesciences.com



	Part Number	Description	Material	Cap Size (mm)	Autoclavable	Qty/Pack
PUREGRIP® Carboys	1551033-FLS	PUREGRIP® Glass Carboys, Round, Clear, 5L	BORO	83	Yes	1
	1551038-FLS	PUREGRIP® Glass Carboys, Round, Clear, 10L	BORO	83	Yes	1
	1551040-FLS	PUREGRIP® Glass Carboys, Round, Clear, 20L	BORO	120	Yes	1
	15R-E121-FLS	PUREGRIP® Glass Carboys, Round, Clear, 50L	BORO	120	Yes	1
Versacap [®] 83B	205-1003-RLS	Versacap®, Closed	PP	83	Yes	1
	205-1004-RLS	Versacap®, Open with closed adapter	PP	83	Yes	1
	205-1007-RLS	Versacap®, Open with adapter, 1/4" HB x 2 and Vent	PP	83	Yes	1
	205-1008-RLS	Versacap®, Open with adapter, 1/2" HB x 2 and Vent	PP	83	Yes	1
	205-3001-RLS	Adapter, Closed	PP	83	Yes	1
	205-3002-RLS	Adapter, Quick Connect, 1/4" HB x 2	PP	83	Yes	1
	205-3003-RLS	Adapter, Quick Connect, 3/8" HB x 2	PP	83	Yes	1
	205-3004-RLS	Adapter, Quick Connect, 1/4" HB x 3	PP	83	Yes	1
	205-3005-RLS	Adapter, Quick Connect, 3/8" HB x 3	PP	83	Yes	1
	205-3006-RLS	Adapter, Quick Connect, 1/4" HB x 4	PP	83	Yes	1
	205-3007-RLS	Adapter, Quick Connect, 3/8" HB x 4	PP	83	Yes	1
	205-5001-RLS	Adapter, Molded, 1/4" HB x 2 and Vent	PP	83	Yes	1
	205-5002-RLS	Adapter, Molded, 1/2" HB x 2 and Vent	PP	83	Yes	1
	205-3001-RLS	Adapter, Closed Clearview®	PC	83	Yes	1

^{*} EZBio® Cap Assemblies are compatible with the 83B VersaCap® and can be purchased separately.

Rev. 08/19/2022 Foxx Life Sciences,12 N Wentworth Ave, Londonderry, NH 03053, USA www.foxxlifesciences.com Order Now: 1 (603) 890-FOXX (3699)
PN: DS00013-02 Foxx Life Sciences India, SN-65, Rajabolaram Village, Medchal, Telangana, India, 501401 www.foxxlifesciences.in Order Now: +91 (8418)- 661210 / 15

