


Catalogue Number:
NA##/371-221 and
NU###/371-221

Specification Sheet

NovAseptic® Valve, Valve Body, Shut off 90°, TC Connectable

Product Description

The NovAseptic® valve is specifically designed for aseptic applications and complies with the most stringent cleanability and sterilization requirements. The focus on aseptic design is a significant feature for all valves. The NovAseptic® valve is designed for minimum dead leg, complete drainability and high chemical resistance.

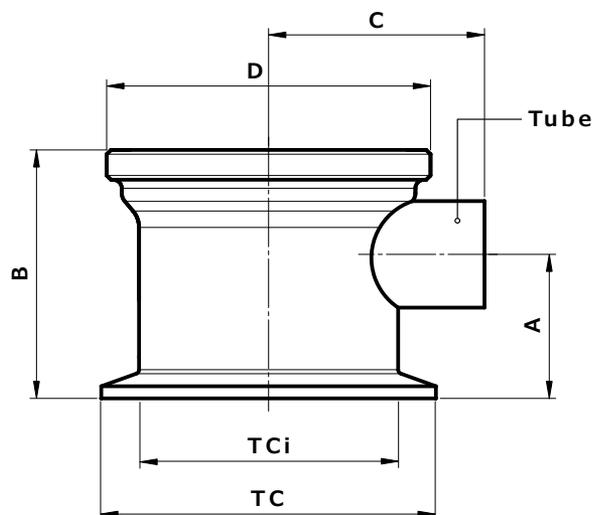

Catalogue No. Structure

NA##/371-221

NU###/371-221

= NA tube size


= NU tube size

Nominal Dimensions in mm (in.)

Catalogue No.	A	B	C	D*	TCi	TC	SMS Tube	US Tube ASTM®
							3008	A269/270
							Tube Size	Tube Size
NA25/371-221							25.0 x 1.2	
NU100/371-221	33.9 (1.335)	58.3 (2.295)	50.0 (1.969)	TR75x2	60.5 (2.382)	77.5 (3.051)		25.4 x 1.65 (1.0 x 0.065)
NA38/371-221							38.0 x 1.2	
NU150/371-221	43.8 (1.724)	75.0 (2.953)	60.0 (2.362)	TR85x2	72.0 (2.835)	91.0 (3.583)		38.1 x 1.65 (1.5 x 0.065)

* Non-standard thread, special thread to fit NovAseptic® Actuators.

Recommended Orientation for full drain ability of NovAseptic® Diaphragm Valves

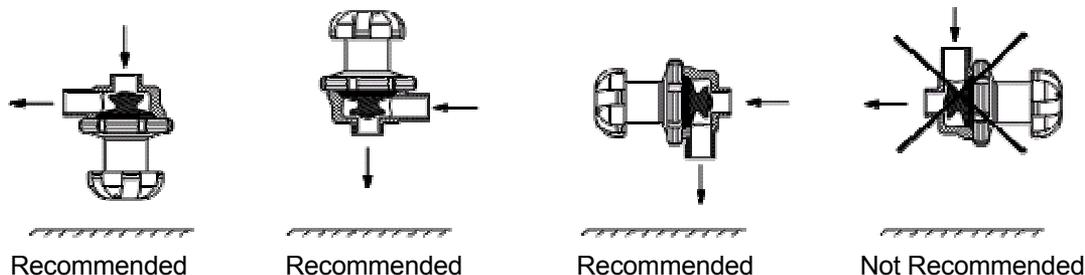

Specifications

Net Volume (valve cavity, diaphragm applied)		
Valve Body	NA25/NU100	NA38/NU150
Volume mL	30	80

Net Weight (approximate)		
Valve Body	NA25/NU100	NA38/NU150
Weight kg (lb)	0.6 (1.3)	1.2 (2.65)

Material				
	Bar		Tubing	
	Stainless Steel in Compliance with			
Material Code	316L	EN 1.4404	TP316L	EN 1.4404
Technical Requirements	ASME® SA-479	EN 10272	ASTM® A269/270	EN 10217-7/EN 10216-5

Specifications	
Surface Roughness	Internal surface (electropolished) $Ra \leq 0.38 \mu\text{m}$ (15 μin) External surface $Ra \leq 1.6 \mu\text{m}$ (63 μin)
Design Temperature, Valve Body	-80 to 200 °C (-112 to 392 °F)
Design Pressure, Valve Body	-1 to 10 bar(g) (-14.5 to 145 psi(g))
Note	The valve body can only be used with NovAseptic® diaphragm, actuator and position indicator (optional). The applied diaphragm and actuator have different design temperature and/or pressure limits. The weakest component in the assembled product determines the maximum operating/design temperature and pressure limits.
Labeling	Each valve body is individually labeled for full traceability and heat No. according to our QA routines.
Packaging	All open ends on the valve body are protected with plastic covers on delivery and packaged in a closed box.
Quality Control	Our quality assurance system guarantees the control and traceability at all stages of the manufacturing.
Regulatory	<ul style="list-style-type: none"> The NovAseptic® valves, used with gas group 2 and liquid group 1, fall under the art. 4 § 3 of the Pressure Equipment Directive (2014/68/EU). They are not CE marked. They have been designed and manufactured following the sound engineering practices. NovAseptic® valves follow the applicable ASME BPE guideline sections.
Options	For non-standard NovAseptic® valve body options, please contact us for further information.

Technical Assistance

For more information, please visit [SigmaAldrich.com](https://www.sigmaaldrich.com) for up-to-date worldwide contact information

