

When downstream quality matters, make sure your upstream tools are the best.

Amicon® Ultra 2 mL Centrifugal Filters for Protein Purification and Concentration

The Amicon® Ultra 2 mL filter provides a fast and easy way to concentrate, desalt and exchange the buffer for both proteins and nucleic acids. Now you can concentrate 2 mL down to 15 – 70 μ L in an Amicon® Ultra centrifugal device.

By incorporating familiar Amicon® Ultra features such as the vertical membrane design and true dead stop, the Amicon® Ultra 2 mL filter delivers high recovery, short spin times and quality results even for the most dilute samples.

Amicon® Ultra 2 mL Advantages

- True engineered dead stop provides consistent recovery without spinning to dryness
 - Reverse spin aids in sample recovery, especially with small, dilute samples
 - Vertical membrane design aligns with the filtrate for less clogging and faster filtering
 - Ultracel® low-binding membrane combines with the innovative vertical housing to provide fast sample processing without sacrificing recovery
 - Fast processing, with 10- to 30-minute spin times
-




Easy to Choose!

Select the Amicon® Ultra Filter that best meets your sample needs. Choose your starting volume, identify the molecular weight of your molecule of interest, and identify final volume expected.

Amicon® Ultra Centrifugal Filters

	Product	Amicon® Ultra-0.5	Amicon® Ultra-2	Amicon® Ultra-4	Amicon® Ultra-15
	Maximum initial sample volume (mL)	0.5	2	4	15
	Final concentrate (retentate) volume (µL)	15 – 20	15 – 70	30 – 70	150 – 300
MWCO	Qty/Pk				
3,000 MWCO	8	UFC500308	UFC200324	UFC800308	UFC900308
	24	UFC500324		UFC800324	UFC900324
	96	UFC500396		UFC800396	UFC900396
	500	UFC5003BK			
10,000 MWCO	8	UFC501008	UFC201024	UFC801008	UFC901008
	24	UFC501024		UFC801024	UFC901024
	96	UFC501096		UFC801096	UFC901096
	500	UFC5010BK			
30,000 MWCO	8	UFC503008	UFC203024	UFC803008	UFC903008
	24	UFC503024		UFC803024	UFC903024
	96	UFC503096		UFC803096	UFC903096
	500	UFC5030BK			
50,000 MWCO	8	UFC505008	UFC205024	UFC805008	UFC905008
	24	UFC505024		UFC805024	UFC905024
	96	UFC505096		UFC805096	UFC905096
	500	UFC5050BK			
100,000 MWCO	8	UFC510008	UFC210024	UFC810008	UFC910008
	24	UFC510024		UFC810024	UFC910024
	96	UFC510096		UFC810096	UFC910096
	500	UFC5100BK			

Search with our online Amicon® Ultra Selector Tool and find the right filter to concentrate your sample. Visit: www.millipore.com/FastEasy

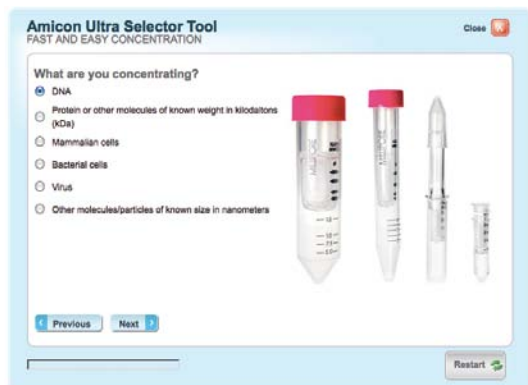

To Place an Order or Receive Technical Assistance

In Europe, please call Customer Service:

France: 0825 045 645

Germany: 01805 045 645

Italy: 848 845 645

Spain: 901 516 645 Option 1

Switzerland: 0848 645 645

United Kingdom: 0870 900 4645

For other countries across Europe, please call: +44 (0) 115 943 0840

Or visit: www.merckmillipore.com/offices

For Technical Service visit:

www.merckmillipore.com/techservice



www.merckmillipore.com

Merck Millipore and the M logo are trademarks of Merck KGaA, Darmstadt, Germany. Amicon and Ultracel are registered trademarks of Merck KGaA, Darmstadt, Germany. Lit. No. DS3336ENEU 3/2012 LS-SBU-12-06120 Printed in the USA. © 2012 EMD Millipore Corporation, Billerica, MA, USA. All rights reserved.

Get Connected!

Join Merck Millipore Bioscience on your favorite social media outlet for the latest updates, news, products, innovations, and contests!

 facebook.com/MerckMilliporeBioscience

 twitter.com/Merck4Bio