




SIGMA QUALITY CONTROL TEST PROCEDURE

Product Information

Enzymatic Assay of ARGININE KINASE (EC 2.7.3.3)

PRINCIPLE:

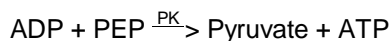

Abbreviations used:

ATP = Adenosine 5'-Triphosphate

AK = Arginine Kinase

ADP = Adenosine 5'-Diphosphate

PEP = Phospho(enol)pyruvate

PK = Pyruvate Kinase

β -NADH = β -Nicotinamide Adenine Dinucleotide, Reduced Form

LDH = L-Lactic Dehydrogenase

β -NAD = β -Nicotinamide Adenine Dinucleotide, Oxidized Form

CONDITIONS: T = 30°C, pH = 8.6, $A_{340\text{nm}}$, Light path = 1 cm

METHOD: Continuous Spectrophotometric Rate Determination

REAGENTS:

- A. 250 mM Glycine Buffer, pH 8.6 at 30°C
(Prepare 100 ml in deionized water using Glycine, Free Base, Sigma Prod. No. G-7126.
Adjust to pH 8.6 at 30°C with 1 M NaOH.)
- B. 100 mM Glycine buffer, with 10 mM 2-Mercaptoethanol (Enz Dil)
(Prepare 50 ml in deionized water using Reagent A and 2-Mercaptoethanol, Sigma Prod. No. M-6250.)
- C. 200 mM Magnesium Sulfate Solution (MgSO_4)
(Prepare 10 ml in deionized water using Magnesium Sulfate, Heptahydrate, Sigma Prod. No. M-1880.)
- D. 2 M Potassium Chloride Solution (KCl)
(Prepare 10 ml in deionized water using Potassium Chloride, Sigma Prod. No. P-4504.)

**Enzymatic Assay of ARGININE KINASE
(EC 2.7.3.3)**

REAGENTS: (continued)

- E. 300 mM Phospho(enol)pyruvate Solution (PEP)
(Prepare 3 ml in deionized water using Phospho(enol)pyruvate, Trisodium Salt, Hydrate, Sigma Prod. No. P-7002.)
- F. 200 mM Adenosine 5'-Triphosphate Solution (ATP)
(Prepare 2 ml in deionized water using Adenosine 5'-Triphosphate, Disodium Salt, Sigma Prod. No. A-5394.)
- G. 500 mM L-Arginine Solution (ARG)
(Prepare 2 ml in deionized water using L-Arginine, Free Base, Sigma Prod. No. A-5006.)
- H. 7.5 mM β -Nicotinamide Adenine Dinucleotide, Reduced Form, Solution (β -NADH)
(Dissolve the contents of one 5 mg vial of β -Nicotinamide Adenine Dinucleotide, Reduced Form, Disodium Salt, Sigma Stock No. 340-105, in the appropriate volume of Reagent A. **PREPARE FRESH.**)
- I. PK/LDH Enzyme Solution (PK/LDH)
(Immediately before use, prepare a solution containing 20 units/ml of Pyruvate Kinase in cold Reagent A using PK/LDH Enzymes Suspension, Sigma Stock No. 40-7.¹)
- J. Arginine Kinase Enzyme Solution (AK)
(Prepare a solution containing 0.2 - 0.4 unit/ml of Arginine Kinase in Reagent B. Let stand at room temperature for 15 minutes before use to obtain full activation.)

PROCEDURE:

Prepare a reaction cocktail by pipetting (in milliliters) the following reagents into a suitable container:

Reagent A (Buffer)	19.50
Reagent C (MgSO ₄)	2.00
Reagent D (KCl)	2.00
Reagent E (PEP)	2.00
Reagent F (ATP)	1.00
Reagent H (β -NADH)	0.50

Mix by swirling and adjust to pH 8.6 at 30°C with 1 M HCl or 1 M NaOH.

**Enzymatic Assay of ARGININE KINASE
(EC 2.7.3.3)**

PROCEDURE: (continued)

Pipette (in milliliters) the following reagents into suitable cuvettes.

	<u>Test</u>	<u>Blank</u>
Reaction Cocktail	2.70	2.70
Reagent I (PK/LDH)	0.10	0.10
Reagent G (ARG)	0.10	0.10

Mix by inversion and equilibrate to 30°C. Monitor the $A_{340\text{nm}}$ until constant, using a suitably thermostatted spectrophotometer. Then add:

Reagent J (AK)	0.10	-----
Reagent B (Enz Dil)	-----	0.10

Immediately mix by inversion and record the decrease in the $A_{340\text{nm}}$ for approximately 5 minutes. Obtain the $\Delta A_{340\text{nm}}$ /minute using the maximum linear rate for both the Test and Blank.

CALCULATIONS:

$$\text{Units/ml enzyme} = \frac{(\Delta A_{340\text{nm}}/\text{min Test} - \Delta A_{340\text{nm}}/\text{min Blank})(3)(\text{df})}{(6.22)(0.1)}$$

3 = Total volume (in milliliters) of assay

df = Dilution factor

6.22 = Millimolar extinction coefficient of β -NADH at 340 nm

0.1 = Volume (in milliliter) of enzyme used

$$\text{Units/mg solid} = \frac{\text{units/ml enzyme}}{\text{mg solid/ml enzyme}}$$

$$\text{Units/mg protein} = \frac{\text{units/ml enzyme}}{\text{mg protein/ml enzyme}}$$

UNIT DEFINITION:

One unit will convert 1.0 μmole of L-arginine and ATP to N^{ω} -phospho-L-arginine and ADP per minute at pH 8.6 at 30°C.

**Enzymatic Assay of ARGININE KINASE
(EC 2.7.3.3)**

FINAL ASSAY CONCENTRATIONS:

In a 3.00 ml reaction mix, the final concentrations are 178 mM glycine, 0.33 mM 2-mercaptoethanol, 13 mM magnesium sulfate, 133 mM potassium chloride, 20 mM phospho(enol) pyruvate, 6.7 mM adenosine 5'-triphosphate, 0.13 mM β -nicotinamide adenine dinucleotide, reduced form, 2 units pyruvate kinase, 3 units lactic dehydrogenase, 17 mM L-arginine, and 0.02 - 0.04 unit arginine kinase.

REFERENCE:

Blethen, S. (1970) *Methods in Enzymology*, XVIIA, 330-335

NOTES:

1. Contains not less than 700 Pyruvate Kinase units and 1000 L-Lactic Dehydrogenase units per ml.
2. Pyruvate Kinase Unit Definition: One unit will convert 1.0 μ mole of phospho(enol)pyruvate to pyruvate per minute at pH 7.6 at 37°C.
3. L-Lactic Dehydrogenase Unit Definition: One unit will reduce 1.0 μ mole of pyruvate to L-lactate per minute of pH 7.5 at 37°C.
4. This assay is based on the cited reference.
5. Where Sigma Product or Stock numbers are specified, equivalent reagents may be substituted.

Sigma warrants that the above procedure information is currently utilized at Sigma and that Sigma products conform to the information in Sigma publications. Purchaser must determine the suitability of the information and products for its particular use. Upon purchase of Sigma products, see reverse side of invoice or packing slip for additional terms and conditions of sale.