## SIGMA-ALDRICH®

#### sigma-aldrich.com

3050 Spruce Street, St. Louis, MO 63103 USA Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757 email: techservice@sial.com sigma-aldrich.com

# **Product Information**

MRCKα, Active human, recombinant GST-tagged, expressed in *Sf*9 cells

Catalog Number **M9447** Lot Number 019K1622 Storage Temperature –70 °C

Synonyms: CDC42BPA; MRCK; PK428; FLJ23347; KIAA0451; DKFZp686L1738; DKFZp686P1738

#### **Product Description**

Myotonic Dystrophy Kinase-Related Cdc42-Binding Kinase  $\alpha$  (MRCK $\alpha$ ) is a Cdc42/Rac/Rho interactive/ binding serine/threonine kinase with multiple functional domains.<sup>1</sup> MRCK are effectors of RhoA and Cdc42 for actin reorganization. MRCK $\alpha$  is a critical regulator of signal transduction pathways in eukaryotic cells that are known principally for their role in regulating the cytoskeleton and they do so by recruiting a variety of downstream effector proteins.<sup>2</sup>

This recombinant product was expressed by baculovirus in *Sf*9 insect cells using an N-terminal GST-tag. The gene accession number is NM 003607. It is supplied in 50 mM Tris-HCl, pH 7.5, with 150 mM NaCl, 0.25 mM DTT, 0.1 mM EGTA, 0.1 mM EDTA, 0.1 mM PMSF, and 25% glycerol.

Molecular mass: ~73 kDa

Purity: ≥70% (SDS-PAGE, see Figure 1)

Specific Activity: 86–116 nmole/min/mg (see Figure 2)

#### **Precautions and Disclaimer**


This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

#### Storage/Stability

The product ships on dry ice and storage at -70 °C is recommended. After opening, aliquot into smaller quantities and store at -70 °C. Avoid repeated handling and multiple freeze/thaw cycles.


### Figure 1.

SDS-PAGE Gel of Lot Number 019K1622: >75% (densitometry)

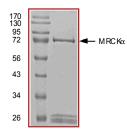

#### Figure 2.

Specific Activity of Lot Number 019K1622: 101 nmole/min/mg

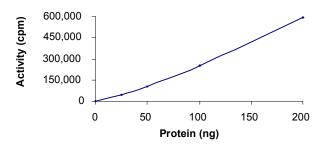

#### Procedure

**Preparation Instructions** 

Kinase Assay Buffer – 25 mM MOPS, pH 7.2, 12.5 mM glycerol 2-phosphate, 25 mM MgCl<sub>2</sub>, 5 mM EGTA, and 2 mM EDTA. Just prior to use, add DTT to a final concentration of 0.25 mM.

Kinase Dilution Buffer – Dilute the Kinase Assay Buffer 5-fold with water.

Kinase Solution – Dilute the Active MRCK $\alpha$  (0.1 µg/µl) with Kinase Dilution Buffer to the desired concentration. Note: The lot-specific specific activity plot may be used as a guideline (see Figure 2). It is recommended that the researcher perform a serial dilution of Active MRCK $\alpha$  kinase for optimal results.

10 mM ATP Stock Solution – Dissolve 55 mg of ATP in 10 ml of Kinase Assay Buffer. Store in 200  $\mu l$  aliquots at –20 °C.

 $\gamma$ -<sup>32</sup>P-ATP Assay Cocktail (250  $\mu$ M) – Combine 5.75 ml of Kinase Assay Buffer, 150  $\mu$ l of 10 mM ATP Stock Solution, 100  $\mu$ l of  $\gamma$ -<sup>32</sup>P-ATP (1 mCi/100  $\mu$ l). Store in 1 ml aliquots at –20 °C.

Substrate Solution – Dissolve the synthetic peptide substrate (KRRRLASLR) in water at a final concentration of 1 mg/ml.

1% phosphoric acid solution – Dilute 10 ml of concentrated phosphoric acid to a final volume of 1 L with water.

#### Kinase Assay

This assay involves the use of the <sup>32</sup>P radioisotope. All institutional guidelines regarding the use of radioisotopes should be followed.

- 1. Thaw the Active MRCK $\alpha$ , Kinase Assay Buffer, Substrate Solution, and Kinase Dilution Buffer on ice. The  $\gamma$ -<sup>32</sup>P-ATP Assay Cocktail may be thawed at room temperature.
- 2. In a pre-cooled microcentrifuge tube, add the following solutions to a volume of 20 μl:
  - 10 µl of Kinase Solution
  - 5 µl of Substrate Solution
  - $5 \mu$  of cold water (4 °C)
- Set up a blank control as outlined in step 2, substituting 5 μl of cold water (4 °C) for the Substrate Solution.
- 4. Initiate each reaction with the addition of 5  $\mu$ l of the  $\gamma$ -<sup>32</sup>P-ATP Assay Cocktail, bringing the final reaction volume to 25  $\mu$ l. Incubate the mixture in a water bath at 30 °C for 15 minutes.
- 5. After the 15 minute incubation, stop the reaction by spotting 20  $\mu$ l of the reaction mixture onto an individually precut strip of phosphocellulose P81 paper.

- 6. Air dry the precut P81 strip and sequentially wash in the 1% phosphoric acid solution with constant gentle stirring. It is recommended the strips be washed a total of 3 times of ~10 minutes each.
- 7. Set up a radioactive control to measure the total  $\gamma^{-32}$ P-ATP counts introduced into the reaction. Spot 5 µl of the  $\gamma^{-32}$ P-ATP Assay Cocktail on a precut P81 strip. Dry the sample for 2 minutes and read the counts. Do not wash this sample.
- 8. Count the radioactivity on the P81 paper in the presence of scintillation fluid in a scintillation counter.
- 9. Determine the corrected cpm by subtracting the blank control value (see step 3) from each sample and calculate the kinase specific activity

#### Calculations:

1. Specific Radioactivity (SR) of ATP (cpm/nmole)

SR =  $cpm of 5 \mu l of \gamma^{-32}P-ATP Assay Cocktail$ nmole of ATP cpm – value from control (step 7) nmole – 1.25 nmole (5  $\mu$ l of 250  $\mu$ M ATP Assay Cocktail)

2. Specific Kinase Activity (SA) (nmole/min/mg)

nmole/min/mg =  $\frac{\Delta \text{cpm} \times (25/20)}{\text{SR} \times \text{E} \times \text{T}}$ 

SR = specific radioactivity of the ATP (cpm/nmole ATP)  $\triangle$ cpm = cpm of the sample – cpm of the blank (step 3) 25 = total reaction volume

- 20 = spot volume
- T = reaction time (minutes)
- E = amount of enzyme (mg)

#### References

- 1. Ivan, T. et al., Genomic organization of human myotonic dystrophy kinase-related Cdc42-binding kinase  $\alpha$  reveals multiple alternative splicing and functional diversity. Gene, **304**, 107-115 (2003).
- Ivan, T. et al., Phosphorylation of a Novel Myosin Binding Subunit of Protein Phosphatase 1 Reveals a Conserved Mechanism in the Regulation of Actin Cytoskeleton. J. Biol. Chem., **276**, 21209-21216 (2001).

BKR,MAM 02/09-1

Sigma brand products are sold through Sigma-Aldrich, Inc.

Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip.