

06966 Atto 655-Biotin


Application

Atto 655 belongs to a new generation of fluorescent labels. The dye is designed for application in the area of life science, e.g. labeling of DNA, RNA or proteins. Characteristic features of the label are strong absorption, good fluorescence quantum yield, excellent thermal and photo-stability, outstanding ozone resistance, very good water solubility, and very little triplet formation. Atto 655 is a zwitterionic dye with a net electrical charge of zero. The fluorescence is efficiently quenched by electron donors like guanine, tryptophan, etc.


The **biotin** derivative can be used as reagent for binding to proteins like avidin and streptavidin.

Product Description

MW	1204 g/mol
λ_{abs}	663 nm
ϵ_{max}	$1.25 \times 10^5 \text{ M}^{-1} \text{ cm}^{-1}$
λ_{fl}	680 nm
η_{fl}	30 %
τ_{fl}	1.8 ns
CF ₂₆₀	0.24
CF ₂₈₀	0.08

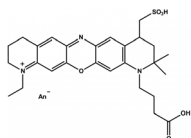

Optical data of the carboxy derivative (in aqueous solution)

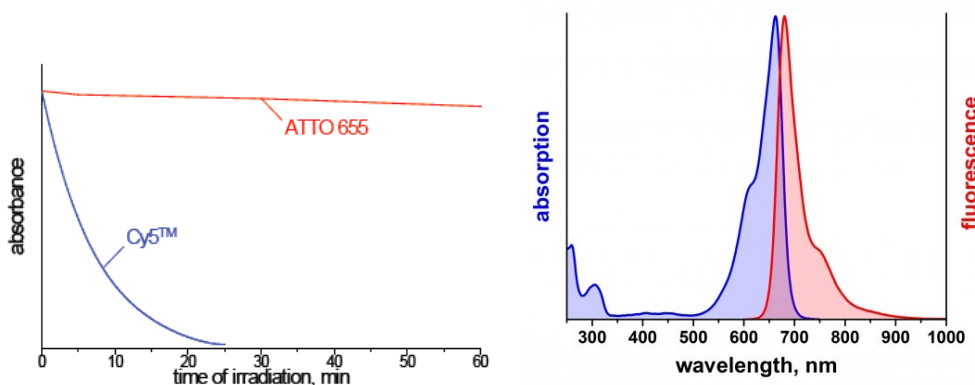

Storage: store at $\leq -20^\circ\text{C}$. Protect from long-term exposure to moisture and light.

Precautions and Disclaimer

This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

The vibrant M and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. Detailed information on trademarks is available via publicly accessible resources.
© 2018 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the US and Canada.

