

# Sodium Hydroxide Solution 0.1M - 58300C

# Liquid Solutions MyWay

58300C is part of Liquid Solutions MyWay – a portfolio of high-quality sterile filtered liquids in prefilled single-use bags that are in stock and ready-to-ship.

- Manufactured in an ISO 9001:2015 certified facility
- A 12-month shelf life from the date of production (with a target of 24 months)
- Manufactured using Emprove® raw materials (complying with Ph. Eur. and most with USP)
- Intended for bioprocessing and further manufacturing uses

#### Formulation 58300C - EMPROVE® raw materials


CAS No.	Description	Amount (per L)
1310-73-2	SODIUM HYDROXIDE SOLUTION 50%	5.281 ml
7732-18-5	WATER FOR INJECTION (USP/EP)	to 1 L

# **Outer packaging**

All packaging sizes will be shipped in a stainless steel drum.

# **Ship testing (ISTA)**

To safeguard delivery, all shipping containers were ISTA ship tested using a filled bag in an appropriate outer container. The 1 L and 20 L 2D bags were tested in accordance with ISTA 1A (International Safe Transit Association) that included a vibration and drop test as well as real-time shipping. The 190 L 3D bags were tested in accordance with ISTA 3E (for palletized loads).

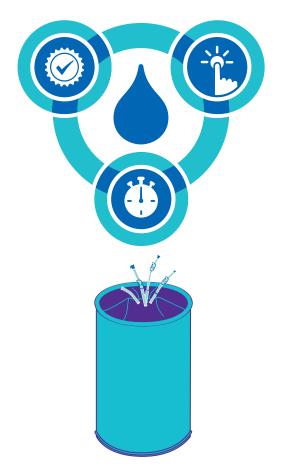

#### **Product**

58300C - Sodium hydroxide solution 0.1M		
58300C-1BC / -1BNC	10 x 1 L Bags	
58300C-20BC / -20BNC	20 L Bag	
58300C-190BC / -190BNC	190 L Bag	




#### 1 L Bag Description

58300C-1BC/	2	
-1BNC	Description	Material
<b>A</b> Bag	1 L Bag	ULDPE
<b>B</b> Filling Line	Filling Line, Sealed End	PVC/Silicone
C Connector 1	1/4 in. x 7/16 in. Tubing, 12 in. Length	TPE
	Coupling Insert 1/4 in. HB, MPC	Polycarbonate
	Sealing Cap, MPC	Polycarbonate
<b>D</b> Connector 2	1/4 in. x 7/16 in. Tubing, 12 in. Length	TPE
	1/4 in. HB x Female Luer Adapter	Polypropylene
	Male Luer Plug	Polycarbonate

# A B C C

#### 20 L Bag Description

58300C-20BC / -20BNC	Description	Material		
<b>A</b> Bag	20 L Bag	ULDPE		
<b>B</b> Filling Line	Filling Line, Sealed End	PVC/Silicone		
C Connector 1	3/8 in. x 5/8 in. Tubing, 24 in. Length Coupling Insert 3/8 in. HB, MPC Sealing Cap MPC	TPE Polycarbonate Polycarbonate		
<b>D</b> Connector 2	1/4 in. x 7/16 in. Tubing, 24 in. Length Coupling Body, 1/4 in. HB, MPC Sealing Plug, MPC	TPE Polypropylene Polycarbonate		
E Connector 3	1/4 in. x 3/8 in. Tubing, 3 in. Length 1/4 in. x 1/8 in. HB Reduction Fitting 1/8 in. x 1/4 in. Tubing, 24 in. Length 1/8 in. HB x Female Luer Adapter (Needleless Injection site for non-hazardous products) Male Luer Plug/Cap	TPE Polypropylene TPE Polycarbonate Polycarbonate Polypropylene/ Polycarbonate		
	Total Tubing Length: approx. 28 in.			

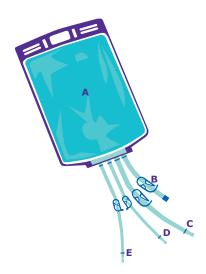

### 190 L Bag Description

58300C-190BC			
/ -190BNC	Description	Material	
<b>A</b> Bag	200 L Bag	ULDPE	
<b>B</b> Filling Line	Filling Line, Sealed End	PVC/Silicone	
C Connector 1	1/2 in. x 3/4 in. Tubing, 36 in. Length	TPE	
	Coupling Body 1/2 in. HB, MPX	Polycarbonate	
	Sealing Plug MPX	Polycarbonate	
<b>D</b> Connector 2	1/2 in. x 3/4 in. Tubing, 3 in. Length	TPE	
	1/2 in. x 3/8 in. HB Reduction Fitting	Polypropylene	
	3/8 in. x 5/8 in. Tubing, 36 in. Length	TPE	
	Coupling Insert, 3/8 in. HB, MPC	Polycarbonate	
	Sealing Cap, MPC	Polycarbonate	
	Total Tubing Length: approx. 40 in.		
E Connector 3	1/2 in. x 3/4 in. Tubing, 3 in. Length	TPE	
	1/2 in. x 1/4 in. HB Reduction Fitting	Polypropylene	
	1/4 in. x 7/16 in. Tubing, 36 in. Length	TPE	
	Coupling Body, 1/4 in. HB, MPC	Polycarbonate	
	Coupling Insert, 1/4 in. HB, MPC	Polycarbonate	
	1/4 in. x 7/16 in. Tubing, 3 in. Length	TPE	
	1/4 in. HB Swabable Valve / Needleless Injection Site	Polycarbonate/Silicone	
	Male Luer Cap	Polycarbonate	
	Total Tubing Length: approx. 43 in.		

